Silk'n®

GLASSPEDI

FTC05

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
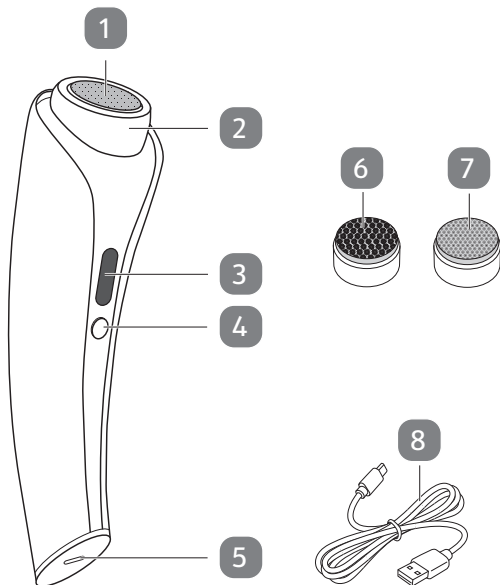
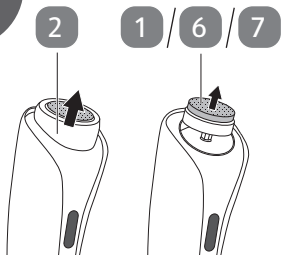
C

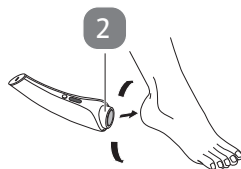Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silkn Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silkn Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silkn Beauty Ltd.

Invention Works B.V.

Galileïstraat 17, 3029 AL Rotterdam, The Netherlands
www.silkn.eu – info@silkn.eu

Table of contents

I	Package contents/device parts	5
2	General information	5
2.1	Reading and storing the user manual	5
2.2	Explanation of symbols	5
3	Safety	6
3.1	Proper use	6
3.2	Safety instructions	6
4	Checking the product and package contents	10
5	Charging	10
6	Operation	11
6.1	Attaching or replacing the grinding disc	11
6.2	Using the product	12
7	Cleaning	12
8	Storage	13
9	Troubleshooting	13
10	Technical data	13
11	Disposal	14
11.1	Disposing of the packaging	14
11.2	Disposing of the product	14
12	Warranty	15
13	Customer service	15

1 Package contents/device parts

- | | |
|-----------------------|------------------------|
| 1 Fine grinding disc* | 5 Charging port |
| 2 Grinding head | 6 Medium grinding disc |
| 3 Display | 7 Coarse grinding disc |
| 4 Operation button | 8 USB charging cable |

*Already attached to the grinding head.

2 General information

2.1 Reading and storing the user manual

This user manual accompanies this GlassPedi (hereafter referred to as the "product"), and contains important information on setup and handling.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product. The user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

WARNING!

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

NOTICE!

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

This symbol identifies voltage and the current rating for direct current.

3 Safety

3.1 Proper use

The product is exclusively designed for removing calluses on the feet. It is only intended for private use and not suitable for commercial purposes.

Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer cannot be held liable for damages or injury incurred through improper or incorrect use.

3.2 Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

- Only use the product with the supplied USB charging cable.
- The product must only be supplied at safety extra-low voltage corresponding to the marking on the product.
- Do not use the product if it is damaged or if the USB charging cable is defective.
- Do not open the housing; instead, have a qualified professional perform repairs. Contact the customer service team for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection or incorrect operation.
- Only parts that comply with the original

device data may be used for repairs. This product contains electrical and mechanical parts which are essential for providing protection against sources of danger.

- Keep the product and USB charging cable away from water or other liquids.
- Always switch the product off when you are not using it, when you intend to clean it, or in the event of a malfunction. Always unplug the USB charging cable when you are not charging the product.
- Never touch the USB charging cable with wet or damp hands.
- Never insert any objects into the housing.
- Keep all parts of the product away from open flames and hot surfaces.
- Lay the USB charging cable so that it does not pose a tripping hazard.
- Do not pull the USB charging cable out of the USB power source by the cord; instead, always pull it out by the cable plug.
- Do not kink the USB charging cable and do not lay it over sharp edges.
- Only use the product indoors.
- Do not use the product while charging.

WARNING! RISK OF INJURY!

- The product can be used by children aged from 8 years and above and persons with

reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the product in a safe way and understand the hazards involved. Children must not play with the product. Cleaning and user maintenance must not be performed by unsupervised children.

- Keep children under the age of 8 away from the product and USB cable.
- Make sure that children do not play with the plastic wrapping. They may get caught in it when playing and suffocate.
- Only use the product after the grinding disc is properly installed.
- Only use the product on your feet to treat calloused skin. Do not use it on other parts of the body and on soft or thin parts of the skin.
- Only use the product for dry treatment of skin (do not use any additives such as oils or lotions).
- Do not continuously treat one spot of your skin for more than 2 to 3 seconds at a time.
- Do not use the product directly on or near open wounds and inflamed areas of the skin or to remove warts.
- Immediately stop using the product if you experience pain or skin irritation.

- Consult a doctor if you are not sure whether the product is suitable for you (e.g. you suffer from skin diseases or diabetes, you are pregnant, or you are using a pacemaker).
- For hygienic reasons, the product should only be used by one person.
- Keep the product away from hair, wigs and other long flexible objects, e.g. shoelaces, threads, textile fibres.
- Under abusive conditions, liquid may be ejected from the battery; avoid contact. Liquid ejected from the battery may cause irritation or burns.
- If battery fluid comes into contact with skin, eyes or mucous membrane, wash and rinse the affected areas immediately with clean water. Consult a physician if necessary.
- The battery must not be disassembled, short-circuited or thrown into a fire.
- Do not try to remove the battery from the product. This product contains batteries that are non-replaceable.

NOTICE! RISK OF DAMAGE!

- Keep the product away from open fire, hot surfaces, high temperatures or the effects of the weather (e.g. rain).
- Do not place the product on the edge of a surface to avoid dropping.

- Only use the product with the original accessories.
- Never fill the housing of the product with liquid.
- Stop using the product if its parts exhibit cracks or breaks or are deformed. Only replace damaged parts with corresponding original spare parts.

4 Checking the product and package contents

NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.
1. Take the product out of the packaging.
 2. Check to make sure that all parts are included (see **Fig. A**).
 3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

5 Charging

WARNING! RISK OF INJURY!

- Charge only with the supplied USB charging cable.
- Only charge the product at temperatures between 0 and 40 °C.

- Fully charge the product before first use and when the battery power is low.
- The remaining battery power level shows in the display **3** (indicated by a number of percentage).
Recharge the product when the battery power is low.
- The product cannot be used during charging.

1. Plug the smaller end of the USB charging cable **8** into the charging port **5**.
2. Plug the larger end of the USB charging cable into the USB port of a compatible USB power source (not included) (e.g. USB wall adapter, 5 V DC, 1 A).
3. If a USB wall adapter is used, plug it into a mains socket.
Charging starts automatically. The current battery power level shows in the display (indicated by a number of percentage).
4. When charging is done, unplug the USB wall adapter from the mains socket and disconnect the USB charging cable from the charging port.

6 Operation

- Before treating an area of your feet for the first time, perform a skin sensitivity test on a spot on your foot for maximum 3 seconds. Start your treatment only if there is no adverse reaction.
- You can purchase new grinding discs on our website at www.silkn.eu.

6.1 Attaching or replacing the grinding disc

The product includes 3 grinding discs **1/6/7** (fine, medium and coarse). Start with the fine grinding disc until you are sure there is no adverse reaction.

1. Pull the cover of the grinding head **2** upwards to remove it (see step 1 in **Fig. B**).
2. Remove the attached grinding disc from the product (see step 2 in **Fig. B**).
3. Gently push the desired grinding disc onto the spindle of the grinding head until it locks into place.
4. Attach the cover of the grinding head back onto the product.

6.2 Using the product

1. Make sure your feet are clean and completely dry.
2. Press and hold the operation button **4** until the product starts operation. The product is now turned on at a low speed.
3. Press the operation button again to change the operation to a high speed if needed.

The indicator light in the display **3** turns on when the product operates at a high speed.

You can press the operation button repeatedly to change between the speed levels during operation.

4. Move the grinding head **2** over the parts of your skin you intend to remove callus. In the process, apply gentle pressure and only keep the grinding disc **1/6/7** on a spot on your skin for 2 to 3 seconds at a time, then check the result (see **Fig. C**).
 5. When done, press and hold the operation button until the product switches off.
- The product has a memory function, remembering the last selected operation speed for next use.

7 Cleaning

WARNING! RISK OF ELECTRIC SHOCK!

- Always unplug the product from the power source before cleaning.

NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners, brushes with metal bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers, and the like.

Clean the product after each use.

1. Make sure the product is turned off and is not being charged.
2. Pull the cover of the grinding head **2** upwards to remove it (see step 1 in **Fig. B**).

3. Remove the attached grinding disc **(1)**/**(6)**/**(7)** from the product (see step 2 in **Fig. B**).
4. Use a soft brush to clean the grinding disc and the inside of the grinding head.
5. Use a damp cloth to wipe off all parts of the product.
6. Wipe dry all parts of the product after cleaning.

8 Storage

All parts must be completely dry before being stored.

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product so that it is not accessible to children.
- For best protection, store the product in its original packaging.

9 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not operate.	• Make sure the product is properly charged. • Make sure the product is not connected to a USB power source. The product does not work while charging.
The product does not charge.	• Make sure the USB charging cable is securely connected between the charging port and the USB power source. • In case a USB wall adapter is used for charging, make sure it is compatible with the product (5 V DC, 1 A).

10 Technical data

Model:	FTC05
Voltage:	5 V DC, 1 A
Power:	2.5 W
Battery capacity:	1,300 mAh (lithium)

Speed settings:	2 (low / high)
Rotation speed:	4,050 rpm (low speed) 5,500 rpm (high speed)
Charging time:	2.5 hours
Use time:	Approx. 300 min (low speed) Approx. 250 min (high speed)
Weight:	164 g
Dimensions (L × W × H):	16.5 × 4.4 × 4.2 cm

II Disposal

II.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

II.2 Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

Old appliances may not be disposed of in the household waste!

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or county.** This ensures that old appliances are recycled

in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

Batteries and rechargeable batteries may not be disposed of with household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances* or not, to a collection point run by the communal authority or to a retailer, so that they can be disposed of in an

environmentally friendly manner.

Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

Note: The battery must be removed by a qualified specialist for separate disposal. Do not attempt to open or disassemble the product in order to remove the battery yourself. Note for qualified specialists: To remove the battery, switch off the product and unplug it from any charging source. Then open the housing of the product and take out the rechargeable battery.

12 Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

13 Customer service

For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: www.silkn.eu.

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Silk'n UK Ltd.