Silk'n®

LUMIGLOW

FACIO




- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DE USUARIO
- IT MANUALE D'USO
- TR KULLANIM KILAVUZU
- PT MANUAL DO UTILIZADOR

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A

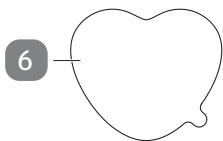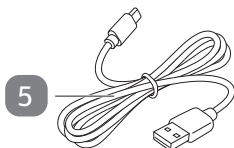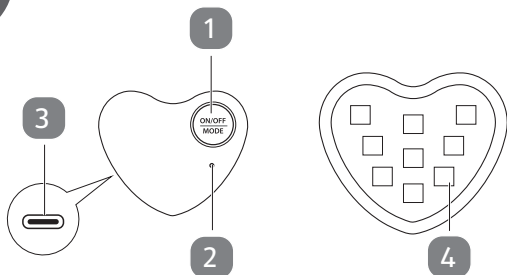Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silk'n Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silk'n Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.

Invention Works B.V.
Galileïstraat 17, 3029 AL Rotterdam, The Netherlands
www.silkn.eu – info@silkn.eu

Table of contents

1	Package contents/device parts _____	5
2	General information _____	5
2.1	Reading and storing the user manual _____	5
2.2	Explanation of symbols _____	5
3	Safety _____	6
3.1	Proper use _____	6
3.2	Safety instructions _____	6
4	Checking the product and package contents _	10
5	Charging _____	10
6	Operation _____	11
6.1	Carrying out the session _____	11
6.2	Selecting the operation modes _____	12
7	Regimen _____	12
8	Cleaning _____	12
9	Storage _____	13
10	Troubleshooting _____	13
11	Technical data _____	14
12	Disposal _____	14
12.1	Disposing of the packaging _____	14
12.2	Disposing of the product _____	14
13	Warranty _____	15
14	Customer service _____	15

I Package contents/device parts

- | | |
|--------------------|----------------------|
| 1 Operation button | 4 LED lights |
| 2 Indicator light | 5 USB charging cable |
| 3 Charging port | 6 Adhesive strip |

2 General information

2.1 Reading and storing the user manual

This user manual accompanies this LumiGlow (hereafter referred to as the “product”), and contains important information on setup and handling.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product. The user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

WARNING!

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

NOTICE!

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the United Kingdom.

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

This symbol identifies voltage and the current rating for direct current.

3 Safety

3.1 Proper use

The product is exclusively designed to help reduce blemishes and make the skin younger looking. It is only intended for private use and not suitable for commercial purposes.

Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer cannot be held liable for damages or injury incurred through improper or incorrect use.

3.2 Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

- Only use the product with the supplied USB charging cable.
- The product must only be supplied at safety extra-low voltage corresponding to the marking on the product.
- Do not use the product if it is damaged or if the USB charging cable is defective.
- Do not open the housing; instead, have a qualified professional perform repairs. Contact the customer service team for this. Liability and warranty claims are waived in the event of repairs performed by the user,

improper connection or incorrect operation.

- Only parts that comply with the original device data may be used for repairs. This product contains electrical and mechanical parts which are essential for providing protection against sources of danger.
- Do not immerse the product or USB charging cable in water or other liquids. Keep them away from water.
- Always switch the product off when you are not using it, when you intend to clean it, or in the event of a malfunction. Always unplug the USB charging cable when you are not charging the product.
- Never touch the USB charging cable with wet or damp hands.
- Never insert any objects into the housing.
- Keep all parts of the product away from open flames and hot surfaces.
- Lay the USB charging cable so that it does not pose a tripping hazard.
- Do not pull the USB charging cable out of the USB power source by the cord; instead, always pull it out by the cable plug.
- Do not kink the USB charging cable and do not lay it over sharp edges.
- Only use the product indoors.
- Do not use the product while charging.

WARNING! RISK OF INJURY!

- The product cannot be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge unless they have been given supervision or instruction concerning the use of the product in a safe way and understand the hazards involved. Children must not play with the product. Cleaning and user maintenance must not be performed by children.
- Keep children away from the product and USB charging cable.
- The product cannot be used by pregnant women, children, people with epilepsy, thyroid patients, people using thyroid medications, people who are photoallergic or use light allergy medication.
- The product cannot be used on skin diseases, wound ulceration, trauma and post-operative wounds.
- Consult your physician if you are uncertain whether you can use the product.
- Stop using the product immediately if you feel any discomfort, develop any abnormal reactions (such as itching or uncomfortable inflammation) or in case your skin becomes warm to the point of discomfort.

- Under abusive conditions, liquid may be ejected from the battery; avoid contact. Liquid ejected from the battery may cause irritation or burns.
- If battery fluid comes into contact with skin, eyes or mucous membrane, wash and rinse the affected areas immediately with clean water. Consult a physician if necessary.
- The battery must not be disassembled, short-circuited or thrown into a fire.
- Do not try to remove the battery from the product. This product contains batteries that are non-replaceable.
- Do not look directly at the light source part of the product.

NOTICE! RISK OF DAMAGE!

- Keep the product away from open fire, hot surfaces, high temperatures or the effects of the weather (e.g. rain).
- Do not stack heavy objects on top of the product.
- Stop using the product if its parts exhibit cracks or breaks or are deformed. Only replace damaged parts with corresponding original spare parts.

4 Checking the product and package contents

NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.
1. Take the product out of the packaging.
 2. Check to make sure that all parts are included (see **Fig. A**).
 3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

5 Charging

WARNING! RISK OF INJURY!

- Charge only with the supplied USB charging cable.
 - Only charge the product at temperatures between 0 and 40 °C.
-  • Fully charge the product before first use and when the battery power is low. The charging time is around 2 hours.
- When the battery power is low, the indicator light **2** flashes red and the product beeps 5 times as a warning. The current treatment will continue but the product needs to be charged before next use.
 - The product cannot be used during charging.
 - If the product is to be stored for an extended period, the charge level should be checked and the battery should be charged regularly.
1. Plug the smaller end of the USB charging cable **5** into the charging port **3**.

2. Plug the larger end of the USB charging cable into the USB port of a compatible USB power source (not included) (e.g. USB wall adapter, 5 V DC, 1 A).
3. If a USB wall adapter is used, plug it into a mains socket.
The indicator light **2** flashes green during charging. Charging is done when the indicator light illuminates green continuously.
4. When charging is done, unplug the USB wall adapter from the mains socket and disconnect the USB charging cable from the charging port.

6 Operation

WARNING! RISK OF INJURY!

- Do not use the product for more than 5 minutes per session and max. once a day. The FI mode should not be used more than 3 times per week.

6.1 Carrying out the session

1. Fully charge the product.
2. Thoroughly clean and dry your face. Make sure it is free from any make-up and dirt as they can reduce the amount of light entering the skin.
3. Attach the adhesive strip **6** on the back of the product and then stick the adhesive strip on the face to carry out the session.
4. Press and hold the operation button **1** to switch on the product.
The product beeps once and the indicator light **2** illuminates white when switched on.
The LED lights **4** turn on. By default, the "FI" operation mode is activated.
5. If desired, press the operation button repeatedly to switch between operation modes.
The product beeps once when the operation mode is changed. See chapter "Selecting the operation modes" for more information.
6. The product beeps twice and automatically switches off after 5 minutes.

You can also press and hold the operation button to switch off the product at any time.

6.2 Selecting the operation modes

- Press the operation button  repeatedly to cycle through the operation modes (F1, F2 and F3).
- The operation timer (5 minutes) is reset whenever you change to another operation mode.
- Description of the operation modes:

	Mode	Best for
F1	Blue + Red (a purple colour is visible)	For anti-ageing and the reduction of blemishes and imperfections to improve the skin's vitality, radiance and clarity.
F2	Red + Near-Infrared (a red colour is visible)	For younger-looking and tighter skin.
F3	Amber + Red + Near-Infrared (a dark orange colour is visible)	For calmer, younger-looking and tighter skin, and a more even skin tone.

7 Regimen

This product is suitable for daily use. However, do not use the product for more than 5 minutes per session and max. once a day. The F1 mode should not be used more than 3 times per week.

8 Cleaning

WARNING! RISK OF ELECTRIC SHOCK!

- Always unplug the product from the power source before cleaning.

NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners (e.g.

vinegar or limescale removers), brushes with metal or nylon bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers and the like.

- Do not immerse the product in water for cleaning. The product is not waterproof.

1. Clean the product with a damp cloth.
2. Dry the product thoroughly with a soft and clean cloth.

9 Storage

All parts must be completely dry before being stored.

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product so that it is not accessible to children.
- For best protection, store the product in its original packaging.

10 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not operate.	<ul style="list-style-type: none"> • Make sure the product is properly charged. • Make sure the product is not connected to a USB power source. The product does not work while charging.
The product does not charge.	<ul style="list-style-type: none"> • Make sure the USB charging cable is securely connected between the charging port and the USB power source. • In case a USB wall adapter is used for charging, make sure it is compatible with the product (5 V DC, 1 A).

II Technical data

Model:	FACIO
Rated voltage:	3.7 V
Battery capacity:	42 mAh lithium
Number of LED lights:	9
Materials:	ABS + Silicone
Weight:	19.5 g
Dimensions (L × W × H):	46.8 × 40.2 × 19.2 mm

12 Disposal

12.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

12.2 Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

Old appliances may not be disposed of in the household waste!

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or county.**

This ensures that old appliances are recycled in a professional manner and also rules out negative

consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

Batteries and rechargeable batteries may not be disposed of with household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances* or not, to a collection point run by the communal authority or to a retailer, so that they can be disposed of in an

environmentally friendly manner.

Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

Note: To remove the rechargeable battery, wrap a towel around the product and use a hammer to smash the housing. Make sure your eyes and hands are protected to avoid injury. Please note, the removal of the rechargeable battery is not reversible.

13 Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

14 Customer service

For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: www.silkn.eu.

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Copyright 2026 © Silkn Beauty Ltd. Toate drepturile rezervate.
Silk'n UK Ltd.