Silk'n®

Flash&Go Pro



H3220

EN INSTRUCTIONS FOR USE

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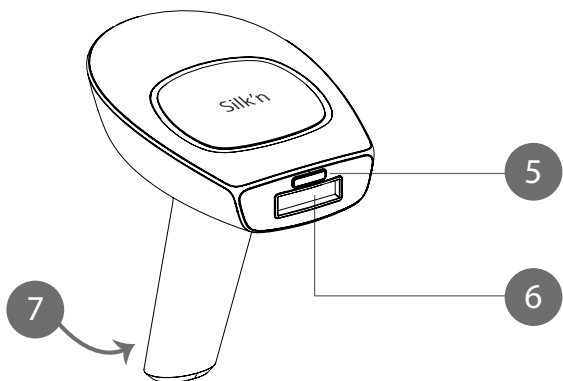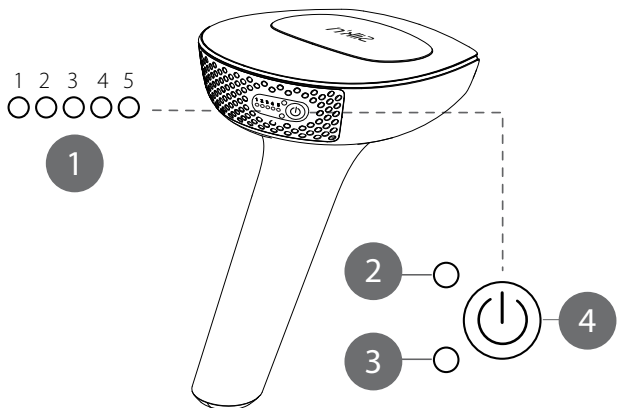⚡ Energy / Énergie / Energia

Skin color Couleur de peau Color de piel		⚡ 1-5 ✓	⚡ 1-5 ✓	⚡ 1-4 ✓	⚡ 1-3 ✓	⚡ 1-2 ✓	⚡ 1-1 ✗
Hair color Couleur de cheveux Color de pelo		⚡ 1-5 ✓	⚡ 1-4 ✓	⚡ 1-3 ✓	⚡ 1-2 ✗	⚡ 1-1 ✗	⚡ 1-0 ✗

Table of contents

1. Contraindications & warnings	2
2. Get to know your Silk'n Flash&Go Pro	6
2.1 Package contents.....	6
2.2 Device description	6
2.3 Intended use.....	7
3. Safety with Silk'n Flash&Go Pro	7
4. Safety tips.....	8
4.1 Exposure to the sun	8
4.2 Exposure of pulse light to the eyes.....	8
4.3 Wounds or infections.....	8
5. Possible side effects.....	8
6. Long-term hair removal the Silk'n way	10
6.1 Expected treatment outcomes	11
6.2 Energy levels	12
6.3 Pulsing and gliding	12
6.4 Device sensors.....	13
6.5 What to expect when treating with Silk'n	14
7. Get started.....	15
7.1 Before your first treatment	15
7.2 Treatments with Silk'n Flash&Go Pro	16
8. After treatment	17
9. Maintenance.....	18
10. Troubleshooting	18
11. Customer service	19
12. Warranty	19
13. Disposal	20
14. Specifications	20
15. Signs and symbols	22
16. Labeling	23
17. Declaration	24
18. Legal notice	28

1. Contraindications & warnings

EN

Silk'n Flash&Go Pro is not designed for everyone. Hair removal by lasers or intense pulsed light sources can cause increased hair growth in some individuals. Based upon currently available data, the highest risk groups for this response are females of Mediterranean, Middle Eastern and South Asian heritage when used on the face and neck. For more information visit www.silkn.com, www.silkn.ca or your local Silk'n domain.

- ⦿ This device must not be used by the following people: children, adults with limited physical, sensory or psychological capacities, or those lacking experience or knowledge in how to use the device. Children should be supervised to ensure that they do not play with the device.
- ⦿ Do not expose treated areas to the sun, and do not use energy levels 2-5 on tanned skin or after recent sun exposure (including artificial tanning machines). Refer to section 4.1 Exposure to the sun.
- ⦿ Use below the cheekbones and not around the eyes.
- ⦿ Do not use on areas with injected fills or toxins.

- ⦿ Do not treat on tattoos or permanent makeup, dark brown or black spots, (such as large freckles, moles, or warts), nipples, genitals, around the anus, or lips.
- ⦿ Do not treat on active eczema, psoriasis, lesions, open wounds or infections (cold sores), abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example).
- ⦿ Abnormal skin conditions caused by diabetes, systemic or other metabolic diseases.
- ⦿ Do not treat on an area that has a history of herpes outbreaks.
- ⦿ Do not use if you are pregnant or nursing.
- ⦿ Do not treat on areas where you may want hair to grow. The results are permanent.
- ⦿ Do not treat the same area of skin more than once per hair removal session.
- ⦿ Stop use immediately if your skin is very red, or if your skin blisters or burns.
- ⦿ Use a lower energy level if you experience any discomfort.
- ⦿ Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol) or acetone to clean the skin before use.

- Keep your device away from water or liquids! Do not use near bathtubs, showers, basins or other vessels containing water. Do not place or store the device where it can fall, be pushed or placed anywhere near water. There is a risk of electrocution. Unplug immediately if the device comes into contact with any fluids, and contact Customer Service.
- Remove all jewelry before use.
- Never leave the device unattended when plugged into an outlet. Keep the AC/DC adapter accessible at all times, so it can be removed from the socket easily in case of emergency.
- Do not attempt to modify, open or repair your Silk'n Flash&Go Pro device. You may expose yourself to dangerous electrical components or pulsed light energy. You may cause serious harm to yourself, and you will void your warranty.

Do not use Silk'n Flash&Go Pro if:

- The adapter, attachments or accessories for this device are not recommended by Silkn Beauty Ltd.
- The device has a damaged cord or adapter. Keep the power cord away from heated surfaces.
- You see or smell smoke when it is in use.
- The skin color sensor, treatment surface, or device casing is cracked or coming apart.

- The device is not working or appears to be damaged in any other way.

In such cases, contact Silk'n Customer Service.

Consult your physician before use if you have:

- a history of keloidal scar formation.
- used medication that makes the skin more sensitive to light, including non-steroidal anti-inflammatory agents, (e.g., aspirins, ibuprofens, acetaminophen), tetracyclines, phenothiazines, thiazide, diuretics, sulfonamides, DTIC, fluorouracil, vinblastine, griseofulvin, Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A®, Accutane® and/or topical retinoids.
- been treated, or receiving current treatment, with Alpha-Hydroxi Acids (AHAs), Beta-Hydroxi Acids (BHAs), Retin-A®, topical retinoids or azelaic acid.
- been treated with Accutane® (isotretinoin) within the past 6 months.
- epilepsy.
- a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- a history of skin cancer or areas of potential skin malignancies.
- received radiation therapy or chemotherapy treatments within the past 3 months, or have areas of potential skin malignancies.
- an active implant, such as a pacemaker, incontinence device, insulin pump, etc.
- used physician-prescribed drugs on your skin within the past 6 months.

- any other condition which, in the opinion of your physician or dermatologist, would make it unsafe for you to be treated.

2. Get to know your Silk'n Flash&Go Pro

Silk'n Flash&Go Pro is a light-based device for permanent hair reduction designed for home-use. It is a powerful electrical device and should be used with special attention to safety. Be sure to read all warnings, contraindications and safety precautions before use, and follow them strictly. We recommend you re-familiarize yourself with the instructions before each treatment.

2.1 Package contents

Your Silk'n Flash&Go Pro package includes:

- Hair removal device
- User manual
- Storage case
- Cleaning cloth
- Adapter with cable

2.2 Device description

Refer to the diagram at the front of this manual to see the location of each feature of your Silk'n Flash&Go Pro.

- Energy level indicator lights (1-5) (#1)
- Skin color sensor light (#2)
- Ready light (#3)
- Control button (#4)
- Skin color sensor (#5)
- Treatment surface (#6)
- Adapter socket (#7)

2.3 Intended use

Silk'n Flash&Go Pro is intended to remove unwanted hair and effect long-term, or permanent hair reduction. Permanent hair reduction is defined as a long-term stable reduction in the number of hairs re-growing after a treatment regime. You may use Silk'n Flash&Go Pro to remove unwanted body hair. Ideal body areas include the underarms, bikini line, arms, legs, face (below the cheekbones), back, shoulders and chest.

Do not use this device for any purpose other than as described in this manual.

3. Safety with Silk'n Flash&Go Pro

With Silk'n Flash&Go Pro, safety comes first

Home Pulsed Light (HPL) technology in Silk'n Flash&Go Pro provides superior safety with lower energy levels. HPL achieves long-term hair removal results at a fraction of the energy level used in other light-based hair removal equipment. The low energy reduces the risk of harm or complications, and contributes to your overall safety.

Silk'n Flash&Go Pro has additional features to protect your:

- **Eyes:** Light pulses are only active when the treatment surface (#6) is in full contact with the skin.
- **Skin:** Your skin is never in direct contact with the optical lens through which light pulses are emitted.

4. Safety tips

This section has information to help you use your device safely.

4.1 Exposure to the sun

- The lowest energy level (1) is designed for safe use even after sun exposure.
- If you use energy levels 2-5, avoid further exposure to the sun for 4 weeks before and 2 weeks after your treatment. This is because higher quantities of melanin expose the skin to greater risk when using any method of light-based hair removal. This applies to all skin types, even those that don't seem to tan quickly.
- We don't advise it, but if you have to go out in the sun after a treatment, be sure to apply sunscreen SPF30 or higher to the treated area. Do this for 2 weeks after a treatment.

4.2 Exposure of pulse light to the eyes

- Although you can use the device to treat facial hair, take extra care to avoid contact with the eyes, and use only below the cheekbones.

4.3 Wounds or infections

- Wait for the affected area to heal before further use.

5. Possible side effects

If you follow the instructions fully and pay attention to all warnings and contraindications, complications and side effects are uncommon. However every cosmetic procedure, including those designed for home use, involves some degree of risk.

Therefore, it is important that you understand and accept the risks and complications that may occur with pulsed light hair removal systems designed for home use.

Minor skin discomfort

Home pulsed light hair removal is generally well-tolerated. However, some users feel mild discomfort during use. This is best described as a mild stinging sensation on the treated skin areas. The stinging sensation may last during or after treatment.

This means that you should use Silk'n Flash&Go Pro at a lower energy level, and if that doesn't help, you should stop using the device.

Skin redness

Your skin may become red immediately after using the device or within 24 hours. Redness generally clears up within 24 hours.

See your doctor if redness does not go away within 2 to 3 days.

Make sure you use the energy level that is comfortable for you and suitable for your skin color. Refer to section 7.1 Before your first treatment.

Increased sensitivity of the skin

The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.

Excessive redness and swelling

In rare cases, treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days and should be treated with frequent applications of ice. Gentle cleansing is OK, but avoid exposure to the sun.

Blisters or burns

Cool the affected area and apply anti-burn cream. Wait until the area has healed completely before another treatment, and use a lower energy level. **If the skin does not return to normal within 3 days, consult your physician.**

Scarring

Although very rare, scarring may occur. Usually when scarring occurs, it is in the form of a flat and white lesion on the skin (hypotrophic). However, it can be large and red (hypertrophic) or large and extend beyond the margins of the injury itself (keloid). Subsequent aesthetic treatments may be required to improve the appearance of the scar.

Pigment changes

Silk'n Flash&Go Pro targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless, there is a risk of temporary hyper-pigmentation (increased pigment or brown discoloration) or hypo-pigmentation (whitening) to the surrounding skin. This risk of changes in skin pigmentation is higher for people with darker skin tones.

Usually discoloration or changes to skin pigment are temporary and hyperpigmentation or hypo-pigmentation rarely occur.

Infection

Infection of the skin is exceedingly rare but is still a possible risk following a skin burn or wound caused when using Silk'n Flash&Go Pro.

Bruising

Very rarely, treatment may cause blue-purple bruising that can last 5 to 10 days. As the bruising fades, there may be a rust-brown discoloration of the skin (hyper-pigmentation).

6. Long-term hair removal the Silk'n way

The efficiency of hair removal varies from person to person according to body area, hair color, and how Silk'n Flash&Go Pro is used. A typical full hair growth cycle may take 18-24 months. During this time, multiple treatment sessions are required to achieve permanent hair removal.

A typical Silk'n Flash&Go Pro hair removal plan during a full hair growth cycle includes:

- **Treatments 1-4:** Plan two weeks apart.
- **Treatments 5-7:** Plan four weeks apart.
- **Treatments 8+:** Treat as needed, until desired results are achieved.

6.1 Expected treatment outcomes

Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.

Skin Type ¹	Body Hair ²	Light Fluence [J/cm ²]	Part of the Body	Average Treatment Time	Effectiveness of Hair Removal
1-5	Dark blonde to brown	2.5-4	Lower or upper legs		60%
			Arms	Session 10-12, approx. 4 weeks apart.	50%
			Bikini line		60%
			Underarms		55%
1-5	Brown to black	2.5-4	Lower or upper legs	Sessions 8-10 approx. 4 weeks apart.	70%
			Arms		60%
			Bikini line		70%
			Underarms		65%

¹ Skin types are based on the Fitzpatrick scale of skin colors. Refer to the skin color chart in the diagram at the front of the manual.


² Not effective on white body hair.

6.2 Energy levels

The energy level determines the intensity of the light pulse delivered to your skin. As the energy level increases, so do hair removal results, as well as the risk of possible side effects and complications.

- The device has 5 levels of energy from 1 (the lowest) to 5 (the highest), represented by 5 energy level indicator lights (#1).
 - When the device is turned on, the energy level is set to level 1 automatically. Only one energy level indicator light (#1) will be on.
 - To change the energy level, press the control button (#4) repeatedly to increase the energy level from 1 to 5. Once at level 5, press the control button (#4) again to reduce the energy level from 5 back to 1.
- i** To determine the right energy level for your treatment area, refer to section 7.1 Before your first treatment.

6.3 Pulsing and gliding

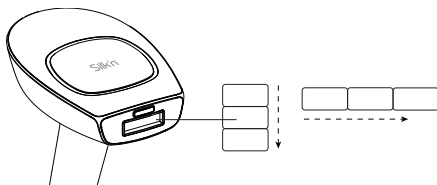Pulsing is the best technique for working at higher energy levels.

- Press the treatment surface (#6) against your skin. Work in rows (either across or down as shown in the image). Start at one end of each row and progress to the other end.

- The treatment surface (#6) creates temporary pressure marks on the treated area so you can see areas you have treated.
- After each pulse, move the treatment surface (#6) to the next spot. Keep doing this pulsing action for the area you are treating.
- Avoid gaps or overlaps as you use the device.

Gliding is the best technique for working at lower energy levels.

- Press the treatment surface (#6) firmly against your skin and slowly glide over your skin (either across or down as shown in the image). Make sure that the treatment surface (#6) is always touching your skin.
- At lower energy levels, the device takes only 1-1.5 seconds to recharge between pulses. This generates a nonstop sequence of flashes.
- Glide at a speed so that you avoid gaps or overlaps.
- ⚠ Do not treat the same area of skin more than once per hair removal session.
- ⚠ Stop use immediately if you have any discomfort.

6.4 Device sensors

Your Silk'n Flash&Go Pro has three sensors. for touch, movement, and color.

- **Touch and Movement:** The device will only emit repeated pulses when it is in firm contact with your skin, and you keep moving it across your skin.

- **Color:** Light-based hair removal on darker skin tones may result in adverse effects such as burns, blisters, and skin color changes (hyper- or hypo-pigmentation). Your Silk'n Flash&Go Pro has a skin color sensor (#5) that detects the color of your skin automatically. The skin color sensor checks your skin color at the start of, and occasionally during, each treatment.
- i If Silk'n Flash&Go Pro detects that your skin tone is dark, and the energy level is set to 3 or higher, the device will auto-adjust and set the energy level back to 2.

6.5 What to expect when treating with Silk'n

For many people, using Silk'n Flash&Go Pro may be their first experience with a light-based device. During a treatment session it is normal to experience or feel:

- **Noise:** The cooling fan makes noise similar to a hair dryer. This is normal.
 - **A popping sound:** With each flash of HPL™ light.
 - **A sensation of warmth and tingling:** During each light pulse.
 - **Some redness:** It is not uncommon to see some very mild, pink-like skin coloration during or just after a treatment. This is usually most noticeable around the hairs themselves.
- ⚠ Stop use immediately if you see your skin is very red.
- ⚠ Contact your physician immediately, if you experience any complication.

7. Get started

To get started:

1. Remove the device and other components from box.
2. Plug the adapter cord into the device socket (#7) and plug the adapter into an electrical outlet. All energy level indicator lights (#1) will turn on, and then go off.

7.1 Before your first treatment

 Remove all jewelry before use.

You must decide which energy level is right for your skin and the area you want to treat. The device has 5 levels of energy from 1 (the lowest) to 5 (the highest). For example, you might treat one area at energy level 1, and another at level 4. For information on energy levels, refer to section [6.2 Energy levels](#).

Spot Test

To choose the right energy level, test each treatment area separately (refer to the skin color chart at the front of this manual).

1. Press the control button (#4) to turn on the device. The level 1 energy level indicator light (#1) will go on.
2. Place the treatment surface (#6) on your skin to deliver one pulse. If you did not experience abnormal discomfort after 24 hours, increase the energy level on your next test spot. To increase the energy level, remove the device from your skin and press the control button (#4) again. This will increase the energy level.

3. Deliver one pulse on a different spot at energy level 2. Use the energy levels indicated on the skin color chart. Wait 24 hours. If there are no adverse effects, use a different test spot and continue testing at increased energy levels to find the highest level that is comfortable for you, without adverse effects.
4. When you find the highest energy level that is most comfortable for that treatment area, you can start your treatment.

Use a lower energy level if you experience any discomfort. Be sure to do a spot test for each part of the body you wish to treat.

7.2 Treatments with Silk'n Flash&Go Pro

Follow these steps each time you treat an area.

1. Before treating an area, the skin should be shaved, clean, dry and free of any powders, antiperspirants or deodorants.
Do not wax, pluck or tweeze the hair out.
 - Be sure to cover up moles, warts or other marks on the skin that should not be treated.
2. Press the control button (#4). The device will turn on, the level 1 energy indicator light (#1) will be on, and a fan sound will start.
3. Approximately 1 second after pressing the control button (#4), the ready light (#3) will turn on. The device is then ready to use at the lowest energy level.
4. **Set the energy level.** Press the control button (#4) repeatedly to increase the energy level until you have the level you want. The energy level indicator lights (#1) indicate the energy level setting (refer to section [6.2 Energy levels](#)).

5. Decide if you wish to treat the area by Gliding or Pulsing (refer to section [6.3 Pulsing and gliding](#)).
6. Apply the treatment surface (**#6**) to the skin. Make sure that the skin is spread evenly and smoothly. When the treatment surface (**#6**) is in full contact with the skin, the ready light (**#3**) will start blinking slowly.
7. The device will first determine the color of your skin automatically, and it will flash a pulse of light onto your skin. You will feel a mild sensation of warmth and tingling, and you may hear a popping sound.
8. If you are pulsing, remove the device from your skin. The device will immediately recharge for the next pulse. After a few seconds (depending on the energy level) the ready light (**#3**) will turn on again, and start blinking slowly when pressed against the skin.
9. If you are pulsing, repeat the process. If you are gliding, simply continue. Continue to pulse or glide until you have treated the area fully.

Stop use immediately, if your skin blisters or burns.

8. After treatment

When you have finished your treatment session:

1. Turn the device off by giving the control button (**#4**) an extended press. Remember the last energy level setting you used, because it will not be restored when you turn the device on again.
2. Remove the adapter from the electrical outlet.

9. Maintenance

It is important to clean your Silk'n Flash&Go Pro device after each treatment. To clean the device:

- Unplug the device.
 - To clean the device including the treatment surface (**#6**), wipe with a clean, dry cloth. Never use acetone or other abrasive fluids, because they will damage your device.
 - After cleaning, we recommend you store the device in its original box.
-  Never immerse your Silk'n Flash&Go Pro or any of its parts in water.

10. Troubleshooting

Do not attempt to open or repair your device. You may expose yourself to dangerous electrical components or pulsed light energy and may cause serious harm to yourself. It will also void your warranty.

Problem	Checks
My device does not start.	Make sure the adapter is properly connected to the device and that it is plugged into an electrical outlet on the wall.
There is no light pulse when the device touches my skin.	Check that the device is connected and the ready light (#3) is on. Make sure that the treatment surface (#6) is evenly and firmly pressed to the skin and that you are moving the device across your skin.

Problem	Checks
One or more lights blink.	The device has malfunctioned.
The device has stopped working and is in standby mode.	Press the control button again to start the device.

For problems or queries, please contact Silk'n Customer Service.

11. Customer service

For more information about Silk'n products, visit your regional Silk'n website: www.silkn.com or www.silkn.ca. If the device is damaged, defective, in need of repair or if you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from these websites.

Service number: 1-800-250-5053 (toll-free)

Email: customercarena@silkn.com

12. Warranty

This product is covered by a 1-year warranty. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service centre. Our customer service is always pleased to help you!

- i** The manufacturer accepts no liability for damage or injuries caused by improper or incorrect use.

13. Disposal

EN

Help protect the environment. Do not dispose with household waste and be sure to follow your country's laws for the disposal of electrical items.

14. Specifications

Model no.	H3220
Spot size	0.9cm x 3cm [2.7cm ²]
Technology	Home Pulsed Light
Maximum energy level	Max 3-5J/cm ²
Wavelength	475-1200nm
Operation & safety	Skin color sensor ensures use on appropriate skin types.
Pulse duration	500-800 micro seconds
Pulse interval	0.9 sec minimum
Repetition rate	1.1 Hz maximum
Package size	220x220x135 mm
System weight	275 gm
Transport & storage between uses and storage condition	Temperature: minus 40° - 70°C Relative humidity : 10% - 90% Atmospheric pressure: 500 - 1060 hPa
Operating conditions	Temperature: 10° - 35°C Relative humidity: 30% - 75% Atmospheric pressure: 700 - 1060 hPa

Time from minimum storage temperature between uses to minimum operating temperature with ambient temperature of 20°C degrees	30 minutes
Time from maximum storage temperature between uses to maximum operating temperature with ambient temperature of 20°C degrees	30 minutes
Mode of operation	Continuous
Protection against electrical shock	Class II Medical equipment type BF applied part
Adapter models and rating	KSAPV0361200300D5
	Input: 100-240VAC, 50-60 Hz, 0.9A. Output: 12V, 3.0A
The device was tested and complies with	IEC 60601-1 (Safety) IEC 60601-1-2 (EMC) IEC 60601-1-11 (Environment) IEC60601-2-57 (Safety of non-laser source)
Expected service life of the equipment	5 years


15. Signs and symbols

	Class II. The adapter is protected against electric shock by means of double insulation.
IP30	The device is protected against the ingress of objects larger than 2.5mm diameter, such as fingers, to protect against electric shock.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment.
	Keep dry.
	Follow the operating instructions.
	Degree of protection against electric shock, type BF applied part.
	Medical device.
	Serial number.
	Polarity.
	Date of manufacture.
	Manufacturer.

16. Labeling

Protection rating against the ingress of water: Ordinary.
This device is not suitable for use in the presence of flammable anaesthetic mixtures, or with air oxygen or nitrous oxide.

EN

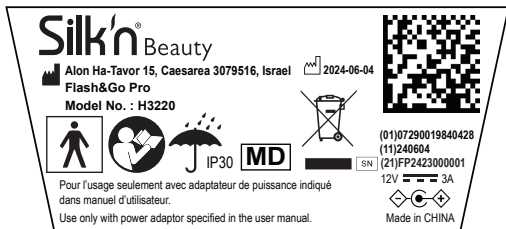The image is for illustration purposes only.

17. Declaration

EN

The use of accessories other than those specified or sold by Home Skinovations Ltd as replacement parts may have the consequence of increasing the emissions or decreasing the immunity of the device. The device is intended for use in the electromagnetic environment specified in tables 1, 2, 4 and 6.

Table 1: Guidance and Manufacturer's Declaration: Electromagnetic Emissions

The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.

Emissions Test	Compliance	Electromagnetic Environment Guidance
RF emissions CISPR 11	Group 1	The device uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment.
RF emissions CISPR 11	Class B	Complies
Harmonic emissions IEC 61000-3-2	Not Applicable	The device is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Complies	

Table 2: Guidance and Manufacturer's Declaration: Electromagnetic Immunity

The device is intended for use in the electromagnetic environment specified below. The customer or the user of the devices should assure that it is used in such an environment.

Immunity Test	IEC 60601 Test Level	Compliance Level	Electromagnetic Environment Guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 8 kV contact	± 8 kV contact	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at
	± 15 kV air	± 15 kV air	
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines	± 2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
	± 1 kV for input/output lines	NA	
Surge IEC 61000-4-5	2 kV 0, 90, 180 & 270 degree phase angles	2 kV 0, 90, 180 & 270 degree phase angles	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	100% drop, 0.5 periods (0°, 45°, 90°, 135°, 180°, 225°, 270°, 315°)	100% drop, 0.5 periods (0°, 45°, 90°, 135°, 180°, 225°, 270°, 315°)	Mains power quality should be that of a typical commercial or hospital environment. If the user of the device requires continued operation during power mains interruptions, it is recommended that the device be powered from an uninterruptible power supply or a battery.
	100% dip, 1 period	100% dip, 1 period	
	30% dip, 25/30 periods	30% dip, 25/30 periods	
	100% drop, 5 seconds	100% drop, 5 seconds	
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.
Note	<i>UT is the a.c. mains voltage prior to application of the test level.</i>		

Table 3: Guidance and Manufacturer's Declaration: Electromagnetic Immunity

The device is intended for use in the electromagnetic environment specified below. The customer or the user of the device should assure that it is used in such an environment.

Immunity Test	IEC 60601-1-2 Test Level	Compliance Level	Electromagnetic Environment Guidance
Conducted RF IEC 61000-4-6	6 V/rms 150 kHz to 80 MHz	6 V/rms 150 kHz to 80 MHz	Portable and mobile RF communications equipment should be used no closer to any part of the device, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter. Recommended separation distance. $d = \left[\frac{3.5}{f_1} \right] \sqrt{P}$
Radiated RF	10 V/m	10 V/m	$d = \left[\frac{3.5}{f_1} \right] \sqrt{P}$ 80 MHz to 800 MHz
IEC 61000-4-3	80 MHz to 2.7 GHz	80 MHz to 2.7 GHz	$d = \left[\frac{7}{f_1} \right] \sqrt{P}$ 800 MHz to 2.5 GHz Where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m). Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey, should be less than the compliance level in each frequency range. Interference may occur in the vicinity of equipment marked with the following symbol: 

Note 1: At 80 MHz and 800 MHz, the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

A Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the device is used exceeds the applicable RF compliance level above, the device should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the device.

B Over the frequency range 150 kHz to 80 MHz, field strengths should be less than [V1] V/m.

Table 4: Recommended Separation Distances between Portable and Mobile RF Communications Equipment and the device

The device is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the device can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the device as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter		
	m		
	150 kHz to 80 MHz $d = [\frac{3,5}{V_1}] \sqrt{P}$	80 MHz to 800 MHz $d = [\frac{3,5}{E_1}] \sqrt{P}$	800 MHz to 2,5 GHz $d = [\frac{7}{E_1}] \sqrt{P}$
0.01	0.12	0.12	0.23
0.1	0.37	0.37	0.74
1	1.17	1.17	2.33
10	3.7	3.7	7.38
100	11.7	11.7	23.33

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

Note 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

Note 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

18. Legal notice

EN

Copyright 2025 © Silkn Beauty Ltd. All rights reserved.

Silkn Beauty Ltd. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Silkn Beauty Ltd. is believed to be accurate and reliable at the time of publication. However, Silkn Beauty Ltd. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Silkn Beauty Ltd.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Silkn Beauty Ltd. Data is subject to change without notification. Silkn Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silkn Beauty Ltd. Specifications are subject to change without notice.

SILK'N and the SILK'N logo are registered trademarks of Silkn Beauty Ltd. Galileïstraat 17, 3029AL Rotterdam, The Netherlands, www.silkn.com, www.silk.ca, info@silkn.com.

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Copyright 2025 © Silkn Beauty Ltd. Tous droits réservés.
Copyright 2025 © Silkn Beauty Ltd. Todos los derechos reservados.
Silk'n Beauty Ltd.