Silk'n®

Motion Premium

H3220

- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DEL USUARIO
- PT MANUAL DO UTILIZADOR
- GR ΕΓΧΕΙΡΙΔΙΟ ΧΡΗΣΤΗ
- TR KULLANIM KILAVUZU
- IT MANUALE PER L'UTENTE
- NO BRUKSANVISNING
- SV ANVÄNDARHANDBOK
- DA BRUGERMANUAL

Motion Premium


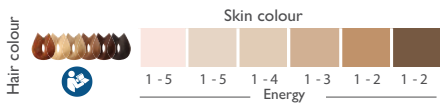
Table of contents

1. Warnings	2
2. Contraindications	4
3. Safety features	6
3.1 Skin protection	6
3.2 Eye protection	7
3.3 Exposure to the sun	7
4. Get to know your device	8
4.1 Device description	8
4.2 Intended use	8
4.3 Package contents	8
4.4 Device parts	9
4.5 Device controls	9
4.6 Energy levels	10
4.8 Pulsing and gliding	11
5. Treatments	12
5.1 Treatment plan	12
5.2 What to expect with the device	13
5.3 Possible side effects	14
6. Get started!	16
6.1 Connecting the device	16
6.2 Pre-treatment preparation	17
6.3 Perform a spot test	17
6.4 Treatment steps	18
6.5 After treatment	19
7. Cleaning, maintenance & storage	20
8. Troubleshooting	20
9. Customer service	21
10. Warranty	22
11. Disposal	22
12. Specifications	23
13. Signs & symbols	24
14. Labeling	25
15. Technical results	26
16. Legal notice	27

1. Warnings

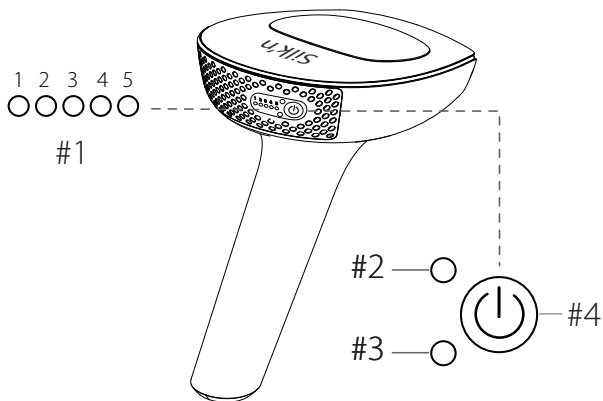This user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws.

Before first use, please read this user manual in full, and keep it for future use. We recommend you refamiliarise yourself with the instructions before each treatment.

- Always keep the device or any of its parts dry. Ensure that your hands are dry before connecting the adapter to the electricity supply.
- Do not use the device in the bathroom.
- Never use the adapter if it has a damaged cable or plug, or isn't working properly, or has been submerged in water. If that happens, unplug the adapter immediately before touching it. Failure to unplug the adapter may result in an electric shock!

- ⦿ Do not operate the device in environments where there is a risk of explosion. The device is not suitable for use in the presence of flammable anaesthetic mixture with air, oxygen or nitrous oxide. Keep the device and all its parts away from open flames and hot surfaces.
- ⦿ Keep the adapter accessible at all times, so it can be removed from the electricity supply easily in case of emergency.
- ⦿ Always switch the device off and disconnect it from the electricity supply when you are not using the device, when you intend to clean it, or in the event of a fault.
- ⦿ Never pull the adapter out of the socket by the cable. Always remove it by taking hold of the adapter.
- ⦿ Do not kink the adapter cable or lay it over sharp edges. Place it in such a way that it does not pose a tripping hazard.
- ⦿ Possible hazardous optical radiation is emitted from this device. Do not look at the operating lamp. Eye injury may result.
- ⦿ Do not use the device where there is monitoring equipment (e.g., cardiac monitors, ECG alarms), which may not operate properly when the device is in use.

- ⦿ Do not use if the parts or accessories are not supplied by the manufacturer as defined in this manual, or if the device is damaged in any way, does not appear to be working properly, or if you see or smell smoke.
- ⦿ Do not attempt to modify, open or repair the device. You may expose yourself to an electric shock. It will also void the warranty.
- ⦿ Children must not play with the device or any of its parts. They may injure themselves or swallow small parts and choke.

2. Contraindications

- ⦿ The device is suitable for adults aged 16 years and above.
- ⦿ The device must not be used by the following people: children, adults with limited physical, sensory or psychological capacities, or those lacking experience or knowledge in how to use the device.
- ⦿ The device is not designed for everyone. Hair removal by intense pulsed light sources can cause increased hair growth in some individuals, although it is extremely rare (less than 1% of the population). Potential areas for women are the face and neck, and for men, the back, shoulders and upper arms, although it could also affect other areas. The highest risk groups are people from the Mediterranean, Middle East and South Asia.
- ⦿ The device is not suitable for use on children or animals.

 **Do not use the device if you:**

- have naturally dark skin (refer to the skin colour chart at the front of this manual). Treating dark skin can result in adverse effects such as burns, blisters, and skin colour changes (hyper- or hypopigmentation).
- have naturally white, grey, light blond, or red body hair.
- are pregnant or nursing.
- have a history of keloidal scar formation.
- have suffered or suffer from epilepsy.
- have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- have a history of skin cancer, a current condition of any cancer, or areas of potential skin malignancies.
- have received drugs from your physician for a skin condition within the past 6 months.
- are allergic to metal, or are sensitive to contact with chrome.
- have any other condition which in your physician's opinion would make it unsafe for you to use the device.
- cannot sense heat properly because of an illness, or other health related issue.

 **Do not use on these parts of the body:**

- on tanned skin.
- around or near the eyes. Use below the cheekbones.
- on tattoos or permanent makeup, dark brown or black spots, (such as large freckles, moles, or warts), nipples, genitals, around the anus, or lips.
- areas with injected fillers or toxins.

- a part of the body that has a history of herpes outbreaks.
- skin that has active eczema, psoriasis, lesions, open wounds or infections (cold sores), abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example).
- areas that are prone to rosacea.
- on varicose veins.
- areas that are painful.
- areas where you may want hair to grow later. The results are permanent.

3. Safety features

With Silk'n safety comes first

Home Pulsed Light™ (HPL) technology provides superior safety with lower energy levels. HPL achieves long-term hair removal reduction results at a fraction of the energy level used in other light-based hair removal equipment. The low energy reduces the risk of harm or complications, and contributes to your overall safety.

3.1 Skin protection

Light-based hair removal is not suitable for tanned skin. The device comes with a built-in skin colour sensor (**#5**) that is designed to measure the complexion of the applied surface and enables energy levels according to skin complexion.

Furthermore, the treatment surface (**#6**), through which pulses of light are delivered, is recessed inside the device. This enables The device to protect your skin by avoiding direct contact between the treatment surface (**#6**) and the skin.

3.2 Eye protection

The device has a built-in safety feature for eye protection. It has been designed so that a light pulse cannot be emitted when the treatment surface (#6) is facing open air. The safety feature enables treatment only when the treatment surface (#6) is in contact with the skin.

3.3 Exposure to the sun

- Sun exposure includes constant unprotected exposure to direct sunlight for over 15 minutes, or constant unprotected exposure to diffuse sunlight for over 1 hour.
- Do not use on tanned skin after recent sun exposure, including tanning machines. Using the device on skin that has been recently exposed to sunlight can result in adverse effects such as burns, blisters, and skin colour changes (hyper- or hypopigmentation).
- Avoid further exposure to the sun for **4 weeks before and 2 weeks after** your treatment. This is because higher quantities of melanin expose the skin to greater risk when using any method of light-based hair removal. This applies to all skin types, even those that don't seem to tan quickly.

4. Get to know your device

4.1 Device description

The device is a light-based device for permanent hair reduction designed for home use. It is a powerful electrical device and should be used with special attention to safety.

4.2 Intended use

The IPL device is an over-the-counter device intended to remove unwanted hair and affect long term, or permanent hair reduction. Permanent hair reduction is defined as a long-term stable reduction in the number of hairs regrowing after a treatment regime. You may use the device to remove unwanted body hair. Ideal body areas include the underarms, bikini line, arms, legs, face (below the cheekbones), back, shoulders and chest.

The device is a medical appliance and is intended exclusively for cosmetic treatments. It is for private use, and is not suitable for commercial or therapeutic use. Do not use for any purpose other than as described in this manual.

4.3 Package contents

- Hair removal device
- Adapter with cable
- Cleaning cloth
- Storage case
- User manual

4.4 Device parts

Refer to the diagram on page (i) at the front of this manual to see the location of each feature.

- Energy level indicator lights (1-5) **(#1)**
- Skin colour sensor light **(#2)**
- Ready light **(#3)**
- Control button **(#4)**
- Skin colour sensor **(#5)**
- Treatment surface **(#6)**
- Adapter socket **(#7)**

4.5 Device controls

Button	Function / Mode
Control button (#4) 	<ul style="list-style-type: none"> • Long press: turns the device on or off. • Short press: each press increases the energy level. When at energy level 5, the next press moves back to energy level 1.
Energy level indicator lights x5 (#1) 1 2 3 4 5 ○ ○ ○ ○ ○	<ul style="list-style-type: none"> • Green constant: the energy level is active. • The energy level 1 indicator light is green constant when the device is turned on. • The number of energy level indicator lights on corresponds to the energy level selected.
Ready light (#3) 	<ul style="list-style-type: none"> • Green constant: when the device is turned on. • Light off: the device is getting ready for the next pulse. • Green blinks: cannot detect the skin colour. or there is no movement across the skin.
Skin colour sensor light (#2) 	<ul style="list-style-type: none"> • Green blinks: skin colour is too dark for safe use.

4.6 Energy levels


- The energy level determines the intensity of the light pulse delivered to your skin. As the energy level increases, so do the results, as well as the risk of possible side effects and complications (see section [5.3 Possible side effects](#)).
- The device has 5 levels of energy from 1 (the lowest) to 5 (the highest), represented by 5 energy level indicator lights **(#1)**. To change levels, short press the control button **(#4)**.
- To determine the right energy level for your treatment area, refer to section [6.3 Perform a spot test](#).

4.7 Device sensors

The device has sensors for touch and colour.

- **Touch:** The device will only emit pulses when it is in firm contact with your skin.
 - **Colour:** Light-based hair removal on darker skin tones may result in adverse effects such as burns, blisters, and skin colour changes (hyper- or hypopigmentation). The device has a skin colour sensor **(#5)** that detects the colour of your skin automatically. The skin colour sensor **(#5)** checks your skin colour at the start of, and occasionally during, each treatment.
- i** If the device detects that your skin tone is dark, and the energy level is set to 3 or higher, the device will auto-adjust and set the energy level back to 2.

4.8 Pulsing and gliding

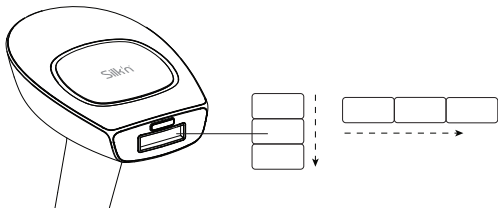Whether you decide to pulse or glide:

- Avoid gaps or overlaps as you use the device. Overlaps increase the likelihood of adverse affects.
- Do not treat the same area of skin more than once per treatment session.

Pulsing is the best technique for treating specific spots or small areas.

- Press the treatment surface **(#6)** against your skin. Work in rows (either across or down as shown in the image). Start at one end of each row and progress to the other end.
- The treatment surface **(#6)** creates temporary pressure marks on the treated area so you can see areas you have treated.
- After each pulse, move the treatment surface **(#6)** to the next spot. Keep doing this pulsing action for the area you are treating.

Gliding is the best technique for treating larger parts of the body.

- Press the treatment surface **(#6)** firmly against your skin and slowly glide over your skin (either across or down as shown in the image). Make sure that the treatment surface **(#6)** is always touching your skin.
- At lower energy levels, the device takes only 1-1.5 seconds to recharge between pulses. This generates a nonstop sequence of flashes.

5. Treatments

5.1 Treatment plan

The efficiency of hair removal varies from person to person according to body area, hair colour, and how the device is used. A typical full hair growth cycle may take 18-24 months. During this time, multiple treatment sessions are required to achieve permanent hair reduction.

A typical hair removal plan during a full hair growth cycle includes:

- **Treatments 1-4:** Plan two weeks apart.
- **Treatments 5-7:** Plan four weeks apart.
- **Treatments 8+:** Treat as needed, until desired results are achieved.

Ideal body areas include underarms, bikini area, arms, legs, back, shoulders and chest (excluding nipples). For the face, use below the cheekbones but not around or near the eyes or lips.

Skin Type ¹	Body Hair ²	Light Fluence [J/cm ²]	Part of the Body	Average Treatment Time	Effectiveness of Hair Reduction
1-6	Dark blond to brown	2.5-4	Lower or upper legs	Sessions 10-12 approx. 4 weeks apart	60%
			Arms		50%
			Bikini line		60%
			Underarms		55%
1-6	Brown to black	2.5-4	Lower or upper legs	Sessions 8-10 approx. 4 weeks apart	70%
			Arms		60%
			Bikini line		70%
			Underarms		65%

1 Individual response depends on hair type as well as biological factors that may affect hair growth patterns. Some users may respond quicker or slower than the average number of treatments.

2 One cannot expect permanent hair reduction in a single or even double treatment sessions. The duration of the resting period for hair follicles also depends on the body area.

5.2 What to expect with the device

During a treatment session it is normal to experience and feel:

- **A fan noise:** The cooling fan makes a noise similar to that of a hair dryer.
- **A pop sound with each pulse:** When a pulse of HPL light is activated, it is normal to hear a subtle pop sound simultaneously with the flash of light.

- **A sensation of warmth and tingling:** During each light pulse.
 - **Some mild red or pink colour:** During and just after your treatment session it is not uncommon to see some very mild, pink-like colour of the skin. This is usually most noticeable around the hairs themselves.
- ⚠ Stop use immediately if you see full redness of the skin, blistering or burns.

5.3 Possible side effects

When used according to the instructions, side effects and complications associated with the use of the device are uncommon. However, every cosmetic procedure, including those designed for home use, involves some degree of risk. Therefore, it is important that you understand and accept the risks and complications that can occur with pulsed light hair removal systems designed for home use.

There is a minor risk that stacking or overlapping of multiple pulses on the same skin spot may have some adverse effects.

Minor skin discomfort

Home pulsed light hair removal is generally well-tolerated. However, some users feel mild discomfort during use. This is best described as a mild stinging sensation on the treated skin areas. The stinging sensation may last during or after treatment.

If you feel pain and a notable skin irritation during treatment, reduce the energy level. If that doesn't help, stop the treatment.

Skin redness

Your skin may become red immediately after using the device or within 24 hours. Redness generally clears up within 24 hours.

See your doctor if redness does not go away within 2 to 3 days.

Make sure you use the energy level that is comfortable for you and suitable for your skin colour.

Excessive redness and swelling

In rare cases, treated skin may become very red and swollen. This is more common in sensitive areas of the body. The redness and swelling should subside within 2 to 7 days and should be treated with frequent applications of ice. Gentle cleansing is OK, but avoid exposure to the sun.

Increased sensitivity of the skin

The skin of the treated area is more sensitive so you may encounter dryness or flaking of the skin.

Blisters or burns; skin wounds

It is very unusual but blister, burns or wounds to the skin may occur following a treatment. These can require a few weeks to heal, and in very rare cases, may leave a noticeable scar.

Cool the affected area and apply anti-burn cream. Wait until the area has healed completely before another treatment, and use a lower energy level. **If the skin does not return to normal within 3 days, consult your physician.**

Infection or bruising

Infection of the skin is exceedingly rare but is still a possible risk following a skin burn or wound caused by the device.

Very rarely, the device may cause blue-purple bruising that can last 5 to 10 days. As the bruising fades, there may be a rust-brown discolouration of the skin (hyperpigmentation).

Scarring

Although very rare, scarring may occur. Usually, the scar is in the form of a flat and white lesion on the skin (hypotrophic). However, it can be large and red (hypertrophic) or large and extend beyond the margins of the injury itself (keloid). Subsequent aesthetic treatments may be required to improve the appearance of the scar.

Pigment changes

The device targets the hair shaft, in particular the pigmented cells in the hair follicle and the hair follicle itself. Nevertheless there is a risk of temporary hyper-pigmentation (increased pigment or brown discolouration) or hypopigmentation (whitening) to the surrounding skin. Risk of changes in skin pigmentation is higher for people with darker skin tones or heavily tanned skin. Usually discolouration or changes to skin pigment are temporary, and permanent hyper-pigmentation or hypopigmentation are rare.

 Stop use immediately if any of these side effects are more severe.

 See your physician if the side effects do not go away within 2 to 3 days.

6. Get started!

6.1 Connecting the device

1. Remove the device and other components from box.
2. Plug the adapter into the device socket **(#7)** and plug the adapter into the electricity supply. All energy level indicator lights **(#1)** will turn on, and then go off.

6.2 Pre-treatment preparation

 Remove all jewellery on or near the treatment area before use.

1. Before treating an area, the skin should be shaved, clean, dry and free of any powders, antiperspirants or deodorants.
Do not wax, pluck or tweeze the hair out.

 Never use flammable liquids such as alcohol (including perfumes, sanitizers, or other applications containing alcohol), or acetone to clean the skin before use.

 Do not use the device if the treatment surface **(#6)** is broken or missing.

2. Be sure to cover up moles, warts or other marks on the skin that should not be treated. Use a material that does not absorb the light such as a white cloth or white medical tape.

6.3 Perform a spot test

You must decide which energy level is right for your skin, and the area you want to treat. For example, you might treat one area at energy level 1, and another at level 4.

To choose the right energy level, test each treatment area separately (refer to the skin colour chart at the front of this manual).

Repeat this spot test for each part of the body you wish to treat.

1. Press the control button **(#4)** to turn on the device. The level 1 energy level indicator light **(#1)** will go on.

2. Place the treatment surface **(#6)** on your skin to deliver one pulse. If you do not experience abnormal discomfort press the control button **(#4)** again. This will increase the energy level. Deliver one pulse on a different spot at energy level 2.

3. Continue testing at increased energy levels. Use the energy levels indicated on the skin colour chart. When you find the highest energy level that is comfortable for that treatment area, you can start your treatment.
- i** Use a lower energy level if you experience discomfort.

You are now ready to perform a treatment.

6.4 Treatment steps

1. Press the control button **(#4)**. The device will turn on, the level 1 energy indicator light **(#1)** will be on, and a fan sound will start.
2. After pressing the control button **(#4)**, the ready light **(#3)** will turn on. The device is then ready to use at the lowest energy level. Always start at energy level 1 for your first treatment.
3. **Set the energy level.** Press the control button **(#4)** repeatedly to increase the energy level until you have the level you want. The energy level indicator lights **(#1)** indicate the energy level setting (see section [4.6 Energy levels](#)). The number of energy level indicator lights **(#1)** that are on, indicates the energy level selected.
4. Decide if you wish to treat the area by **Pulsing** or **Gliding** (refer to section [4.8 Pulsing and gliding](#)).
5. Apply the treatment surface **(#6)** to the skin. Make sure that the skin is spread evenly and smoothly. When the treatment surface **(#6)** is in full contact with the skin, the ready light **(#3)** will start blinking slowly.
6. The device will first determine the colour of your skin automatically and it will flash a pulse of light onto your skin. You will feel a mild sensation of warmth and tingling, and you may hear a popping sound.

7. If you are pulsing, remove the device from your skin. The device will immediately recharge for the next pulse. After a few seconds (depending on the energy level) the ready light (#3) will turn on again, and start blinking slowly when pressed against the skin.
8. If you are pulsing, repeat the process. If you are gliding, simply continue. Continue to pulse or glide until you have treated the area fully.

i Avoid gaps or overlaps.

i The device will only emit repeated pulses when it is in firm contact with your skin, and you keep moving it across your skin.

 Stop use immediately, if your skin blisters or burns, if you feel very uncomfortable, or if the device gets too hot.

6.5 After treatment

When you have finished your treatment session:

1. Long press the control button (#4) to switch off the device. Remember the last energy level setting you used, because it will not be restored when you turn the device on again.
2. Disconnect the adapter from the electricity supply, then clean the device and store safely (see section 7. Cleaning, maintenance & storage).
3. After your treatment session we recommend you apply a body lotion to the treated area to help your skin recover from the treatment faster.

7. Cleaning, maintenance & storage

It is important to clean the device after each treatment.

- Disconnect the device from the electricity supply.
- To clean the device including the treatment surface (**#6**), wipe with a clean, dry cloth and 70% alcohol. Never use acetone or other abrasive fluids, because they will damage your device.
- Use a wipe to dry the device.
- After cleaning, we recommend you store the device in its original box.

8. Troubleshooting

Do not attempt to open or repair your Silk'n device. You may expose yourself to dangerous electrical components or pulsed light energy and cause serious harm to yourself. It will also void your warranty. For unresolved problems or queries, please contact customer service.

Problem	Checks
My device does not start.	Make sure the adapter is properly connected to the device and that it is plugged into an electrical outlet on the wall.
There is no light pulse when the device touches my skin.	Check that the device is connected and the ready light (#3) is on. Make sure that the treatment surface (#6) is evenly and firmly pressed to the skin. It is also possible that your skin colour is too dark.

Problem	Checks
The skin colour sensor light (#2) and the ready light (#3) blink.	The skin colour is too dark.
One or more lights blink continuously.	The device has malfunctioned.
The device has stopped working and is in standby mode.	Press the control button (#4) again to start the device.

9. Customer service

See our website for more information: www.silkn.eu. This manual is also available as a PDF download from: www.silkn.eu. Contact Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

Country	Service number	Email
Belgium	0800-29316	serviceconsommateurbe@silkn.eu
France	0891-655557	serviceconsommateurfr@silkn.eu
Germany	089 51 23 44 23	kundenservicede@silkn.eu
Italy	+31 (0)180-330 550	servizioclientiit@silkn.eu
Netherlands	0900-2502217	klantenservicenl@silkn.eu
Spain	900 823 302	servicioalcliente@silkn.eu
UK	0906-2130009	customer-careuk@silkn.eu
Other countries	+31 (0)180-330 550	info@silkn.eu

10. Warranty

EN

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service centre. Our customer service is always pleased to help you!

The manufacturer accepts no liability for damage or injuries caused by improper or incorrect use.

11. Disposal

If your device reaches its End of Life stage, be sure to follow your country's laws for the disposal of electrical items.

Packaging: Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.

Device: Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. This ensures that the device is recycled professionally, and prevents the release of harmful substances* into the environment. Make sure that the device is in an uncharged state before disposal.

* Products labelled with Cd = cadmium, Hg = mercury, Pb = lead.

12. Specifications

Model no.	H3220
Spot size	0.9cm x 3cm [2.7cm ²]
Technology	Home Pulsed Light™
Maximum energy level	Max 4J/cm ²
Wavelength	475-1200nm
Operation & safety	Skin colour sensor ensures use on appropriate skin types.
System weight	275g
Transport & storage between uses and storage condition	Temperature: -40 to +70°C Relative humidity: 10 to 90% Atmospheric pressure: 500 to 1060hPa
Operating conditions	Temperature: 10 to 35°C Relative humidity: 30 to 75% Atmospheric pressure: 700 to 1060hPa
Time from minimum storage temperature between uses to minimum operating temperature with ambient temperature of 20°C	30 minutes
Time from maximum storage temperature between uses to maximum operating temperature with ambient temperature of 20°C	30 minutes
Mode of operation	Continuous
Adapter models and rating	KSAPV0361200300D5
Input Output	100-240Vac, 50-60Hz, 0.9A 12Vdc, 3.0A


Pulse light duration Repetition rate	500-800 μ s 0.7Hz maximum
Protection against electric shock	Class II medical equipment type BF applied part
Expected service life of the device	5 years
Pulse parameters (seconds)	Energy level 1: 0.9 Energy level 2: 1.1 Energy level 3: 1.4 Energy level 4: 1.7 Energy level 5: 2.0

13. Signs & symbols

	Important information about using or maintaining the device.
	A statement or event that indicates a possible or impending danger, problem, or other unpleasant situation when using the device.
	Products marked with this symbol fulfil the requirements of the EU directives.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.
	Follow the instructions for use.
	Degree of protection against electric shock: type BF applied part.

	Keep dry.
	Date of manufacture.
	Manufacturer.
	EU Authorized Representative.
	Product serial number.
IP30	The device is protected against the ingress of objects larger than 2.5 mm diameter, such as fingers, to protect against electric shock.
	Class II. The adapter is protected against electric shock by means of double insulation.

14. Labeling

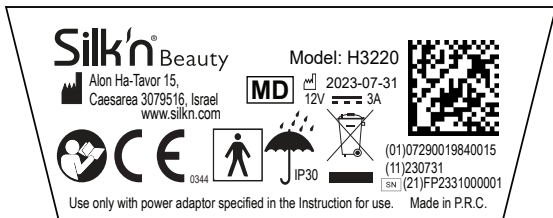15. Technical results

EMC requirements

IEC 60601-1-2 Edition 4.0 (2014)

Environment of intended uses:

Home Healthcare Facility Environment

Summary of Test Results

Test	Standard	Class/ Severity level	Test result
Emission (IEC 60601-1-2 section 7)			
Radiated emission Freq. range: 30 - 1000 MHz	CISPR 11	Group 1 Class B	Complies
Immunity (IEC 60601-1-2 section 8)			
Immunity from Electrostatic discharge (ESD)	IEC 61000-4-2	8 kV contact discharges & 15 kV air discharges	Complies
Immunity from radiated electromagnetic fields	IEC 61000-4-3	3.0; 10.0 V/m; 80 MHz + 2.7 GHz, 80% AM, 1 kHz	Complies
Immunity from Proximity field from wireless communications equipment	IEC 61000-4-3	List of frequencies, from 9 V/m up to 28 V/m, PM (18 Hz or 217 Hz), FM 1 kHz	Complies
Immunity from power frequency magnetic field	IEC 61000-4-8	30 A/m @ 50 Hz	Complies

16. Legal notice

Copyright 2024 © Silkn Beauty Ltd. All rights reserved.

Silkn Beauty Ltd reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Silkn Beauty Ltd is believed to be accurate and reliable at the time of publication. However, Silkn Beauty Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Silkn Beauty Ltd. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Silkn Beauty Ltd. Data is subject to change without notification. Silkn Beauty Ltd has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silkn Beauty Ltd. Specifications are subject to change without notice. Silk'n and the Silk'n logo are registered trademarks of Silkn Beauty Ltd, Alon Ha-Tavor 15, Caesarea, 3079516, Israel, www.silkn.com, info@silkn.com.

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Copyright 2024 © Silkn Beauty Ltd. Alle Rechte vorbehalten.
Copyright 2024 © Silk'n Beauty Ltd. Tous droits réservés.
Copyright 2024 © Silkn Beauty Ltd. Alle rechten voorbehouden.
Copyright 2024 © Silkn Beauty Ltd. Todos los derechos reservados.
Copyright 2024 © Silkn Beauty Ltd. Todos os direitos reservados.
Copyright 2024 © Silk'n Beauty Ltd. Tutti i diritti riservati.