# Silk'n®

EN





FR

ES

## VACUPEDI NV8615D

- EN USER MANUAL
- FR MODE D'EMPLOI
- ES MANUAL DE USUARIO
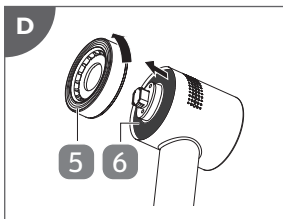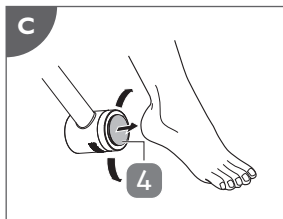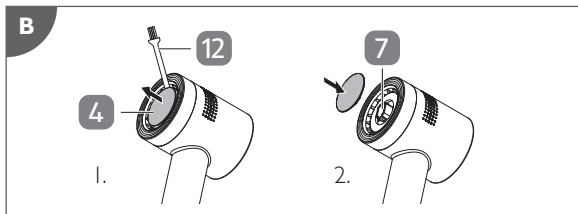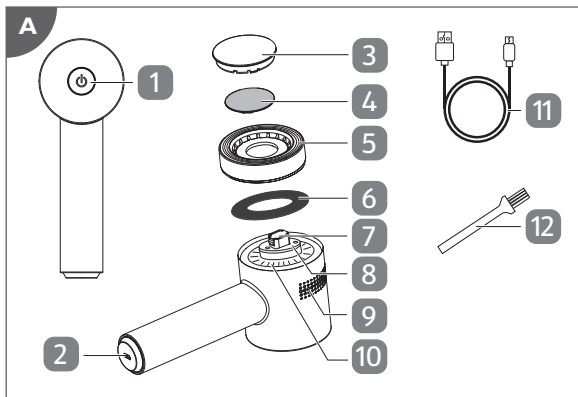
## Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silkn Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silkn Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silkn Beauty Ltd.

Invention Works B.V.  
Galileïstraat 17, 3029 AL Rotterdam, The Netherlands  
[www.silkn.com](http://www.silkn.com), [www.silkn.ca](http://www.silkn.ca) – [info@silkn.com](mailto:info@silkn.com)

# Table of contents

<b>1</b>	<b>Safety instructions</b>	<b>5</b>
<b>2</b>	<b>OPERATING INSTRUCTIONS</b>	<b>7</b>
2.1	Package contents/device parts	7
2.2	Checking the product and package contents	7
2.3	General information	8
2.4	Reading and storing the user manual	8
2.5	Explanation of symbols	8
2.6	Proper use	9
2.7	Operation	10
2.7.1	Attaching or replacing the grinding disc	10
2.7.2	Starting and ending the treatment	10
2.8	Charging	11
<b>3</b>	<b>Cleaning</b>	<b>12</b>
<b>4</b>	<b>Storage</b>	<b>13</b>
<b>5</b>	<b>Troubleshooting</b>	<b>13</b>
<b>6</b>	<b>Technical data</b>	<b>14</b>
<b>7</b>	<b>Disposal</b>	<b>14</b>
7.1	Disposing of the packaging	14
7.2	Disposing of the product	15
<b>8</b>	<b>Warranty</b>	<b>16</b>
<b>9</b>	<b>Customer service</b>	<b>16</b>

## I Safety instructions

# IMPORTANT SAFEGUARDS

When using electrical products, especially when children are present, basic safety precautions should always be followed, including the following:

### **READ ALL INSTRUCTIONS BEFORE USING** **DANGER – To reduce the risk of electrocution:**

1. Always unplug this product immediately after use.
2. Do not use while bathing.
3. Do not place or store product where it can fall or be pulled into a tub or sink.
4. Do not place in or drop into water or other liquid.
5. Do not reach for a product that has fallen into water. Unplug immediately.

### **WARNING – To reduce the risk of burns, electrocution, fire, or injury to persons.**

1. A product should never be left unattended when plugged in.
2. Close supervision is necessary when this product is used by, on, or near children or invalids.
3. Use this product only for its intended use as described in this manual. Do not use attachments not recommended by the manufacturer.
4. Never operate this product if it has a damaged cord or plug, if it is not working properly, if it has been dropped or damaged, or dropped into water. Return the product to a service center for examination and repair.
5. Keep the cord away from heated surfaces.

6. Never block the air openings of the product or place it on a soft surface, such as a bed or couch, where the air openings may be blocked. Keep the air openings free of lint, hair, and the like.
7. Never use while sleeping or drowsy.
8. Never drop or insert any object into any opening or hose.
9. Do not use outdoors or operate where aerosol (spray) products are being used or where oxygen is being administered.
10. Only use the product with the supplied USB cable.
11. Never touch the USB cable with wet or damp hands.
12. Always unplug the USB cable when you are not charging the product.
13. Keep the product and USB cable away from open flames and hot surfaces.
14. Lay the USB cable so that it does not pose a tripping hazard.
15. Do not kink the USB cable or wrap it around the product. Do not lay the USB cable over sharp edges.

## **SAVE THESE INSTRUCTIONS**

## 2 OPERATING INSTRUCTIONS

### 2.1 Package contents/device parts (see Fig. A)

- |   |                     |
|---|---------------------|
| 1 Control button (⏻)                              | 7 Spindle           |
| 2 Charging port                                   | 8 LED lights        |
| 3 Cover   | 9 Ventilation holes |
| 4 Grinding discs<br>(Fine / Medium / Coarse), 3×* | 10 Dust chamber     |
| 5 Fixing ring                                     | 11 USB cable        |
| 6 Filter, 3×*                                     | 12 Cleaning brush   |

\* One piece is already attached to the product.

### 2.2 Checking the product and package contents

#### NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.
1. Take the product out of the packaging.
  2. Check to make sure that all parts are included (see **Fig. A**).
  3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

## 2.3 General information

**WARNING:** This product is not intended to be opened and/or repaired by the user and contains a battery that is not user-replaceable. Service/Repairs of the product must be performed only by an authorized person or at an authorized service center.

## 2.4 Reading and storing the user manual

This user manual accompanies the callus remover – VacuPedi (hereafter referred to as the “product”), and contains important information on setup and handling.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product.

Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

## 2.5 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

### **WARNING!**

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

### **NOTICE!**

This signal word warns of possible damage to property.

This signal symbol identifies useful additional information on handling and use.

Read the user manual.

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

This symbol identifies voltage and the current rating for direct current.

## 2.6 Proper use

The product is exclusively designed for removing calluses on the feet. It is only intended for private use and not suitable for commercial purposes. Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer or vendor cannot be held liable for damages or injury incurred through improper or incorrect use.

### **WARNING! RISK OF INJURY!**

- Only use the product on your feet to treat calloused skin. Do not use it on other parts of the body and on soft or thin parts of the skin.
- Only use the product for dry treatment of skin (do not use any additives such as oils or lotions).
- Do not continuously treat one spot of your skin for more than 2 to 3 seconds at a time.
- Do not use the product directly on or near open wounds and inflamed areas of the skin or to remove warts.
- Immediately stop using the product if you experience pain or skin irritation.

- Consult a doctor if you are not sure whether the product is suitable for you (e.g. you suffer from skin diseases or diabetes, you are pregnant, or you are using a pacemaker).
- For hygienic reasons, the product should only be used by one person.
- Keep the product away from hair, wigs and other long flexible objects, e.g. shoelaces, threads, textile fibres.

## 2.7 Operation

- Clean the product before first use as described in the chapter “Cleaning”.
- Before treating an area on your feet for the first time, perform a skin sensitivity test on a spot on your foot for maximum 3 seconds. Start your treatment only if there is no adverse reaction.

### 2.7.1 Attaching or replacing the grinding disc

The product includes 3 grinding discs **4**: Fine, Medium and Coarse. With particularly rough skin, we recommend you to start with the Coarse grinding disc.

1. Make sure the product is turned off.
2. Remove the cover **3** if it is attached to the product.
3. Remove the attached grinding disc from the product. If necessary, use the handle of the cleaning brush **12** to aid the removal (see step 1 in **Fig. B**).
4. Attach the desired grinding disc onto the spindle **7** (see step 2 in **Fig. B**).

You will hear a click sound when the grinding disc is secured into place.

### 2.7.2 Starting and ending the treatment

1. Make sure your feet are clean and completely dry.
2. Press the control button (⏻) **1** once.  
The product is now turned on at a low speed. The LED lights **8** will light up blue.

3. Press the control button (⏻) again if you want a higher speed setting.  
The product is now turned on at a high speed. The LED lights will blink blue.
4. Move the rotating grinding disc (4) over the parts of your skin you intend to treat. In the process, apply gentle pressure and only keep the grinding disc on a spot on your skin for two to three seconds at a time, then check the result (see **Fig. C**).  
Do not cover the ventilation holes (9) during operation to avoid overheating.  
The dust created during treatment will be sucked into the dust chamber (10).
5. When treatment is done, press the control button (⏻) a third time to turn the product off. The LED lights will also turn off.

## 2.8 Charging

- The product cannot be used during charging.
  - Only charge the product at temperatures between 0 and 40 °C.
  - Fully charge the product before first use and when the battery power is low.
  - It takes approx. 2.5 hours for a full charge.
  - When the battery power is low, the LED lights become red and start to blink during operation. Recharge the product.
1. Plug the smaller end of the USB cable (11) into the charging port (2).
  2. Plug the larger end of the USB cable into the USB port of a compatible USB power source (not included) (e.g. USB wall adapter, 5V DC).
  3. Plug the USB wall adapter into a mains socket for charging.

During charging, the LED lights **8** will become blue and start to blink slowly.

4. The LED lights will become solid blue when the product is fully charged. Unplug the USB wall adapter from the wall socket and disconnect the USB cable from the product.

### 3 Cleaning

#### NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners, brushes with metal bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers, and the like.

Clean the product after each use.

1. Remove the attached grinding disc **4** from the product. If necessary, use the handle of the cleaning brush **12** to aid the removal (see step 1 in **Fig. B**).
2. Turn the fixing ring **5** counterclockwise to detach it from the product (see **Fig. D**).
3. Remove the filter **6** from the dust chamber **10** (see **Fig. D**).

4. Clean these parts accordingly:

**Main body:** Wipe it off with a damp cloth. Then wipe it completely dry with a dry cloth.

**Dust chamber and grinding discs:** Use the cleaning brush to remove the dust accumulated on them.

**Filter, fixing ring and cover **3**:** Clean them under running water. Wipe the fixing ring and cover completely dry with a dry cloth and let the filter air dry

thoroughly. Replace the filter if it is too dusty and the grinding disc if it is worn out.

5. Reassemble the filter, fixing ring, grinding disc and cover after cleaning. Please note, that all parts need to be completely dry before reassembling.

You can purchase new grinding discs and filters on our website at [www.silkn.com](http://www.silkn.com) or [www.silkn.ca](http://www.silkn.ca).

## 4 Storage

All parts must be completely dry before being stored.

- Store all parts in a clean and dry place, avoiding heat and direct sunlight.
- Attach the cover **3** to protect the grinding disc **4** for storage.

## 5 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not start.	Make sure the product is properly charged.
The product does not charge.	Make sure the USB cable <b>11</b> is securely connected to the charging port <b>2</b> and the USB power source.

## 6 Technical data

### Product

Model:	NV8615D
Rated input:	5V DC
Wattage:	3.5 W
Battery type:	1× 18650 Lithium-Ion battery (rechargeable)
Battery voltage/capacity:	3.7 V DC / 2,000 mAh
Rotation speed settings:	2 (low/high)
Rotation speed:	1,350 rpm (low) 1,550 rpm (high)
Charging time:	2.5 hours
Use time:	Approx. 1 hour 50 min per full charge
Weight:	266 g
Dimensions (L × W × H):	180 × 64 × 89.5 mm

## 7 Disposal

### 7.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

## 7.2 Disposing of the product

### **Old appliances may not be disposed of in household waste!**

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or county.**

This ensures that old appliances are recycled in a professional manner and also rules out negative consequences for the environment.

### **Batteries and rechargeable batteries may not be disposed of with household waste!**

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances\* or not, to a collection point run by the communal authority or to a retailer, so that they can be disposed of in an environmentally friendly manner. Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

\* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

**Note:** To remove the rechargeable battery, dismantle the fan unit and remove the battery terminals from the circuit board. Then remove the cover from the end of the handle to remove the battery. Make sure the battery is not damaged during the removal. Please note, the cover will be damaged and the removal of the battery is not reversible.

## 8 Warranty

This product is covered by a 1-year warranty. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

## 9 Customer service

For more information about Silk'n products, visit your regional Silk'n website: [www.silkn.com](http://www.silkn.com) or [www.silkn.ca](http://www.silkn.ca). If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from the websites.

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© 2024 Silk'n Beauty Ltd. All rights reserved. Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.