# Silk'n®

## BRIGHT LUX

### FAC03

EN

DE

FR

NL

ES

IT

TR

PT

DA

NO

SV

CS

PL



LT

HU

RO

- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DE USUARIO
- IT MANUALE D'USO
- TR KULLANIM KILAVUZU
- PT MANUAL DO UTILIZADOR

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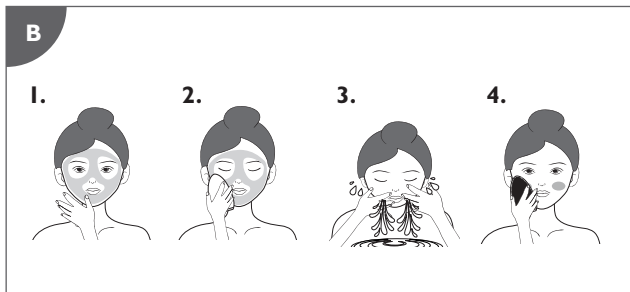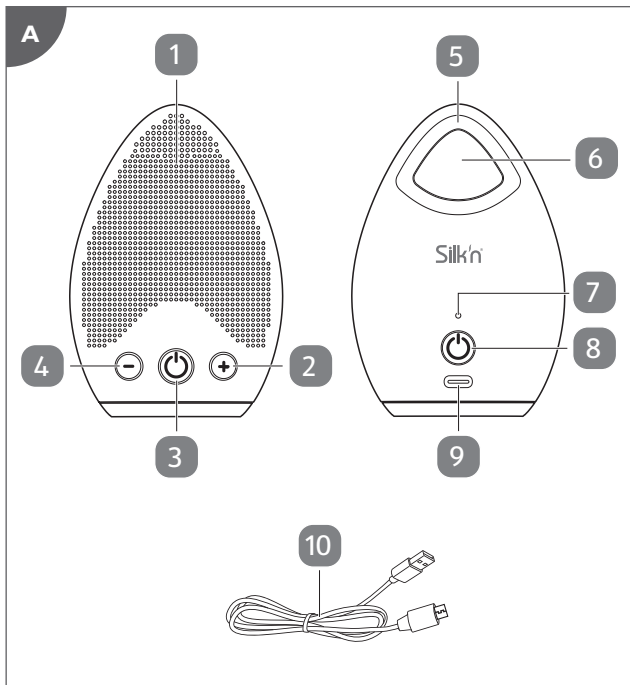## **Legal notice**

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silk'n Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silk'n Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.

Invention Works B.V.  
Galileïstraat 17, 3029 AL Rotterdam, The Netherlands  
[www.silkn.eu](http://www.silkn.eu) – [info@silkn.eu](mailto:info@silkn.eu)

# Table of contents

<b>1</b>	<b>Package contents/device parts</b>	<b>5</b>
<b>2</b>	<b>General information</b>	<b>5</b>
2.1	Reading and storing the user manual	5
2.2	Explanation of symbols	5
<b>3</b>	<b>Safety</b>	<b>6</b>
3.1	Proper use	6
3.2	Safety instructions	6
<b>4</b>	<b>Checking the product and package contents</b>	<b>10</b>
<b>5</b>	<b>Charging</b>	<b>11</b>
<b>6</b>	<b>Facial cleansing</b>	<b>11</b>
6.1	Tips on facial cleansing	11
6.2	Switching on/off and selecting an operation mode	12
6.3	Selecting a speed setting	12
6.4	How to use	12
<b>7</b>	<b>LED light, Infrared and cooling mode</b>	<b>13</b>
7.1	Tips on LED light, Infrared and cooling mode	13
7.2	Switching on/off and selecting an LED mode	14
7.3	How to use	14
7.4	Regimen	15
<b>8</b>	<b>Travel lock</b>	<b>15</b>
<b>9</b>	<b>Cleaning</b>	<b>15</b>
<b>10</b>	<b>Storage</b>	<b>16</b>
<b>11</b>	<b>Troubleshooting</b>	<b>16</b>
<b>12</b>	<b>Technical data</b>	<b>17</b>
<b>13</b>	<b>Disposal</b>	<b>17</b>
13.1	Disposing of the packaging	17
13.2	Disposing of the product	17
<b>14</b>	<b>Warranty</b>	<b>18</b>
<b>15</b>	<b>Customer service</b>	<b>18</b>

# I Package contents/device parts

- |                                |                                |
|--------------------------------|--------------------------------|
| 1 Facial brush                 | 6 Cooling surface              |
| 2 Increase button (+)          | 7 Cooling mode indicator light |
| 3 Cleansing button (⏻)         | 8 LED button (⏻)               |
| 4 Decrease button (-)          | 9 Charging port                |
| 5 LED light + Infrared surface | 10 USB cable                   |

## 2 General information

### 2.1 Reading and storing the user manual

This user manual accompanies this Bright Lux (hereafter referred to as the "product"), and contains important information on setup and handling.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product. The user manual is based on the standards and rules in force in the European Union and United Kingdom. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

### 2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

#### **WARNING!**

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

#### **NOTICE!**

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the United Kingdom.

This symbol identifies voltage and the current rating for direct current.

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

**IPX6**

The product is protected against powerful water jets.

## 3 Safety

### 3.1 Proper use

The product is exclusively designed for facial cleansing and skin improvement with LED light technology, Infrared and cooling. It is only intended for private use and not suitable for commercial purposes.

Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer cannot be held liable for damages or injury incurred through improper or incorrect use.

### 3.2 Safety instructions

#### **WARNING! RISK OF ELECTRIC SHOCK!**

- Only use the product with the supplied USB cable.
- Only connect the USB cable of the product to a suitable mains adaptor (IEC/EN 60335 certified SELV power supply) connected to an easily accessible socket so that you can

quickly disconnect it from the mains supply in the event of a problem.

- Do not use the product if it is damaged or if the USB cable is defective.
- Do not open the housing; instead, have a qualified professional perform repairs. Contact the customer service team for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection or incorrect operation.
- Only parts that comply with the original device data may be used for repairs. This product contains electrical and mechanical parts which are essential for providing protection against sources of danger.
- Do not immerse the USB cable in water or other liquids. Keep it away from water.
- Make sure the charging port is dry before charging.
- Always switch the product off when you are not using it, when you intend to clean it, or in the event of a malfunction. Always unplug the USB cable when you are not charging the product.
- Never touch the USB cable with wet or damp hands.
- Never insert any objects into the housing.
- Keep all parts of the product away from open flames and hot surfaces.

- Lay the USB cable so that it does not pose a tripping hazard.
- Do not pull the USB cable out of the mains adaptor by the cord; instead, always pull it out by the cable plug.
- Do not kink the USB cable and do not lay it over sharp edges.
- Only charge the product at temperatures between 0 and 40 °C.
- Only use the product indoors.

** WARNING! RISK OF INJURY!**

- The product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the product in a safe way and understand the hazards involved. Children must not play with the product. Cleaning and user maintenance must not be performed by children.
- Keep children under the age of 8 away from the product and USB cable.
- The product cannot be used by pregnant women, children, people with epilepsy, thyroid patients, people using thyroid medications, people who are photoallergic or use light allergy medication.
- Do not use the product on the eye lids,

inside the eye orbits, around areas that contain inflammatory acne, or areas that have rashes or wounds.

- The product cannot be used on skin diseases, wound ulceration, trauma and post-operative wounds.
- Do not focus the product on one spot on your skin during use.
- Consult your physician if you are uncertain whether you can use the product.
- Stop using the product immediately if you feel any discomfort, develop any abnormal reactions (such as itching or uncomfortable inflammation) or in case your skin becomes warm to the point of discomfort.
- Under abusive conditions, liquid may be ejected from the battery; avoid contact. Liquid ejected from the battery may cause irritation or burns.
- If battery fluid comes into contact with skin, eyes or mucous membrane, wash and rinse the affected areas immediately with clean water. Consult a physician if necessary.
- The batteries must not be disassembled, short-circuited or thrown into a fire.

## **NOTICE! RISK OF DAMAGE!**

- Do not use clay-based, silicone-based or grainy cleansers, exfoliators or scrubs with

the product as they may damage the soft silicone touch-points of the product.

- Avoid leaving the product in direct sunlight and never expose it to extreme heat or boiling water.
- Clean the product after each use to avoid the accumulation of grease and residues of face wash products.
- Stop using the product if its parts exhibit cracks or breaks or are deformed. Only replace damaged parts with corresponding original spare parts.
- Do not operate the product if it has been dropped or is damaged in any way. Have the product checked and/or repaired by qualified technicians if necessary.

#### **4 Checking the product and package contents**

##### **NOTICE! RISK OF DAMAGE!**

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.
  1. Take the product out of the packaging.
  2. Check to make sure that all parts are included (see **Fig. A**).
  3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

## 5 Charging

- Fully charge the product before first use and when the battery power is low.
- When the battery power is low, the product emits 2 beep sounds and the cleansing button flashes 10 times before switching off automatically. Recharge the product.
- The product cannot be used during charging.
- If the product is to be stored for an extended period, the battery level should be checked and the product should be charged regularly.

1. Plug the smaller end of the USB cable **10** into the charging port **9** of the product.
2. Plug the larger end of the USB cable into the USB port of a compatible USB power source (e.g. USB power adaptor).
3. Plug the USB power adaptor into a mains socket for charging. Charging starts automatically. The cleansing button ( **3**) flashes white during charging and stays lit when charging is done.
4. When charging is done, unplug the USB power adaptor from the mains socket and disconnect the USB cable from the charging port of the product.

## 6 Facial cleansing

### WARNING! RISK OF INJURY!

- Do not use the facial brush for more than 2 minutes per cleansing session and more than 2 times per day.

### 6.1 Tips on facial cleansing

- Clean the facial brush **11** before first use as described in the chapter "Cleaning".
- Before cleansing an area of the face for the first time, perform a skin sensitivity test on a spot. Start your cleansing session only if there is no adverse reaction.
- The facial brush is equipped with 2 operation modes: vibrating mode and pulsating mode. The vibrating mode uses uniform vibrations for deep, gentle cleansing. The pulsating mode uses alternating strong and weak vibrations for a massage-like soothing effect. There are also 7 speed settings for each mode.

## 6.2 Switching on/off and selecting an operation mode

1. Press the cleansing button (⏻) **3** once.  
The facial brush **1** is switched on in the vibrating mode.  
The cleansing button illuminates white during operation.
2. Press the cleansing button (⏻) a second time.  
The device is switched to the pulsating mode.
3. Press the cleansing button (⏻) a third time.  
The device is switched off.

## 6.3 Selecting a speed setting

- In the selected operation mode (vibrating or pulsating), press the increase button (+) **2** or decrease button (-) **4** to select a speed setting (1–7). The default level is 4.

## 6.4 How to use

1. Dampen your face and the facial brush **1** with lukewarm water. Apply your regular face wash product on your face (see step 1 in **Fig. B**).
2. Switch on the facial brush, select the vibrating mode and a speed setting (see instructions in the previous chapters).
3. Clean your face by moving the facial brush over the face in circular movements (see step 2 in **Fig. B**).
  - Use the thicker studs for areas that need precise cleansing (e.g. around the nose, eyes and chin) and use the thinner studs for the larger areas (e.g. forehead, cheeks and neck). We recommend cleansing each area for no more than 20 seconds.
  - Once the facial brush starts, you hear a beep sound every 30 seconds, hinting to you about the time that has passed. The device automatically switches off after 2 minutes.  
Note that whenever you change the operation mode or speed setting during the cleansing session, the timer will be reset to zero.
  - If you want to end the cleansing session earlier, press and hold the cleansing button (⏻) **3** until the facial brush switches off.
4. Thoroughly rinse your face and the facial brush with water and then wipe them dry with a clean, dry towel (see step 3 in **Fig. B**).
5. Apply your serum or day/night cream on your face.
6. Turn on the facial brush and select the pulsating mode and a speed setting (see instructions in the previous chapters).

7. Massage the face by moving the facial brush over the face in circular movements to stimulate blood circulation and promote the absorption of your skin care product (see step 4 in **Fig. B**).
  - The facial brush automatically switches off after 2 minutes.
  - If you want to end the cleansing session earlier, press and hold the cleansing button (⏻) until the facial brush switches off.
8. Thoroughly clean the facial brush after use (see chapter "Cleaning").

## 7 LED light, Infrared and cooling mode

### WARNING! RISK OF INJURY!

- Do not use the LED light + Infrared surface and cooling surface for more than 15 minutes per LED session.

#### 7.1 Tips on LED light, Infrared and cooling mode

- Clean the LED light + Infrared surface (5) and the cooling surface (6) before first use as described in the chapter "Cleaning".
- The LED light + Infrared surface is designed to recover your skin's glow and elasticity. You may also select the cooling mode, which helps to reduce puffy eyes for a brighter appearance.
- For skin improvement, you can choose from 5 modes with various benefits:

Mode	Best for
Red LED light + Infrared	Reduces fine lines and wrinkles, stimulates collagen and helps to improve blood circulation
Blue LED light + Infrared	Reduces acne-related redness and imperfections for clearer skin
Yellow LED light + Infrared	Improves skin hydration and soothes dry and sensitive skin
Purple LED light + Infrared	Anti-ageing and reduction of skin imperfections to improve the skin's vitality and radiance

Mode	Best for
Red LED light + Infrared + cooling	Reduces puffy eyes for a brighter appearance

## 7.2 Switching on/off and selecting an LED mode

1. Press the LED button (⏻) **8** once.  
The LED light + Infrared surface **5** is switched on in the **red LED** mode.
2. Press the LED button (⏻) a second time.  
The device is switched to the **blue LED** mode.
3. Press the LED button (⏻) a third time.  
The device is switched to the **yellow LED** mode.
4. Press the LED button (⏻) a fourth time.  
The device is switched to the **purple LED** mode.
5. Press the LED button (⏻) a fifth time.  
The device is switched to the **red LED + cooling mode**.  
In this mode, the LED light + Infrared surface illuminates red and the cooling surface **6** cools down to around 6 °C. The cooling mode indicator light **7** also turns on.
6. Press the LED button (⏻) a sixth time.  
The device is switched off.

**Tip:** You can directly select the red LED + cooling mode without pressing the LED button several times. To do so, press and hold the LED button (⏻) in any LED mode until the red LED + cooling mode is activated.

## 7.3 How to use

1. Thoroughly clean and dry your face. Make sure it is free from any make-up or dirt.
2. Turn on the LED light + Infrared surface **5** and select an LED mode (see instructions in the previous chapter).
3. Move the LED light + Infrared surface over your face in circular movements.
  - Once the LED light + Infrared surface starts, you hear a beep sound every 30 seconds, hinting to you about the time that has passed. The device automatically switches off after 10 minutes (in any LED mode) or 2 minutes (red LED + cooling mode). Note that whenever you change the mode during the LED session, the timer will be reset to zero.

- If you want to end the LED session earlier, press the LED button (⏻) **8** repeatedly until the device switches off.
- 4. After the LED session, apply a serum and/or day/night cream on your face.

## 7.4 Regimen

If you are using the LED light + Infrared surface **5** for the first time, follow the schedule below:

First week	Use once every 3 days, max. 10 minutes at a time
Second week	Use once every 2 days, max. 10 minutes at a time
Third week	Use once every 2 days, max. 10–15 minutes at a time
Fourth week	Use once every day, max. 10–15 minutes at a time

As soon as results show, adjust the frequency to once or twice per week.

## 8 Travel lock

- Press and hold the increase button (+) **2** and decrease button (–) **4** at the same time until you hear 2 beep sounds. The buttons are now locked to prevent accidental activation during travel.
- Repeat the step above to unlock the buttons.

## 9 Cleaning

### **WARNING! RISK OF ELECTRIC SHOCK!**

- Always unplug the product from the power source before cleaning.

### **NOTICE! RISK OF DAMAGE!**

- Do not use any aggressive cleaners, brushes with metal bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers, and the like.

- Do not use hot water for cleaning.

1. Clean the product with lukewarm water and soap.
2. Thoroughly rinse the product with water after cleaning and wipe it dry with a dry cloth.

## 10 Storage

All parts must be completely dry before being stored.

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product so that it is not accessible to children.

## 11 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not operate.	<ul style="list-style-type: none"> <li>• Make sure the product is properly charged.</li> <li>• Make sure the product is not connected to a USB power source. The product does not work while charging.</li> </ul>
The product does not charge.	<ul style="list-style-type: none"> <li>• Make sure the USB cable <b>10</b> is securely connected to the charging port <b>9</b> and the USB power source.</li> </ul>

## 12 Technical data

Model:	FAC03
Rated input:	5 V === 1 A
Working voltage:	3.7 V
Battery capacity:	1,500 mAh, 5.55 Wh (lithium)
Number of LED lights:	5
Charging time:	2–3 hours
Weight:	197.7 g
IP code:	IPX6
Dimensions (L × W × H):	110.5 × 75 × 44 mm

## 13 Disposal

### 13.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

### 13.2 Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

#### **Old appliances may not be disposed of in the household waste!**

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or county**. This ensures that old appliances are recycled in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

## **Batteries and rechargeable batteries may not be disposed of with household waste!**

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances\* or not, to a collection point run by the communal authority or to a retailer, so that they can be disposed of in an environmentally friendly manner.

Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

\* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

## **14 Warranty**

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

## **15 Customer service**

For more information about the products of Silk'n visit your regional Silk'n website: [www.silkn.eu](http://www.silkn.eu). If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: [www.silkn.eu](http://www.silkn.eu).

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Copyright 2023 © Silk'n Beauty Ltd. Toate drepturile rezervate.
© 2023 Silk'n Beauty Ltd. All rights reserved. Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.
Silk'n UK Ltd.