# Silk'n™

Beautiful Technology


## FaceTite

H2111 / H2112

- DE BEDIENUNGSANLEITUNG
- EN USER MANUAL
- FR MODE D'EMPLOI
- BG РЪКОВОДСТВО ЗА УПОТРЕБА
- HU HASZNÁLATI UTASÍTÁS
- RO INSTRUCȚIUNI DE UTILIZARE
- CS NÁVOD K POUŽITÍ
- SL PRIROČNIK ZA UPORABNIKA
- PL INSTRUKCJA OBSŁUGI
- RU РУКОВОДСТВО ПОЛЬЗОВАТЕЛЯ

# Silk'n FaceTite

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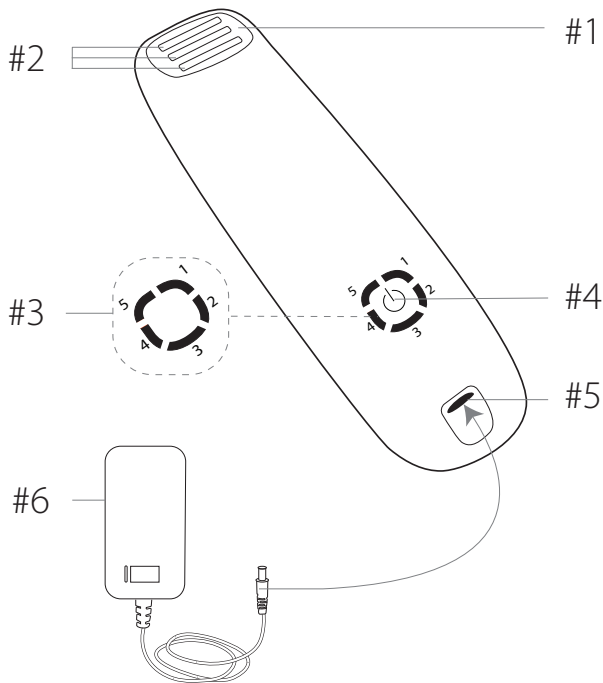# Table of contents

<b>1. Warnings &amp; contraindications.....</b>	<b>22</b>
<b>2. Signs &amp; symbols .....</b>	<b>26</b>
<b>3. Safety features .....</b>	<b>27</b>
<b>4. Device controls &amp; settings.....</b>	<b>27</b>
<b>5. Get to know Silk'n FaceTite .....</b>	<b>28</b>
5.1 Package contents.....	28
5.2 Device parts.....	28
5.3 Device description .....	29
5.4 Intended use.....	29
5.5 Reading & storing the user manual.....	29
<b>6. Silk'n FaceTite treatments .....</b>	<b>30</b>
6.1 Treatment areas.....	30
6.2 Treatment plan.....	31
<b>7. What to expect with Silk'n FaceTite .....</b>	<b>32</b>
<b>8. Get started!.....</b>	<b>33</b>
8.1 Before first use.....	33
8.2 Treatment steps .....	33
<b>9. Cleaning &amp; maintenance .....</b>	<b>35</b>
<b>10. Storage.....</b>	<b>35</b>
<b>11. Troubleshooting .....</b>	<b>35</b>
<b>12. Consumables .....</b>	<b>36</b>
<b>13. Customer service.....</b>	<b>36</b>
<b>14. Warranty .....</b>	<b>37</b>
<b>15. Disposal .....</b>	<b>37</b>
<b>16. Specifications .....</b>	<b>38</b>
<b>17. Legal notice .....</b>	<b>39</b>

## 1. Warnings & contraindications

---

A statement or event that indicates a possible or impending danger, problem, or other unpleasant situation when using the device.

### **Electricity & safety warnings**

- Only authorized personnel are permitted to perform repairs. Do not attempt to repair or open sealed parts of your device. You may expose yourself to dangerous components and cause serious harm to yourself. It will also void your warranty.
- Do not place the device or adapter where it can fall or be pulled into a tub, sink, water, or any other liquid. The device and adapter are not waterproof.
- Do not reach for the device or adapter if they have fallen into water. Unplug immediately before retrieving either of them. Failure to unplug the device and/or the adapter may result in an electric shock!
- Never use the adapter if it has a damaged cable or plug, it isn't working properly, or has been submerged in water. If that happens, unplug the adapter immediately.

- ⦿ The device is only to be used with the adapter provided.
- ⦿ Do not use if the parts or accessories are not supplied by the manufacturer as defined in this manual, or if the device is damaged in any way, does not appear to be working properly, or if you see or smell smoke. In such cases, stop using the device, and contact customer service.
- ⦿ Do not clean the device or any of its parts in the dishwasher.
- ⦿ Children shall not play with the device. Cleaning and user maintenance shall not be made by children.
- ⦿ Children must not play with the packaging. They may injure themselves with the packaging material, or swallow small parts and choke.
- ⦿ Do not dispose of the device with household waste! The law requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices.

## Contraindications

This device including all its parts must not be used by the following people: children below the age of 18, adults with limited physical, sensory or psychological capacities, those lacking experience or knowledge in how to use the device in a safe way, or those who do not understand the hazards involved.

### Do not use if you:

- have a pacemaker, internal defibrillator, or any other active electrical implant anywhere in the body.
- have poorly controlled endocrine disorders, such as diabetes and thyroid conditions.
- have abnormal skin conditions caused by diabetes, systemic or other metabolic diseases.
- have a current history or condition of cancer, or pre-malignant moles.
- have had any surgery in the treated area within 3 months prior to treatment.
- suffer from epilepsy.

- have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, or lupus, etc.
- are pregnant or nursing.

### **Do not treat areas that have:**

- a permanent implant such as metal plates and screws, silicone implants, or areas with injected fills or toxins.
- active eczema, psoriasis, lesions, open wounds, open infections (cold sores) or abnormal skin conditions caused by systemic or metabolic diseases (diabetes for example).
- traumatized or injured skin.
- had recent exposure to the sun without protection (wait 1 day); for tanned or sunburnt skin, wait 3 days.

### **Consult your physician if you have:**

- couperose or rosacea.
- a history of diseases stimulated by heat, such as recurrent herpes simplex in the treatment area. If so, you may treat only following a prophylactic regime.

- ⦿ received radiation therapy or chemotherapy treatments within the past 3 months, or have areas of potential skin malignancies.
- ⦿ used steroids within 3 months prior to treatment.
- ⦿ used physician-prescribed drugs on your skin within the past 6 months.
- ⦿ any other condition which, in the opinion of your physician or dermatologist, would make it unsafe for you to be treated.

## 2. Signs & symbols

	Important information about using or maintaining the device.
	Products marked with this symbol fulfil the requirements of the EU directives.
	Class II. The adapter is protected against electric shock by means of double insulation.
IP32	The device is protected against the ingress of objects such as fingers, and against drops of water.
	Do not use the device in the vicinity of bathtubs, showers, washbasins or other vessels containing water.
	Do not dispose with household waste. Dispose in accordance with your country's laws and help protect the environment. In accordance with the WEEE Directive.

### 3. Safety features

Silk'n FaceTite protects your skin from overheating and has these safety features:

- **Thermal sensor and heat stabilizer:** These measure the skin temperature, stabilize the heat energy and prevent overheating. When the skin temperature reaches 42°C - 43°C, the stabilizer will automatically stop emitting heat, until the skin temperature lowers enough for the stabilizer to resume heating.
- **Movement sensor:** Silk'n FaceTite will only emit RF heat when you are moving the device across your skin. If you keep the device on one part of your face for too long, the device will stop emitting heat automatically. This is a safety device.

### 4. Device controls & settings

Button	Action	Function / Mode
N/A	The device is connected to the electricity supply via the <b>adapter</b> (#6).	Standby mode  All 5 <b>energy level indicator lights</b> (#3) cycle continuously.
Control button (#4)  	Press x 1	Energy level 1 - 1 light on
	Press x 2	Energy level 2 - 2 lights on
	Press x 3	Energy level 3 - 3 lights on
	Press x 4	Energy level 4 - 4 lights on
	Press x 5	Energy level 5 - 5 lights on
	Press x 6	Goes back to level 1
	Long press	Goes back to standby mode

Button	Action	Function / Mode
Energy level indicator lights (#3) 	Automatically lit to show the energy level selected.	There is one light for each energy level. Level 1 is the lowest energy level, and level 5 is the highest.

## 5. Get to know Silk'n FaceTite

### 5.1 Package contents

- Silk'n FaceTite device
- Adapter
- Silk'n Slider Gel
- Storage case
- User manual
- Warranty card

### 5.2 Device parts

The device has a **treatment surface** (#1), **electrodes** (#2), and a **control button** (#4) with surrounding **energy level indicator lights** (#3).

Refer to the illustration on page (i) at the front of this manual to see the location of each feature.

- Treatment surface (#1)
- 3 Electrodes (#2)
- 5 Energy level indicator lights (#3)
- Control button (#4)
- Device socket (#5)
- Adapter with cable (#6)

### 5.3 Device description

Silk'n FaceTite is an anti-aging face treatment device which utilizes HT™ technology. This technology combines three kinds of energy: Infrared (IR) Heat Energy that affects the upper layer of skin for improved skin texture, LED Light Energy that reaches into the dermis for rejuvenation and collagen restoration, and Bi-Polar Radio Frequency (RF) Energy that penetrates deep into the sub-dermal fat for deep tissue tightening and wrinkle reduction.

### 5.4 Intended use

Silk'n FaceTite is intended for improving skin texture, tightening the skin, restoring collagen, and reducing wrinkles.

The device is for self-treatment only. It is intended for private household use and is not suitable for commercial purposes. Only use the device as described in this user manual. Any other use is considered improper, and may result in damage to property or persons.

The manufacturer or vendor cannot be held liable for damages or injury incurred through improper or incorrect use. For additional information see our website [www.silk'n.eu](http://www.silk'n.eu).

### 5.5 Reading & storing the user manual

This user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws.


Before first use, please read this user manual in full, and keep it for future use. We recommend you refamiliarise yourself with the instructions before each treatment.

## 6. Silk'n FaceTite treatments

---

### 6.1 Treatment areas

There are 6 treatment areas for the face and neck.

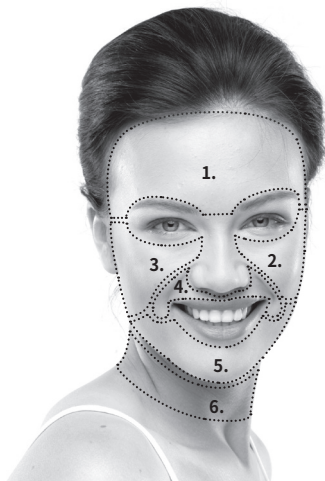Treatment Area	Purpose
Forehead, excluding the eye brows	Tightening the skin
Left cheek, including the area next to and below the eye	Treating fine lines and building volume
Right cheek, including the area next to and below the eye	
Upper lip, including the nasolabial fold	Tightening and eliminating wrinkles
Chin, including the jawline	Tightening and V-Shaping
Neck	Smoothing the skin


## 6.2 Treatment plan

Treat each area separately for 15 minutes. You can treat as many areas as you like in one session.

For maximum effectiveness, we recommend the following:

- Make sure you have enough gel on the **treatment surface** (#1) at all times during the treatment.
- Do not leave the device on one spot of skin. Keep moving the device covering the entire treatment area. This will prevent overheating of the treatment area, and reduce the possibility of any adverse effect.
- Make sure the 3 **electrodes** (#2) are always in full contact with the skin.
- If the energy level feels uncomfortable, reduce it to a level that is comfortable for you.
- Only change energy levels when the **treatment surface** (#1) is **not** in contact with the skin.
- Treat each area for at least 10 minutes at energy level 5.

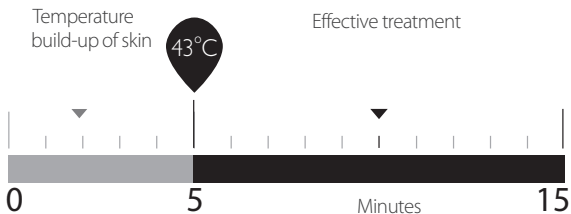- Use the device for 2 sessions a week for 10 weeks. You will begin to see visible results after 5 weeks / 10 treatments.

- Treat once every 2 weeks after that, to maintain results.

## 7. What to expect with Silk'n FaceTite

If you follow the instructions fully and pay attention to all warnings and contraindications, complications and side effects are uncommon. However, every cosmetic treatment, including those designed for home use, involves some degree of risk. Therefore it is important that you understand and accept the risks and complications that may occur.

- Excessive edema (swelling) and erythema (redness) are rare reactions to the treatment that, in most cases, subside within a few days. In case these signs appear, immediately use a damp cool cloth to cool the area for at least 10 minutes and apply an anti-burn cream.
- In rare cases the user may experience dilation of the superficial blood capillaries. This may result in a temporary redness for a few days. Also, some of the capillaries may present red spots that will also reduce with time. If these signs appear, we recommend that you cool the area with chilled water immediately.
- ❗ Consult your physician if you experience any persistent symptoms including pain, skin redness, swelling or a burning sensation.

## 8. Get started!

---

### 8.1 Before first use

- Take the device out of the packaging, and check that all parts are included.
- If any items are missing or appear to be damaged, contact customer service.
- Before first use, remove the protective film that covers the electrodes.

### 8.2 Treatment steps

1. Remove the device and the other components from the box.
2. Remove the protective film from the **treatment surface** (#1).
3. Plug the cable of the **adapter** (#6) into the **device socket** (#5). Then connect the **adapter** (#6) into the electricity supply. All 5 **energy level indicator lights** (#3) will cycle continuously. The device is now in standby mode.
4. Before you start, make sure that your skin is clean and dry, and that you have removed all jewellery.
5. Press the **control button** (#4). The device starts at energy level 1 automatically, and one **energy level indicator light** (#3) will be on.
6. Apply a small amount of Silk'n Slider Gel to the **treatment surface** (#1).

Silk'n Slider Gel (or any other water-based clear gel) facilitates smooth movement of the **treatment surface** (#1) over the skin, and enables conductivity. **To ensure safe and efficient treatment, always start by applying Silk'n Slider Gel and apply more gel if needed.**

- EN
- Using energy level 1, place the device on the area you want to treat and make sure that it is in full contact with the skin. You will feel a mild sensation of warmth.
  - Move the device gently in linear massage-like motions across the area you wish to treat. Progress slowly and carefully. Keep moving the device ensuring that the 3 **electrodes** (#2) stay in good contact with the skin as you treat the area.
  - Increase the energy level gradually, by pressing the **control button** (#4). Using the highest energy level will get the best results, but it's important that you set the energy level that is most comfortable for you.
  - With the device constantly in motion, treat the area for at least 15 minutes in total. You will feel the skin temperature increase as you increase the energy level. During the treatment, your skin will get a little red. Don't worry; this is normal.

Stop immediately if you feel any discomfort.

- After your treatment, long press the **control button** (#4) to put the device into standby mode.
- Remove any residual Silk'n Slider Gel from the skin, then rinse the skin with lukewarm water. Dry the skin with a towel, and apply a facial cream to sooth and moisturize it.

## 9. Cleaning & maintenance

---

Clean Silk'n FaceTite after each treatment.

Never immerse the device or any of its parts in water!

- Put the device in standby mode, then unplug it.
- Use absorbent paper (such as toilet paper or wipes) to remove the gel residue off the **treatment surface** (#1). Make sure you remove all gel residue to prevent the gel from solidifying or damaging the device.
- Use a clean damp cloth to clean the device and the **treatment surface** (#1). You may use cleaning fluids with up to 70% alcohol, but never use acetone or other abrasive liquids, because they will damage your device.
- Finally, wipe the device and the **treatment surface** (#1) dry with a clean dry cloth.

## 10. Storage

---

After cleaning we recommend you keep the device in the storage case. Make sure that the device is completely dry before storage.

## 11. Troubleshooting

---

Do not attempt to open or repair your device. You may expose yourself to dangerous electrical components and may cause serious harm to yourself. It will also void the warranty.

Problem	Checks
The device does not start.	Make sure the <b>adapter</b> (#6) is properly connected to the device, that it is plugged into the electricity supply, and that you have pressed the <b>control button</b> (#4).
I can't feel any heat from the device even though I pressed the <b>control button</b> (#4).	Heat will only be emitted if the <b>treatment surface</b> (#1) is in full contact with the skin, and you have applied sufficient gel before and during your treatment. This is a safety feature. Gel is essential for the device to work.
One or more lights blink continuously.	The device has malfunctioned. Unplug the device and restart. If the problem continues, contact customer service.
The device has stopped working and is in standby mode.	Press the <b>control button</b> (#4) again to restart the device.

For problems or queries, please contact customer service.

## 12. Consumables

---

You can purchase more Silk'n Slider Gel from our website: [www.silkn.eu](http://www.silkn.eu).

## 13. Customer service

---

For more information see our website [www.silkn.eu](http://www.silkn.eu). This manual is also available as a PDF download from the website. Contact Silk'n Customer Service to report unexpected operations or events, or any other problem with your device. Our contact details are:

Country	Service number	Email
Belgium	0900-25006	serviceconsommateurbe@silkn.eu
France	0891-655557	serviceconsommateurfr@silkn.eu
Germany	089 51 23 44 23	kundenservicede@silkn.eu
Italy		servizioclientiit@silkn.eu
Netherlands	0900-2502217	klantenservicenl@silkn.eu
Spain	900 823 302	servicioalcliente@silkn.eu
UK	0906-2130009	customerareuk@silkn.eu
Other countries	+31 (0)180-330 550	info@silkn.eu

## 14. Warranty

---

Please refer to the separate Product Warranty booklet for complete information.

## 15. Disposal

---

**Packaging:** Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.

**Device:** Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. This ensures that the device is recycled professionally, and prevents the release of harmful substances\* into the environment. Make sure that the device is in an uncharged state before disposal.

\* Products labelled with Cd = cadmium, Hg = mercury, Pb = lead.

## 16. Specifications

Model		H2111 / H2112
Technology		Home Tightening™ (HT™)
Wavelength red light		620nm±20nm
Wavelength IR		850nm±20nm
Radio frequency		1MHz ±30%; 12W max
Operation and safety		Heat stabilizer to prevent skin overheating
System weight		125g
Temperature	Operating	10 to 35°C
	Storage	-40 to +70°C
Relative humidity	Operating	30 to 75%rH
	Storage	10 to 90%rH
Atmospheric pressure	Operating	700 to 1060hPa
	Storage	500 to 1060hPa
Adapter model		KSA24A1200150HE
- Input		100-240V; 50-60Hz; 0.5A
- Output		12Vdc; 1.5A

## 17. Legal notice

---

Copyright 2019 © Home Skinovations Ltd. All rights reserved.

Home Skinovations Ltd reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Home Skinovations Ltd is believed to be accurate and reliable at the time of publication. However, Home Skinovations Ltd assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Home Skinovations Ltd. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Home Skinovations Ltd. Data is subject to change without notification. Home Skinovations Ltd has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Home Skinovations Ltd. Specifications are subject to change without notice. Silk'n and the Silk'n logo are registered trademarks of Home Skinovations Ltd, Tavor Building, Shaar Yokneam P.O.Box 533, Yokneam 2069206, ISRAEL.

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Copyright 2019 © Home Skinovations Ltd. Alle Rechte vorbehalten.
Copyright 2019 © Home Skinovations Ltd. Minden jog fenntartva.
Copyright 2019 © Home Skinovations Ltd. Toate drepturile rezervate.
Copyright 2019 © Home Skinovations Ltd. Všechna práva vyhrazena.
Copyright 2019 © Home Skinovations Ltd. Wszelkie prawa zastrzeżone.
Invention Works B.V.