Silk'n®

VACUPEDI

FTC04

- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DE USUARIO
- PT MANUAL DO UTILIZADOR

- GR ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ
- TR KULLANIM KILAVUZU
- IT MANUALE D'USO
- NO BRUKERHÅNDBOK
- SV ANVÄNDARMANUAL
- DA BRUGERMANUAL

EN

DE

FR

NL

ES

PT

GR





TR

IT

NO

SV

DA

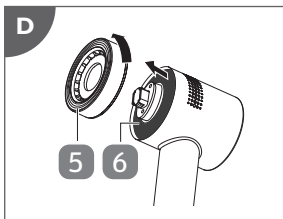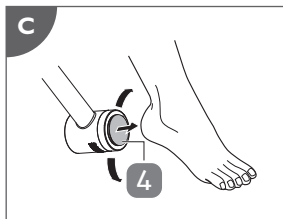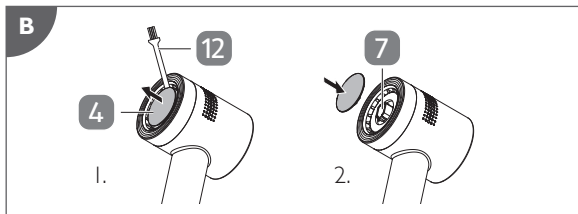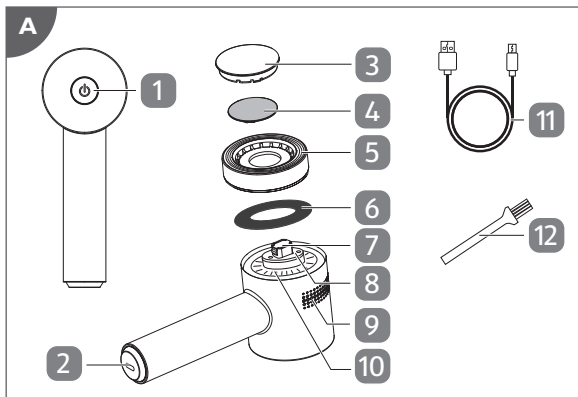Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silkn Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silkn Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silkn Beauty Ltd.

Invention Works B.V.
Galileïstraat 17, 3029 AL Rotterdam, The Netherlands
www.silkn.eu – info@silkn.eu

Table of contents

1	Package contents/device parts	5
2	General information	5
2.1	Reading and storing the user manual	5
2.2	Explanation of symbols	6
3	Safety	7
3.1	Proper use	7
3.2	Safety instructions	7
4	Checking the product and package contents	11
5	Charging	12
6	Operation	13
6.1	Attaching or replacing the grinding disc	13
6.2	Starting and ending the treatment	13
7	Cleaning	14
8	Storage	15
9	Troubleshooting	15
10	Technical data	16
11	Disposal	16
11.1	Disposing of the packaging	16
11.2	Disposing of the product	17
12	Warranty	18
13	Customer service	18

I Package contents/device parts

- | | |
|---|---------------------|
| 1 Control button (⏻) | 7 Spindle |
| 2 Charging port | 8 LED lights |
| 3 Cover | 9 Ventilation holes |
| 4 Grinding discs
(Fine / Medium / Coarse), 3×* | 10 Dust chamber |
| 5 Fixing ring | 11 USB cable |
| 6 Filter, 3×* | 12 Cleaning brush |

* One piece is already attached to the product.

2 General information

2.1 Reading and storing the user manual

This user manual accompanies the callus remover – VacuPedi (hereafter referred to as the “product”), and contains important information on setup and handling.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product.

The user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

WARNING!

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

NOTICE!

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the United Kingdom.

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

This symbol identifies voltage and the current rating for direct current.

3 Safety

3.1 Proper use

The product is exclusively designed for removing calluses on the feet. It is only intended for private use and not suitable for commercial purposes. Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer or vendor cannot be held liable for damages or injury incurred through improper or incorrect use.

3.2 Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

-  Keep the product and USB cable away from water or liquids. Do not use near bathtubs, showers, basins or other vessels containing water. They are not waterproof.
- Only use the product indoors.
- Only connect the USB cable of the product to a suitable mains adaptor (IEC/EN 60335 certified SELV power supply), connected to an easily accessible socket so that you can quickly disconnect it from the mains

supply in the event of a problem.

- Only use the product with the supplied USB cable.
- Do not use the product or USB cable if they are damaged or defective. The manufacturer or customer service must replace the damaged or defective parts to avoid risks.
- Do not open the housing and never insert any objects into the housing.
- Always switch the product off when you are not using it, when you intend to clean it, or in the event of a malfunction. Always unplug the USB cable when you are not charging the product.
- Never touch the USB cable with wet or damp hands.
- Never move, pull, or carry the product by the USB cable.
- Keep the product and USB cable away from open flames and hot surfaces.
- Lay the USB cable so that it does not pose a tripping hazard.

- Do not kink the USB cable or wrap it around the product. Do not lay the USB cable over sharp edges.

 **WARNING! RISK OF INJURY!**

- The product can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the product in a safe way and understand the hazards involved. Children must not play with the product. Cleaning and user maintenance must not be performed by unsupervised children.
- Keep children under the age of 8 away from the product and USB cable.
- Make sure that children do not play with the plastic wrapping. They may get caught in it when playing and suffocate.
- Only use the product on your feet to treat calloused skin. Do not use it on

other parts of the body and on soft or thin parts of the skin.

- Only use the product for dry treatment of skin (do not use any additives such as oils or lotions).
- Do not continuously treat one spot of your skin for more than 2 to 3 seconds at a time.
- Do not use the product directly on or near open wounds and inflamed areas of the skin or to remove warts.
- Immediately stop using the product if you experience pain or skin irritation.
- Consult a doctor if you are not sure whether the product is suitable for you (e.g. you suffer from skin diseases or diabetes, you are pregnant, or you are using a pacemaker).
- For hygienic reasons, the product should only be used by one person.
- Keep the product away from hair, wigs and other long flexible objects, e.g. shoelaces, threads, textile fibres.
- Only charge the product at temperatures between 0 and 40 °C.

EN NOTICE! RISK OF DAMAGE!

- Avoid leaving the product in direct sunlight and never expose it to extreme heat.
- Do not place the product on the edge of a surface to avoid dropping.
- Only use the product with the original accessories.
- Never fill the housing of the product with liquid.
- Do not cover the ventilation holes during operation to avoid overheating.

4 Checking the product and package contents

NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.

1. Take the product out of the packaging.
2. Check to make sure that all parts are included (see **Fig. A**).
3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

5 Charging

- Fully charge the product before first use and when the battery power is low.
 - When the battery power is low, the LED lights become red and start to blink during operation. Recharge the product.
 - It takes approx. 2.5 hours for a full charge.
 - The product cannot be used during charging.
1. Plug the smaller end of the USB cable **11** into the charging port **2**.
 2. Plug the larger end of the USB cable into the USB port of a compatible USB power source (not included) (e.g. USB wall adapter, 5V DC).
 3. Plug the USB wall adapter into a mains socket for charging.
During charging, the LED lights **8** will become blue and start to blink slowly.
 4. The LED lights will become solid blue when the product is fully charged. Unplug the USB wall adapter from the wall socket and disconnect the USB cable from the product.

6 Operation

- Clean the product before first use as described in the chapter “Cleaning”.
- Before treating an area on your feet for the first time, perform a skin sensitivity test on a spot on your foot for maximum 3 seconds. Start your treatment only if there is no adverse reaction.

6.1 Attaching or replacing the grinding disc

The product includes 3 grinding discs (4): Fine, Medium and Coarse. With particularly rough skin, we recommend you to start with the Coarse grinding disc.

1. Make sure the product is turned off.
2. Remove the cover (3) if it is attached to the product.
3. Remove the attached grinding disc from the product. If necessary, use the handle of the cleaning brush (12) to aid the removal (see step 1 in **Fig. B**).
4. Attach the desired grinding disc onto the spindle (7) (see step 2 in **Fig. B**).
You will hear a click sound when the grinding disc is secured into place.

6.2 Starting and ending the treatment

1. Make sure your feet are clean and completely dry.
2. Press the control button (⏻) (1) once.
The product is now turned on at a low speed. The LED lights (8) will light up blue.
3. Press the control button (⏻) again if you want a higher speed setting.
The product is now turned on at a high speed. The LED lights will blink blue.
4. Move the rotating grinding disc (4) over the parts of your skin you intend to treat. In the process, apply gentle pressure and only keep the grinding disc on a

spot on your skin for two to three seconds at a time, then check the result (see **Fig. C**).

Do not cover the ventilation holes (9) during operation to avoid overheating.

The dust created during treatment will be sucked into the dust chamber (10).

5. When treatment is done, press the control button (⏻) a third time to turn the product off. The LED lights will also turn off.

7 Cleaning

NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners, brushes with metal bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers, and the like.

Clean the product after each use.

1. Remove the attached grinding disc (4) from the product. If necessary, use the handle of the cleaning brush (12) to aid the removal (see step 1 in **Fig. B**).
2. Turn the fixing ring (5) counterclockwise to detach it from the product (see **Fig. D**).
3. Remove the filter (6) from the dust chamber (10) (see **Fig. D**).
4. Clean these parts accordingly:
Main body: Wipe it off with a damp cloth. Then wipe it completely dry with a dry cloth.
Dust chamber and grinding discs: Use the cleaning brush to remove the dust accumulated on them.

Filter, fixing ring and cover 3: Clean them under running water. Wipe the fixing ring and cover completely dry with a dry cloth and let the filter air dry thoroughly. Replace the filter if it is too dusty and the grinding disc if it is worn out.

5. Reassemble the filter, fixing ring, grinding disc and cover after cleaning. Please note, that all parts need to be completely dry before reassembling.

You can purchase new grinding discs and filters on our website at “www.silkn.eu”.

8 Storage

All parts must be completely dry before being stored.

- Store all parts in a clean and dry place, avoiding heat and direct sunlight.
- Attach the cover 3 to protect the grinding disc 4 for storage.

9 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not start.	Make sure the product is properly charged.
The product does not charge.	Make sure the USB cable is securely connected to the product and the USB port of a compatible USB power source.

10 Technical data

Product

Model:	FTC04
Rated voltage:	3.7 V
Wattage:	7 W
Battery capacity:	2,000 mAh
Speed settings:	2 (low/high)
Rotation speed:	1,350 rpm (low) 1,550 rpm (high)
Charging time:	2.5 hours
Use time:	Approx. 1 hour 50 min per full charge
Weight:	241 g
Dimensions (L × W × H):	180 × 64 × 89.5 mm

11 Disposal

11.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

11.2 Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

Old appliances may not be disposed of in household waste!

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or county.** This ensures that old appliances are recycled in a professional

manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

Batteries and rechargeable batteries may not be disposed of with household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances* or not, to a collection point run by the communal authority or to a retailer, so that

they can be disposed of in an environmentally friendly manner. Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

Note: To remove the rechargeable battery, dismantle the fan unit and remove the battery terminals from the circuit board. Then remove the cover from the end of the handle to remove the battery. Make sure the battery is not damaged during the removal. Please note, the cover will be damaged and the removal of the battery is not reversible.

12 Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

13 Customer service

For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: www.silkn.eu.

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Copyright 2024 © Silk'n Beauty Ltd. Todos os direitos reservados.
Copyright 2024 © Silk'n Beauty Ltd. Με την επιφύλαξη παντός δικαιώματος.
Copyright 2024 © Silkn Beauty Ltd. Alle rettigheder forbeholdes.
Silk'n UK Ltd.