Silk'n®

LED EMS FACE MASK



FAC06

- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DE USUARIO
- IT MANUALE D'USO
- TR KULLANIM KILAVUZU
- PT MANUAL DO UTILIZADOR

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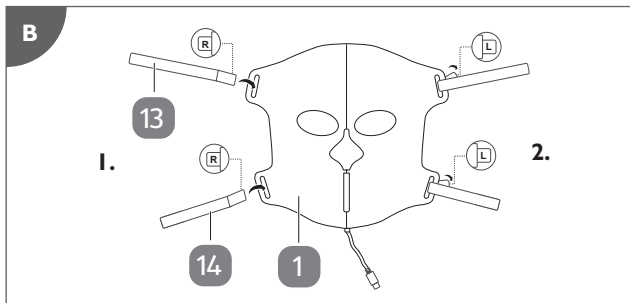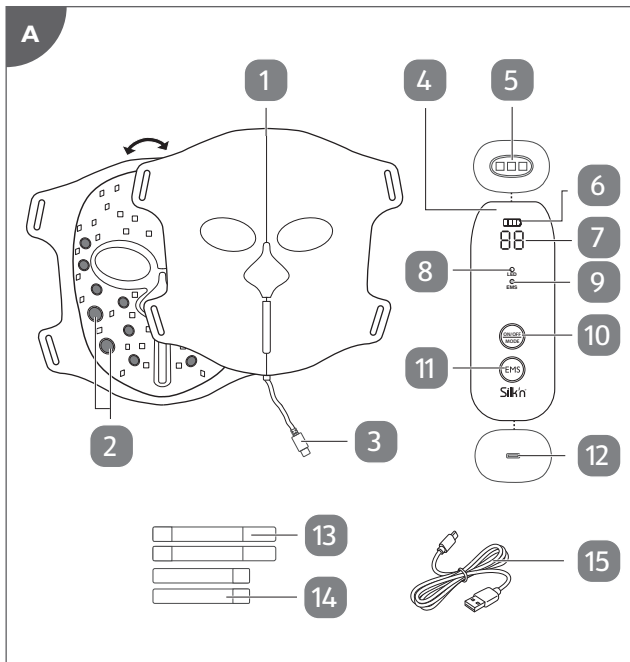Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silk'n Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silk'n Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.

Invention Works B.V.
Galileïstraat 17, 3029 AL Rotterdam, The Netherlands
www.silkn.eu – info@silkn.eu

Table of contents

I	Package contents/device parts _____	5
2	General information _____	5
2.1	Reading and storing the user manual _____	5
2.2	Explanation of symbols _____	5
3	Safety _____	6
3.1	Proper use _____	6
3.2	Safety instructions _____	6
4	Checking the product and package contents _	10
5	Charging _____	10
6	Assembly _____	11
7	Operation _____	11
7.1	Carrying out the session _____	11
7.2	Selecting the operation modes _____	12
7.3	EMS function _____	13
8	Regimen _____	13
9	Cleaning _____	14
10	Storage _____	15
11	Troubleshooting _____	15
12	Technical data _____	15
13	Disposal _____	16
13.1	Disposing of the packaging _____	16
13.2	Disposing of the product _____	16
14	Warranty _____	17
15	Customer service _____	17

1 Package contents/device parts

- | | |
|---------------------------------|----------------------------|
| 1 Face mask | 9 EMS mode indicator light |
| 2 Electrodes | 10 Operation button |
| 3 Plug | 11 EMS button |
| 4 Remote control | 12 Charging port |
| 5 Blue LED lights | 13 Long elastic strap, 2× |
| 6 Battery level indicator light | 14 Short elastic strap, 2× |
| 7 Display | 15 USB charging cable |
| 8 LED mode indicator light | |

2 General information

2.1 Reading and storing the user manual

This user manual accompanies this LED EMS Face Mask (hereafter referred to as the “product”), and contains important information on setup and handling.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product. The user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

WARNING!

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

NOTICE!

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

This symbol identifies voltage and the current rating for direct current.

3 Safety

3.1 Proper use

The product is exclusively designed to recover your skin's glow, elasticity and texture, and to activate the facial muscles through EMS (electrical muscle stimulation). It also helps to reduce the appearance of fine lines, wrinkles and other skin imperfections. It is only intended for private use and not suitable for commercial purposes.

Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer cannot be held liable for damages or injury incurred through improper or incorrect use.

3.2 Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

- Only use the product indoors.
- The product must only be supplied at safety extra-low voltage corresponding to the marking on the product.
- Do not open the housing and never insert any objects into the housing.

- Only parts that comply with the original device data may be used for repairs. This product contains electrical and mechanical parts which are essential for providing protection against sources of danger.
- Only use the product with the supplied USB charging cable.
- Do not use the product or USB charging cable if they are damaged or defective. The manufacturer or customer service must replace the damaged or defective parts to avoid risks.
- Always switch the product off when you are not using it, when you intend to clean it, or in the event of a malfunction. Always unplug the USB charging cable when you are not charging the product.
- Never touch the USB charging cable with wet or damp hands.
- Never move, pull, or carry the product by the USB charging cable.
- Keep the product and USB charging cable away from open flames and hot surfaces.
- Lay the USB charging cable so that it does not pose a tripping hazard.
- Do not kink the USB charging cable or wrap it around the product. Do not lay the USB charging cable over sharp edges.

- Do not immerse the product or USB charging cable in water or other liquids. Keep them away from water.
- Do not use the product while charging.
- Do not try to remove the battery from the product. This product contains batteries that are non-replaceable.

⚠ WARNING! RISK OF INJURY!

- The product cannot be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge unless they have been given supervision or instruction concerning the use of the product in a safe way and understand the hazards involved. Children must not play with the product. Cleaning and user maintenance must not be performed by children.
- The product cannot be used by pregnant women, children, people with epilepsy, thyroid patients, people using thyroid medications, people who are photoallergic or use light allergy medication.
- The product cannot be used on skin diseases, wound ulceration, trauma and post-operative wounds.
- Do not use the product if you have any metal implants in the face/mouth.

- Do not use the Blue LED light or EMS mode for more than 10 minutes per session, and max. 3 sessions per week.
- Consult your physician if you are uncertain whether you can use the product.
- Stop using the product immediately if you feel any discomfort, develop any abnormal reactions (such as itching or uncomfortable inflammation) or in case your skin becomes warm to the point of discomfort.
- Do not look directly at the light source part of the product.
- Keep children away from the product and USB charging cable.

NOTICE! RISK OF DAMAGE!

- Keep the product away from open fire, hot surfaces, high temperatures or the effects of the weather (e.g. rain).
- Do not fold or roll up the product.
- Do not stack heavy objects on top of the product.
- Stop using the product if its parts exhibit cracks or breaks or are deformed. Only replace damaged parts with corresponding original spare parts.

4 Checking the product and package contents

NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.
1. Take the product out of the packaging.
 2. Check to make sure that all parts are included (see **Fig. A**).
 3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

5 Charging

WARNING! RISK OF INJURY!

- Charge only with the supplied USB charging cable.
 - Only charge the product at temperatures between 0 and 40 °C.
-  • Fully charge the product before first use and when the battery power is low. The charging time is around 3.5 hours.
- When the battery power is low, the battery level indicator light **6** flashes 5 times before switching off automatically. Recharge the product.
 - The product cannot be used during charging.
 - If the product is to be stored for an extended period, the charge level should be checked and the battery should be charged regularly.
1. Plug the smaller end of the USB charging cable **15** into the charging port **12** of the remote control **4**.

2. Plug the larger end of the USB charging cable into the USB port of a compatible USB power source (not included) (e.g. USB wall adapter, 5 V DC, 1 A).
3. Plug the USB wall adapter into a mains socket for charging. Check the battery level indicator light **6** for charging progress. The charging is done when all sections of the battery level indicator light illuminate continuously.
4. When charging is done, unplug the USB wall adapter from the mains socket and disconnect the USB charging cable from the charging port of the remote control.

6 Assembly

1. Open the hook-and-loop fasteners on the elastic straps.
2. Guide the long elastic straps **13** through the top openings and the short elastic straps **14** through the bottom openings of the face mask **1** (see step 1 in **Fig. B**).
3. Close the hook-and-loop fasteners to secure the elastic straps on the face mask (see step 2 in **Fig. B**).

7 Operation

WARNING! RISK OF INJURY!

- Do not use the Blue LED light or EMS mode for more than 10 minutes per session, and max. 3 sessions per week.

7.1 Carrying out the session

1. Fully charge the product and properly attach the elastic straps **13**/**14** (see chapter "Charging" and "Assembly").
2. Thoroughly clean and dry your face and hands. Make sure they are free from any make-up or dirt.
3. Sit or lie down and place the face mask **1** over your face.
4. Use the hook-and-loop fasteners on the ends of the elastic straps to fit the face mask comfortably on your face. Make sure that the electrodes **2** are in good contact with the skin.
5. Connect the plug **3** of the face mask into the charging port **12**.
6. Press and hold the operation button **10** to switch on the product.

By default, the product activates the LED function after switching on. The message “F1” shows briefly in the display **7** before switching to the timer countdown (10 minutes).

- If desired, press the operation button repeatedly to select another operation mode.

See chapter “Selecting the operation modes” for more information.

- The product automatically switches off after the countdown. Remove the face mask and apply a serum and/or day/night cream on your face.

You can also press and hold the operation button to switch off the product at any time.

7.2 Selecting the operation modes

- Press the operation button **10** repeatedly to cycle through the operation modes (F1– F4).
- Description of the operation modes:

	Mode	Best for
F1	LED (Red, Dark red, Near-Infrared & Deep Near-Infrared)	Anti-ageing: reduces fine lines and wrinkles, reconstructs collagen and elastine fibres, helps to improve microcirculations for firmer and tighter skin.
F2	LED (Amber, Dark red, Near-Infrared & Deep Near-Infrared)	Reduces pigmentation and redness, calms the skin and promotes the production of collagen for a more even skin tone and tighter skin.
F3	LED (Red & Dark red)	Balances oil secretion, soothes sensitive skin and reduces redness.
F4	LED (Blue)	Accelerates the fading of acne-related redness for a more even skin tone

- The LED mode indicator light **8** illuminates when the LED function is activated.
- The blue LED lights **5** turn on when F4 mode is selected.
- The selected operation mode shows briefly in the display **7** before switching to the timer countdown (10 minutes). The timer is reset whenever you change to another operation mode.

7.3 EMS function

The EMS (electrical muscle stimulation) function generates micro-currents to stimulate muscular movements, improves blood circulation and stimulates the production of collagen and elastin.

To activate EMS function:

1. Press and hold the EMS button **11** to switch on the EMS function. The default intensity is level 1, indicated briefly as "L1" in the display **7**.

The EMS function only works in LED modes 1, 2, and 3.

The EMS mode indicator light **9** illuminates when the EMS function is selected.

2. Press the EMS button repeatedly to cycle through the EMS power levels (1–8).

To avoid any discomfort, a higher power level only activates around 1.5 seconds after selection.

The display shows the selected EMS intensity level (indicated by "L"), the current operation mode (indicated by "F") and the remaining session time alternately.

- You will feel small electrical pulses that contract the facial muscles. Use the highest energy level that still feels comfortable to you. Lower the energy level once you feel any discomfort.
- The session lasts for the duration of your selected LED session. The product switches off automatically after the countdown. You can also press and hold the operation button to switch off the product at any time.
- If you want to stop using the EMS function during a LED session, press and hold EMS button until the EMS function is turned off.

8 Regimen

WARNING! RISK OF INJURY!

- Do not use the Blue LED light or EMS mode for more than 10 minutes per session, and max. 3 sessions per week. All other modes

can be used on a daily basis starting from the fourth week.

If you are using the product for the first time, follow the schedule below:

First week	Use once every 3 days, max. 10 minutes at a time
Second week	Use once every 2 days, max. 10 minutes at a time
Third week	Use once every 2 days, max. 10 minutes at a time
Fourth week	Use once every day, max. 20 minutes at a time

As soon as results show, adjust the frequency to once or twice per week.

9 Cleaning

WARNING! RISK OF ELECTRIC SHOCK!

- Always unplug the product from the power source before cleaning.

NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners (e.g. vinegar or limescale removers), brushes with metal or nylon bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers and the like.
- Do not use hot water for cleaning.
- Use a damp cloth or a water-based wipe to wipe off the face mask . Then wipe it dry with a dry cloth.

10 Storage

All parts must be completely dry before being stored.

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product so that it is not accessible to children.
- For best protection, store the product in its original packaging.

11 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not operate.	<ul style="list-style-type: none"> • Make sure the product is properly charged. • Make sure the product is not connected to a USB power source. The product does not work while charging.
The product does not charge.	<ul style="list-style-type: none"> • Make sure the USB charging cable is securely connected between the charging port and the USB power source.

12 Technical data

Model:	FAC06
Rated voltage:	3.7 V
Battery capacity:	2,000 mAh lithium
Number of LED lights:	60 (face mask)
	3 (remote control)
Total number of LED light bulbs:	243
Materials:	Food grade silicon, ABS, PC
Weight:	352.4 g
Dimensions (L × W × H):	235 × 223 × 108 mm

13 Disposal

13.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

13.2 Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

Old appliances may not be disposed of in the household waste!

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or county.** This ensures that old appliances are recycled

in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

Batteries and rechargeable batteries may not be disposed of with household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances* or not, to a collection point run by the communal authority or to a retailer, so that they can be disposed of in an

environmentally friendly manner.

Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

Note: To remove the rechargeable battery, wrap a towel around the remote control and use a hammer to smash the housing. Make sure your eyes and hands are protected to avoid injury. Please note, the removal of the rechargeable battery is not reversible.

14 Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

15 Customer service

For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: www.silkn.eu.

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Silk'n UK Ltd.