

# Silk'n®

## DUAL LED MASK FAC04

EN  
DE  
FR  
NL  
ES

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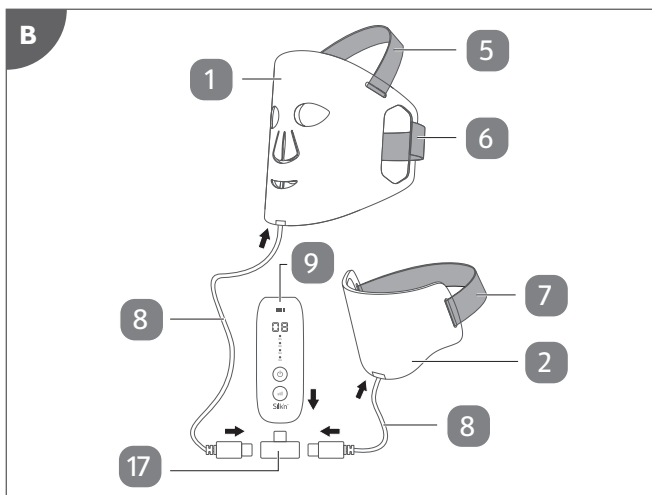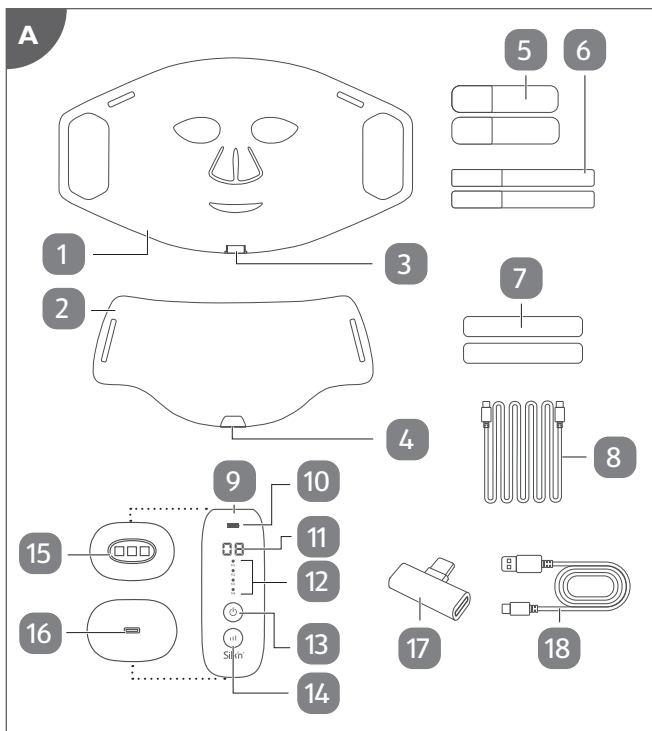## Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silkn Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silkn Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silkn Beauty Ltd.

Invention Works B.V.  
Galileïstraat 17, 3029 AL Rotterdam, The Netherlands  
[www.silkn.eu](http://www.silkn.eu) – [info@silkn.eu](mailto:info@silkn.eu)

# Table of contents

<b>1</b>	<b>Package contents/device parts</b>	<b>5</b>
<b>2</b>	<b>General information</b>	<b>5</b>
2.1	Reading and storing the user manual	5
2.2	Explanation of symbols	6
<b>3</b>	<b>Safety</b>	<b>6</b>
3.1	Proper use	6
3.2	Safety instructions	6
<b>4</b>	<b>Checking the product and package contents</b>	<b>9</b>
<b>5</b>	<b>Charging</b>	<b>10</b>
<b>6</b>	<b>Assembly</b>	<b>10</b>
<b>7</b>	<b>Operation</b>	<b>11</b>
7.1	Preparation	11
7.2	Carrying out the session	12
<b>8</b>	<b>Regimen</b>	<b>12</b>
<b>9</b>	<b>Cleaning</b>	<b>13</b>
<b>10</b>	<b>Storage</b>	<b>13</b>
<b>11</b>	<b>Troubleshooting</b>	<b>13</b>
<b>12</b>	<b>Technical data</b>	<b>14</b>
<b>13</b>	<b>Disposal</b>	<b>14</b>
13.1	Disposing of the packaging	14
13.2	Disposing of the product	14
<b>14</b>	<b>Warranty</b>	<b>15</b>
<b>15</b>	<b>Customer service</b>	<b>15</b>

# I Package contents/device parts

- 1 Face mask
- 2 Neck mask
- 3 Face mask connector
- 4 Neck mask connector
- 5 Small elastic strap, 2×
- 6 Large elastic strap, 2×
- 7 Medium elastic strap, 2×
- 8 Connection cable, 2×
- 9 Controller
- 10 Battery level indicator light
- 11 Numeric display
- 12 Mode indicator light, 4×
- 13 Operation button (⏻)
- 14 Energy level button (⏻)
- 15 Operation lights
- 16 Charging port
- 17 Hub
- 18 USB cable

## 2 General information

### 2.1 Reading and storing the user manual

This user manual accompanies this Dual LED Mask (hereafter referred to as the “product”), and contains important information on setup and handling.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product. The user manual is based on the standards and rules in force in the European Union and United Kingdom. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

## 2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

### **WARNING!**

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

### **NOTICE!**

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

Direct current.

## 3 Safety

### 3.1 Proper use

The product is exclusively designed to recover your skin's glow, elasticity and texture. It is only intended for private use and not suitable for commercial purposes. Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer cannot be held liable for damages or injury incurred through improper or incorrect use.

### 3.2 Safety instructions

#### **WARNING! RISK OF ELECTRIC SHOCK!**

- Only use the product with the supplied USB cable.
- Only connect the USB cable of the product to a suitable mains adaptor (IEC/EN 60335)

certified SELV power supply) connected to an easily accessible socket so that you can quickly disconnect it from the mains supply in the event of a problem.

- Do not use the product if it is damaged or if the USB cable is defective.
- Do not open the housing; instead, have a qualified professional perform repairs. Contact the customer service team for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection or incorrect operation.
- Only parts that comply with the original device data may be used for repairs. This product contains electrical and mechanical parts which are essential for providing protection against sources of danger.
- Do not immerse the product or USB cable in water or other liquids. Keep them away from water.
- Always switch the product off when you are not using it, when you intend to clean it, or in the event of a malfunction. Always unplug the USB cable when you are not charging the product.
- Never touch the USB cable with wet or damp hands.
- Never insert any objects into the housing.
- Keep all parts of the product away from open flames and hot surfaces.
- Lay the USB cable so that it does not pose a tripping hazard.
- Do not pull the USB cable out of the mains adaptor by the cord; instead, always pull it out by the cable plug.
- Do not kink the USB cable and do not lay it over sharp edges.

- Only use the product indoors.
- Do not use the product while charging.

**⚠ WARNING! RISK OF INJURY!**

- The product can be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning the use of the product in a safe way and understand the hazards involved. Children must not play with the product. Cleaning and user maintenance must not be performed by children.
- Keep children away from the product and USB cable.
- The product cannot be used by pregnant women, children, people with epilepsy, thyroid patients, people using thyroid medications, people who are photoallergic or use light allergy medication.
- The product cannot be used on skin diseases, wound ulceration, trauma and post-operative wounds.
- Do not use the product for more than 10 minutes per session.
- Consult your physician if you are uncertain whether you can use the product.
- Stop using the product immediately if you feel any discomfort, develop any abnormal reactions (such as itching or uncomfortable inflammation) or in case your skin becomes warm to the point of discomfort.
- Under abusive conditions, liquid may be ejected from the battery; avoid contact. Liquid ejected from the battery may cause irritation or burns.

- If battery fluid comes into contact with skin, eyes or mucous membrane, wash and rinse the affected areas immediately with clean water. Consult a physician if necessary.
- The battery must not be disassembled, short-circuited or thrown into a fire.

### **NOTICE! RISK OF DAMAGE!**

- Keep the product away from open fire, hot surfaces, high temperatures or the effects of the weather (e.g. rain).
- Stop using the product if its parts exhibit cracks or breaks or are deformed. Only replace damaged parts with corresponding original spare parts.

## **4 Checking the product and package contents**

### **NOTICE! RISK OF DAMAGE!**

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.
1. Take the product out of the packaging.
  2. Check to make sure that all parts are included (see **Fig. A**).
  3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

## 5 Charging

- Fully charge the product before first use and when the battery power is low.
  - When the battery power is low, the battery level indicator lights flash 3 times before the product automatically switches off. Recharge the product.
  - If the product is to be stored for an extended period, the battery level should be checked and the product should be charged regularly.
  - The product cannot be used during charging.
1. Plug the smaller end of the USB cable (18) into the charging port (16) of the controller (9).
  2. Plug the larger end of the USB cable into the USB port of a compatible USB power source (e.g. USB power adaptor).
  3. Plug the USB power adaptor into a mains socket for charging.

Charging starts automatically. Check the battery level indicator lights (10) for the charging progress:

- 1 indicator light = 25 % charged
  - 2 indicator lights = 50 % charged
  - 3 indicator lights = 75 % charged
  - 4 indicator lights = fully charged
4. When charging is done, unplug the USB power adaptor from the mains socket and disconnect the USB cable from the charging port.

## 6 Assembly

Assemble the product as follows (see **Fig. B**):

1. Open the hook-and-loop fasteners on the elastic straps.
2. Guide the hook-and-loop fasteners through the openings on the face mask (1) and neck mask (2):

Small elastic straps (5)	Small top openings of the face mask
Large elastic straps (6)	Large side openings of the face mask
Medium elastic straps (7)	Side openings of the neck mask

3. Close the hook-and-loop fasteners to secure the elastic straps on the masks.

4. Connect one end of a connection cable **8** to the face mask connector **3** and the other end to the hub **17**.
5. Connect one end of the other connection cable to the neck mask connector **4** and the other end to the hub.
6. Attach the hub to the charging port **16** of the controller **9**.

The face mask and neck mask can be used separately. In this case, directly connect the mask to the controller with a connection cable.

## 7 Operation

### WARNING! RISK OF INJURY!

- Be aware that the lights flicker with an increasing speed when turning on the product. For eye safety, first place the masks over your face and neck before turning them on.
- When switching to another mode, the lights flicker again with an increasing speed to indicate the start of the new mode.
- Do not use the product for more than 10 minutes per session.

Clean all parts of the product before first use as described in the chapter "Cleaning".

### 7.1 Preparation

1. Fully charge the product and properly attach the elastic straps (see chapter "Charging" and "Assembly").
2. Thoroughly clean and dry your face and neck. Make sure they are free from any make-up or dirt.
3. Sit or lie down and place the masks over your face and neck.
4. Use the hook-and-loop fasteners on the ends of the elastic straps to fit the face mask **1** and neck mask **2** comfortably on your face and neck.

The small elastic straps **5** should be closed at the top of your head and the large elastic straps **6** should be closed at the back of your head. The medium elastic straps **7** should be closed at the back of your neck.

## 7.2 Carrying out the session

1. Press and hold the operation button (⏻) 13 to switch on the product.

The operation lights 15 illuminate when the product is switched on. By default, the product goes into mode “M1” at energy level 1.

The corresponding mode indicator light 12 illuminates.

The current energy level briefly shows in the numeric display 11 before changing to show the remaining operating time.

2. If necessary, press the operation button (⏻) to select another mode and the energy level button (||) 14 to select another energy level (1–3), respectively.

Mode	LED	Best for
M1	Red + Dark red + Near-Infrared	Reduces fine lines and wrinkles, reconstructs collagen and elastine fibres and helps to improve microcirculations
M2	Dark red + Near-Infrared	For brighter and smoother skin
M3	Red + Dark red	Balances oil secretion and soothes sensitive skin
M4	Blue	Accelerates the fading of acne-related redness for a more even skin tone

3. The session begins and lasts for 10 minutes for modes M1, M2 and M3, and 3 minutes for mode M4. The product automatically switches off after the pre-defined operating time.

You can also press and hold the operation button (⏻) to switch off the product any time.

4. After the session, remove the masks and apply a serum and/or day/night cream on your face and neck.
5. Thoroughly clean the product after use (see chapter “Cleaning”).

## 8 Regimen

If you are using the product for the first time, use the product 3–5 times a week in the first 4 weeks. Afterwards, use the product 1–2 times a week.

## 9 Cleaning

### **!** WARNING! RISK OF ELECTRIC SHOCK!

- Always unplug the product from the power source before cleaning.

### NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners, brushes with metal bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers, and the like.
- Do not use hot water for cleaning.
- Use a damp cloth or a water-based wipe to wipe off the masks **1**/**2**. Then wipe them dry with a dry cloth.

## 10 Storage

All parts must be completely dry before being stored.

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product so that it is not accessible to children.
- For best protection, store the product in its original packaging.

## 11 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not operate.	<ul style="list-style-type: none"> <li>• Make sure the product is properly charged.</li> <li>• Make sure the product is not connected to a USB power source. The product does not work while charging.</li> </ul>
The product does not charge.	<ul style="list-style-type: none"> <li>• Make sure the USB cable <b>18</b> is securely connected between the charging port <b>16</b> and the USB power source.</li> </ul>

## 12 Technical data

Model:	FAC04
Rated input:	5 V === 1 A
Working voltage:	3.7 V
Battery capacity:	2,000 mAh, 7.4 Wh (lithium)
Number of LED lights:	76 (face mask) 36 (neck mask) 3 (controller)
Weight:	200 g (face mask) 80 g (neck mask) 94 g (controller)
Dimensions (L × W × H):	381 × 190.5 × 10 mm (face mask) 237 × 108 × 7.5 mm (neck mask) 133.5 × 50 × 35 mm (controller)
Material:	Silicone (face and neck masks) ABS (controller)
Working temperature:	-10–40 °C

## 13 Disposal

### 13.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

### 13.2 Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

**Old appliances may not be disposed of in the household waste!**

Should the product no longer be capable of being used at some point in time, **dispose of it in**

**accordance with the regulations in force in your city or county.** This ensures that old

appliances are recycled in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

## **Batteries and rechargeable batteries may not be disposed of with household waste!**

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances\* or not, to a collection point run by the communal authority or to a retailer, so that they can be

disposed of in an environmentally friendly manner.

Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

\* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

## **14 Warranty**

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

## **15 Customer service**

For more information about the products of Silk'n visit your regional Silk'n website: [www.silkn.eu](http://www.silkn.eu). If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: [www.silkn.eu](http://www.silkn.eu).

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Copyright 2023 © Silkn Beauty Ltd. Alle rechten voorbehouden.
© 2023 Silkn Beauty Ltd. All rights reserved. Silk'n and the Silk'n logo are registered trademarks of Silkn Beauty Ltd.
Silk'n UK Ltd.