Silk'n®


Beautiful Technology

Pure (Deluxe)

H4101

- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DEL USUARIO
- PT MANUAL DO UTILIZADOR
- GR ΕΓΧΕΙΡΙΔΙΟ ΧΡΗΣΤΗ
- TR KULLANIM KILAVUZU
- IT MANUALE PER L'UTENTE
- NO BRUKSANVISNING
- SV ANVÄNDARHANDBOK
- DA BRUGERMANUAL

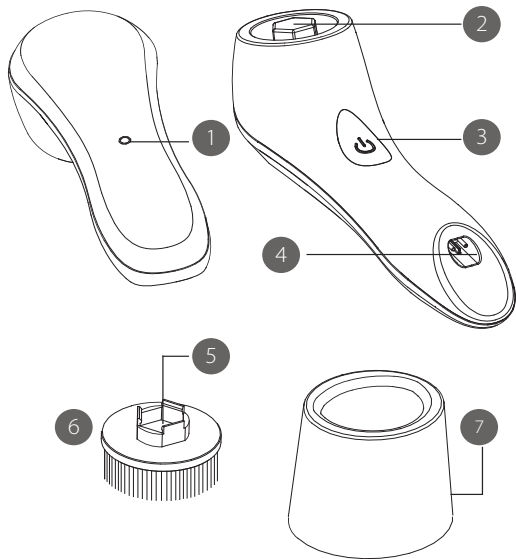Pure Deluxe

Table of contents

1. Warnings & contraindications.....	2
2. Safety with Silk'n Pure (Deluxe).....	7
3. Get to know Silk'n Pure (Deluxe)	7
4. Device controls & settings.....	8
5. Treatment plan.....	9
6. Get started!.....	10
7. Cleaning & maintenance	11
8. Troubleshooting	12
9. Parts & consumables	13
10. Customer service	13
11. Disposal	14
12. Specifications	14
13. Signs & symbols	15
14. Warranty	16
15. Legal notice	16

1. Warnings & contraindications

EN

This user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Please read all warnings and contraindications before use.

- ⚠ Do not use if the parts or accessories are not supplied by the manufacturer as defined in this manual, or if the device is damaged in any way, does not appear to be working properly, or if you see or smell smoke. In such cases, stop using the device and contact Customer Service.
- ⚠ Only authorized personnel are permitted to perform repairs. Do not attempt to repair or open sealed parts of your device. You may expose yourself to dangerous components and cause serious harm to yourself. It will also void your warranty.
- ⚠ The manufacturer accepts no liability for damage or injuries caused by improper or incorrect use.
- ⚠ Keep the adapter away from water. The adapter is not waterproof.

- ⚠ Always ensure that your hands are dry before connecting the adapter to the electricity supply.
- ⚠ Do not place or store the adapter where it can fall or be pulled into a tub, sink, water or any other liquid.
- ⚠ Do not reach for the adapter if it has fallen into water. Unplug immediately before retrieving it. Failure to unplug the adapter may result in an electric shock!
- ⚠ Never use the adapter if it has a damaged cable or plug, it isn't working properly, or has been submerged in water. If that happens, unplug the adapter immediately.
- ⚠ Never pull the adapter out of the socket by the cable. Always remove it by taking hold of the adapter.
- ⚠ Place the adapter cable in such a way that it does not pose a tripping hazard.
- ⚠ Do not kink the adapter cable and do not lay it over sharp edges.
- ⚠ Do not leave the device submerged in bath or shower water for more than 20 minutes.
- ⚠ Do not attach wet brush heads when charging the device.

- ⚠ Do not allow face wash to accumulate on the brush heads.
- ⚠ The device is only to be used with the power supply unit provided.
- ⚠ If a battery has leaked, avoid contact with the skin, eyes and mucous membranes. If necessary, rinse the affected areas with water and seek medical assistance.
- ⚠ Do not clean the device or any of its parts in the dishwasher.
- ⚠ Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. For details, see section 11. Disposal.
- ⚠ This appliance contains batteries that cannot be replaced by users.
- ⚠ This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.
- ⚠ The battery must be removed from the appliance before it is scrapped.

- ⚠ The appliance must be disconnected from the supply mains when removing the battery.
- ⚠ The battery is to be disposed of safely.
- ⚠ Attention should be drawn to the environmental aspects of battery disposal. Don't throw used batteries in a dustbin. Please contact your retailer in order to protect the environment.

Battery Removal

To remove the battery, the device must be broken apart. The device will no longer function after this procedure, nor can it be repaired. Any warranty will be void when the appliance housing is broken and the battery is removed. It is strongly recommended to wear safety protective goggles and gloves.

Silk'n is not liable for any physical or material damages as a result of performing this procedure.

To remove the battery place the device on a hard, strong surface. Use a hammer to hit the appliance on the seam and continue to do so along the whole seam on both sides of the device.

Once the seam is sufficiently broken, use pliers or a screwdriver to take the housing apart. Unscrew the PCB (printed circuit board) and remove the battery by disconnecting the plug from the PCB.

Contraindications

- ⚠ This device can be used by children aged from 8 years and above, and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the device in a safe way and understand the hazards involved.
- ⚠ Children must not play with the device or any of its parts.
- ⚠ Cleaning and maintenance must not be done by children without supervision.
- ⚠ For hygienic reasons, do not share the device with others.
- ⚠ Avoid contact with the eyes.
- ⚠ Do not use a dry brush head on the face or neck, because it may cause skin irritation.

Do not use this device if you have:

- ⚠️ sore or injured skin on the face or neck.
- ⚠️ a pacemaker.

Consult your physician if you have:

- ⚠️ any medical condition or any kind of sensitivity.
- ⚠️ acne or any other skin disorder.
- ⚠️ a rash, redness or itching during or after use.

2. Safety with Silk'n Pure (Deluxe)

Silk'n Pure (Deluxe) has been tested for water ingress to the extent covered by IPX7. This means your device is waterproof and suitable for use in the shower or while taking a bath, but do not leave it submerged in water for more than 20 minutes.

The device also has an overcharge protection. This means that you can safely leave the adapter in the device even if it is connected to the electricity supply, without the risk of damage.

3. Get to know Silk'n Pure (Deluxe)

This user manual refers to both Silk'n Pure and Silk'n Pure Deluxe.

Silk'n Pure (Deluxe) effectively cleanses the skin and removes the outer layer of dry and flaky skin cells. It lifts out impurities from within the pores as well as any build-up of debris from within the folds of fine lines and wrinkles. The device uses vibrations to help smooth the skin's surface, and massages your skin to improve blood circulation.

Silk'n Pure (Deluxe) removes up to 6 times more makeup than manual cleansing, leaving your skin clean, smooth and soft.

Silk'n Pure (Deluxe) is an electrical device and should be used with special attention to safety. **Before use, it is important that you read and understand the instructions on how to use the device including after-use procedures, and follow them strictly.** We recommend you refamiliarise yourself with the instructions in this manual before each treatment.

3.1 Package contents

- Your Silk'n Pure (Deluxe) device
- 2 brush heads (Regular, Soft)
- Face brush stand
- Adapter
- User manual

Silk'n Pure Deluxe includes the following extras:

- Travel case
- Face wash travel bottle
- Face Massage brush head
- Eye Massage brush head

3.2 Device description

Silk'n Pure (Deluxe) is a small, hand-held, electrically powered skin cleanser. Refer to the diagram at the front of this manual to see the location of each feature.

- Charging light (#1)
- Brush attachment point (#2)
- Control button (#3)
- Adapter socket (#4)
- Brush head center (#5)
- Brush head (#6)
- Face brush stand (#7)
- Travel case (#8) *
- Face wash travel bottle (#9) *
- Face Massage brush head (#10) *
- Eye Massage brush head (#11) *

* Silk'n Pure Deluxe

3.3 Intended use

Silk'n Pure (Deluxe) is an over the counter device intended for cleansing your face and neck, and massaging your skin. Do not use for any purpose other than as described in this manual.

4. Device controls & settings

4.1 Charging

Before first use, your device needs to be fully charged. For best results, charge your Silk'n Pure (Deluxe) for 6 hours before using it for the first time, even if there is some residual charge left in the device as a result of quality control procedures.

1. For your safety, make sure the device is completely dry before charging.
2. Connect the **adapter** to the device and the electricity supply.
3. When charging, the **charging light** is red and blinks continuously. When fully charged, the **charging light** is constant red.

 If you are charging for the first time but the **charging light** is constant red, you still need to charge the device for 3 hours.

4. We recommend you unplug the device and store the **adapter** after each full charge until you need to charge the device again.

 We recommend you charge the device fully at least every 3 months to sustain battery life.

 Do not attach wet brush heads when charging the device.

4.2 Indicator lights and function settings

Your Silk'n Pure (Deluxe) is easy to use and has one **control button** and one **indicator light**. Refer to the diagram at the front of the manual to see the location of each feature.


Button	Action	Function / Mode
Control button (#3) 	Press x 1	Vibrates - low speed
	Press x 2	Vibrates - high speed
	Press x 3	Pulsates - low speed
	Press x 4	Pulsates - high speed
	Press x 5	Device switches off
Charging light (#1) 	Light is red and blinks for 10 seconds	Low battery
	Light is red and blinks continuously	Charging
	Light is constant red (visible only when connected to the electricity supply)	Fully charged. The light will go out when the device is disconnected from the adapter , or when you remove the adapter from the electricity supply.

5. Treatment plan

-  For daily use, but not more than twice a day.
-  We recommend cleansing in these areas: forehead, nose, chin, cheeks and neck.
-  The recommended treatment time is 10 seconds per area.
-  Silk'n Pure (Deluxe) is designed to cleanse your skin in just 1 minute. We recommend that you do not to exceed this time.

6. Get started!

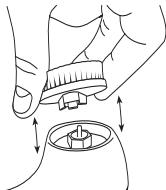1. Select the **brush head** suitable for your skin type. Gently push the **brush head** onto the **brush attachment point**, until you hear a click.
2. Dampen your face with lukewarm water.
3. Moisten the **brush head** with lukewarm water and apply the face wash of your choice onto the brush head.

-  Do not use a dry **brush head** on the face or neck, because it may cause skin irritation.
4. Press the **control button** to turn on the device. The device automatically starts in the low speed vibration mode.
 5. To select the desired mode, keep pressing the control button to cycle through the modes (see section 4.2 Indicator lights and function settings).
 6. Gently press the **brush head** against your skin and move it in an upward circular motion. There's no need to apply any pressure.
-  Try to avoid delicate skin areas such as around the eyes or lips.
7. To finish, rinse your skin with water, and pat dry with a towel.

 The device is set to switch off automatically after 1 minute. However, if you change modes during a treatment, the timer will reset to zero and will start again.

6.1 Silk'n Pure Deluxe accessories

In addition to the handy travel case and face wash travel bottle, Silk'n Pure Deluxe includes two additional treatment accessories. Both are suitable for daily use.

The **Face Massage brush head** improves blood circulation and cleanses the skin mildly, leaving your skin fresher and brighter. Follow the steps in section [6. Get started!](#)

The **Eye Massage brush head** helps reduce puffiness and dark circles around the eyes. To use this accessory, follow these steps:

1. Cleanse your face and make sure the skin is clean and dry.
2. Place the metal part on the area underneath your eyes and then slide from left to right. Do this for 30 seconds to a 1 minute in one area. You can also use it with your favourite eye serum, lotion or cream.

7. Cleaning & maintenance

Use warm water up to a maximum of 45°C / 113°F.

7.1 Cleaning the brush heads

Clean the Regular and/or Soft **brush head(s)** after each use.

- Make sure the device is switched off. Rinse with warm water and mild soap to remove any build-up that may have accumulated. Then shake off any excess water.
- Dry the **brush head** with a clean cloth. You can also dry the **brush head** by switching the device back on, and rubbing it against a clean cloth for 5-10 seconds.
- ❗ For hygienic reasons, we recommend you replace used **brush heads** every 3 months, or earlier if you see that they are worn.

⚠ Do not allow face wash to accumulate on the **brush heads**.

7.2 Cleaning Silk'n Pure Deluxe accessories

- Wash the **Face Massage brush head** with warm water, and remove any build-up that may have accumulated. Then dry it with a clean cloth.
- Wipe the **Eye Massage brush head** with a damp cloth using a mild detergent. Alternatively, wash it with warm water and remove any build-up that may have accumulated. Then dry it with a clean cloth.


7.3 Cleaning the device and the face brush stand

We recommend that you clean the device after each use, and the **face brush stand** once a month. Before cleaning:

- Disconnect the device from the **adapter**.
 - Do not use abrasive cleaners or chemicals to clean any parts.
 - Clean the device or **face brush stand** with a damp cloth using a mild detergent, then wipe it dry with a clean cloth.
- ⚠ Do not clean the device or any of its parts in the dishwasher.

7.4 Replacing the brush head

1. To remove the **brush head**, simply grip the outer plastic base of the **brush head** with one hand, and the device handle with the other, and pull upwards until it snaps loose.
2. To attach a **brush head**, place the center of the **brush head** on the **brush attachment point**.
3. Gently press the **brush head** and push down until it clicks into place.

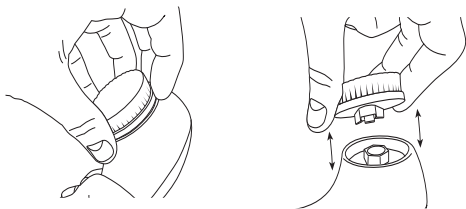7.5 Storage

- You can leave the device in the **face brush stand** and use it as a counter top stand.
- If you don't plan to use the device for some time, we recommend you store it in a dry place.

8. Troubleshooting

Do not attempt to repair or open sealed parts of your device. You may expose yourself to dangerous components and cause serious harm to yourself. It will also void your warranty.

- ⚠ Only authorized personnel are permitted to perform repairs. For problems or queries, please contact Silk'n Customer Service.

Problem	Action
The device does not start.	Make sure the device is properly charged. Attach the adapter to the device and the electricity supply and leave the device to charge for at least 6 hours before the next use.
The device does not charge.	Make sure the adapter is connected correctly to the device and the electricity supply. Unplug the adapter and then reconnect it again to completely reset the device. If the problem persists, contact Customer Service.

Problem	Action
The device shuts off all the time.	Make sure the device is properly charged. Attach the adapter to the device and the electricity supply, and charge the device for at least 6 hours before the next use.

9. Parts & consumables

You can purchase replacement brushes through the retailer where you bought the device, or you can order directly from our website: www.silkn.eu.

10. Customer service

For more information see our website www.silkn.eu. This manual is also available as a PDF download from the website. Contact Silk'n Customer Service to report unexpected operations or events, or any other problem with your device.

Our contact details are:

Country	Service number	Email
Belgium	0800-29316	serviceconsommateurbe@silkn.eu
France	0891-655557	serviceconsommateurfr@silkn.eu
Germany	089 51 23 44 23	kundenservicede@silkn.eu
Italy	+31 (0)180-330 550	servizioclientiit@silkn.eu
Netherlands	0900-2502217	klantenservicenl@silkn.eu
Spain	900 823 302	servicioalcliente@silkn.eu
UK	0906-2130009	customercareuk@silkn.eu
Other countries	+31 (0)180-330 550	info@silkn.eu

11. Disposal

Packaging: Dispose of the packaging separated into single type materials. Dispose of cardboard and carton as waste paper and film via the recyclable material collection service.

Device: Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices. This ensures that the device is recycled professionally, and prevents the release of harmful substances into the environment. Make sure that the device is in an uncharged state before disposal.

12. Specifications

AC/DC Adapter	
Model number:	RSS1002-025050-W2E (EU) / RSS1002-25050-W3U (UK)
Input:	100–240 V~ 50/60 Hz 0.2 A
Output:	5VDC; 500mA
Output power:	2,5 W
Protection class:	II
IP code:	IP20

Published information	Value and precision	Unit
Name or trademark of the manufacturer, commercial register number and address:	Dongguan GaoYi Electronic Co., Ltd, West Sha Yao Road XiSha Road ShiJie Town Dongguan, Guangdong, China	
Model identifier:	RSS1002-025050-W2E (EU) RSS1002-25050-W3U (UK)	
Input voltage:	100–240	V
Input alternating current frequency:	50/60	Hz
Output voltage:	5	V
Output current:	0.5	A
Output power:	2,5	W
Average efficiency during operation:	77.24	%
Power consumption at zero load:	0.04	W

13. Signs & symbols

The following symbols and signal words are used in this user manual, on the device or on the packaging.

EN

	Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.
IPX7	Protected against water immersion for up to 30 minutes at a depth of 1m.
	Do not dispose with household waste! The WEEE directive requires that you dispose of this electrical device (including all accessories and parts) at a designated recycling collection point for electrical devices.
	This symbol identifies electrical devices that belong to the protection class II.
	The word "GS" stands for tested safety. Products marked with this symbol meet the requirements of the German Product Safety Act (ProdSG) (EU adapter).
	Devices marked with this symbol are only suitable for household use (dry indoor rooms).
	The adapter is equipped with a short-circuit proof safety transformer
IP20	The adapter is protected against solid objects up to 12 mm.
	The adapter has a switching power supply.
	The international efficiency mark indicates that the power supply meets the level VI requirements.

	This symbol indicates the polarity of the pins in the adapter.
	The product fulfills national and international safety standards (GB adapter).
	Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the United Kingdom.

14. Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service centre. Our customer service is always pleased to help you!

15. Legal notice

Copyright 2023 © Silk'n Beauty Ltd. All rights reserved.

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V. No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification. Silk'n Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silk'n Beauty Ltd. Specifications are subject to change without notice. Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd., Invention Works B.V. Galileistraat 17, 3029 AL Rotterdam, The Netherlands www.silkn.eu, info@silkn.eu.

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Copyright 2023 © Silkn Beauty Ltd. Alle Rechte vorbehalten.
Copyright 2023 © Silk'n Beauty Ltd. Tous droits réservés.
Copyright 2023 © Silk'n Beauty Ltd. Alle rechten voorbehouden.
Copyright 2023 © Silk'n Beauty Ltd. Todos los derechos reservados.
Copyright 2023 © Silk'n Beauty Ltd. Todos os direitos reservados.
Copyright 2023 © Silk'n Beauty Ltd. Tutti i diritti riservati.