

## FACETITE MINI FAC01

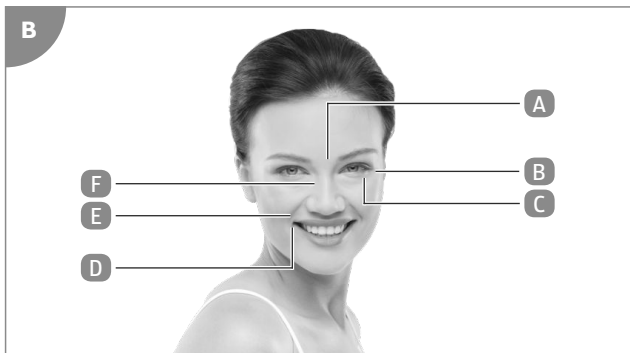- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DE USUARIO
- PT MANUAL DO UTILIZADOR

- GR ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ
- TR KULLANIM KILAVUZU
- IT MANUALE D'USO
- NO BRUKERHÅNDBOK
- SV ANVÄNDARMANUAL
- DA BRUGERMANUAL

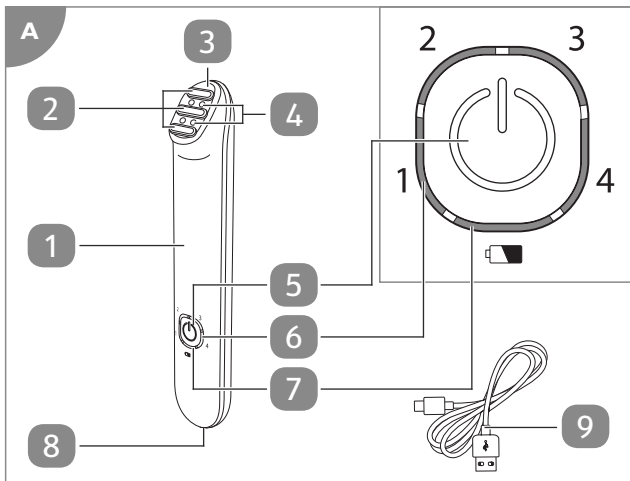Copyright 2023 © Silkn Beauty Ltd. All rights reserved.

## **Legal notice**

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silkn Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silkn Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silkn Beauty Ltd.

Invention Works B.V.  
Galileïstraat 17, 3029 AL Rotterdam, The Netherlands  
[www.silkn.eu](http://www.silkn.eu) – [info@silkn.eu](mailto:info@silkn.eu)

# Table of content

<b>1</b>	<b>Package contents/product parts</b>	<b>5</b>
<b>2</b>	<b>General information</b>	<b>5</b>
2.1	Reading and storing the user manual	5
2.2	Explanation of symbols	5
<b>3</b>	<b>Safety</b>	<b>6</b>
3.1	Proper use	6
3.2	Safety instructions	7
<b>4</b>	<b>Product description</b>	<b>12</b>
<b>5</b>	<b>Check the product and package contents</b>	<b>13</b>
<b>6</b>	<b>Charging</b>	<b>13</b>
<b>7</b>	<b>Operation</b>	<b>14</b>
7.1	Switching the product on/off and selecting an energy level	14
7.2	Treatments	15
7.3	What to expect	15
7.4	Possible side effects	15
7.5	Get started!	16
<b>8</b>	<b>Cleaning</b>	<b>17</b>
<b>9</b>	<b>Storage</b>	<b>18</b>
<b>10</b>	<b>Troubleshooting</b>	<b>18</b>
<b>11</b>	<b>Technical data</b>	<b>20</b>
<b>12</b>	<b>Declaration of conformity</b>	<b>20</b>
<b>13</b>	<b>Disposal</b>	<b>21</b>
13.1	Disposing of packaging	21
13.2	Dispose of old electrical product	21
<b>14</b>	<b>Warranty</b>	<b>21</b>
<b>15</b>	<b>Customer service</b>	<b>22</b>

# I Package contents/product parts

- |   |                   |   |                              |
|---|-------------------|---|------------------------------|
| 1 | FaceTite Mini     | 6 | Energy level indicator light |
| 2 | Electrodes        | 7 | Battery indicator light      |
| 3 | Treatment surface | 8 | Charging port                |
| 4 | LED lights        | 9 | USB cable                    |
| 5 | Control button    |   |                              |

## 2 General information

### 2.1 Reading and storing the user manual

This user manual is a part of the FaceTite Mini (hereinafter called "product"). It contains important information on how to set up and use the product.

This user manual is based on the standards and regulations that are valid within the European Union. Outside the EU, please also note the country-specific directives and laws.

Read the user manual carefully, in particular the safety instructions, before using the product. Failure to follow this user manual may lead to severe injuries or product damage.

Keep this user manual for future reference. If you pass this product on to a third party, you must also supply this user manual.

### 2.2 Explanation of symbols

The following symbols are used in this user manual, on the product or on the packaging.

## WARNING!

This signal symbol/word denotes a hazard with an average risk level that could lead to death or severe injury if it is not avoided.

## CAUTION!

This signal symbol/word denotes a hazard with a low risk level that could lead to mild or moderate injury if it is not avoided.

## NOTICE!

This signal word provides a warning about potential material damage.

This symbol provides you with useful additional information on handling and use.

Declaration of conformity (see chapter “Declaration of conformity”): Products marked with this symbol meet all the applicable Community regulations of the European Economic Area.

Do not use the product in the vicinity of bathtubs, showers, washbasins or other vessels containing water.

## 3 Safety

### 3.1 Proper use

The product is exclusively intended for anti-ageing treatments of the facial skin.

The product is intended for private use only and is not suitable for commercial use. Use the product only as described in this

operating manual. Any other use is improper and may lead to product damage or even personal injury. This product is not a toy. The manufacturer or retailer assumes no liability for damage caused by improper or incorrect use.

### 3.2 Safety instructions

#### **WARNING! RISK OF ELECTRIC SHOCK!**

- Only connect the USB cable if the voltage of the voltage source matches the information shown on the rating plate.
- The product must only be operated with extra low voltage.
- Do not use the product or USB cable if they are damaged or defective. The manufacturer or customer service must replace the damaged or defective parts to avoid risks.
- Only use the product with the supplied USB cable.
- The supplied USB cable can only be connected to the product. Do not use for other purposes.
- Do not open the housing; instead, have a qualified professional perform repairs. Contact the customer service team for this. Liability and warranty claims are waived in the event of repairs performed

by the user, improper connection or incorrect operation.

- Do not use the product while charging.
- Do not immerse the product or USB cable in water or other liquids.
- Keep the USB cable away from water. It is not waterproof.
- Always switch the product off when you are not using it, when you intend to clean it, or in the event of a malfunction. Always unplug the USB cable when you are not charging the product.
- Never touch the USB cable with wet or damp hands.
- Never insert any objects into the housing.
- Never move, pull, or carry the product by the USB cable.
- Keep the product and USB cable away from open flames and hot surfaces.
- Only use the product indoors.

 **WARNING! DANGER TO CHILDREN AND PERSONS WITH REDUCED PHYSICAL, SENSORY OR MENTAL ABILITIES (E.G. PARTIALLY DISABLED, ELDERLY PEOPLE WITH REDUCED PHYSICAL AND MENTAL ABILITIES) OR LACK OF EXPERIENCE AND KNOWLEDGE (E.G. OLDER CHILDREN).**

- Children and persons with reduced physical, sensory or mental abilities (e.g. partially disabled, elderly persons with reduced physical and mental abilities) or lack of experience and knowledge (e.g. older children) are not allowed to operate the product.
- Keep children away from the product and the power cable.
- Children may not play with the product.
- Cleaning and user maintenance must not be carried out by children.

 **CAUTION!** **Do not use if you:**

- have a pacemaker or internal defibrillator, or any electrical implant anywhere in the body.
- suffer from poorly controlled endocrine disorders, such as diabetes and thyroid conditions.
- have a history of skin cancer, or current condition of cancer, or pre-malignant moles.
- have received radiation therapy or chemotherapy treatments within the past 3 months, or have areas of potential skin malignancies.
- have an impaired immune system due to immunosuppressive diseases such as AIDS and HIV, or use immunosuppressive medications.
- have used Isotretinoin (Accutane®) within 6 months prior to treatment.
- have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.
- suffer from epilepsy.
- are pregnant or nursing.

 **Do not use on areas that have:**

- a permanent implant such as metal plates and screws, silicone implants, or areas with injected fills or toxins.
- any active condition in the treatment area, such as sores, psoriasis, eczema, rashes, open lacerations, abrasions, or lesions.
- traumatized or injured skin.
- had recent exposure to the sun without protection (wait 1 day); for tanned or sunburnt skin, wait 3 days.
- had any surgery within 3 months prior to treatment.
- tattoos.
- long or thick hair. Use of the product may result in light bruising.

**i** **Consult your physician before continue to use, if you have side effects, such as itchy skin within 48 hours of first use.**

**i** **Consult your physician if you have:**

- couperose or rosacea.
- a history of diseases stimulated by heat, such as recurrent herpes simplex in the treatment area. If so, you may treat only following a prophylactic regime.
- used physician-prescribed drugs on your skin within the past 6 months.
- a history of skin disorders, keloids, abnormal wound healing, as well as very dry and fragile skin.
- any other condition which, in the opinion of your physician or dermatologist, would make it unsafe for you to be treated.

## **NOTICE! RISK OF DAMAGE!**

- Operate the product on an easily accessible, level, dry, heat resistant and sufficiently stable work surface. Never place the product on the edge of the work surface.
- Never place the product on or near hot surfaces (stove tops, etc.) and never expose the product to high temperatures (heaters, etc.) or to bad weather (rain, etc.).
- Keep the product and the power cable away from open flames and hot surfaces.

- Never immerse the product in water to clean it.
- Never wash the product in the dishwasher. This would destroy it.
- Stop using the product if any of its components are cracked or deformed. Only replace damaged components with matching original spare parts.
- Do not attempt to repair or modify the product yourself. Have repairs carried out by trained personnel.
- Only use the product at ambient temperatures between 10 °C and 35 °C.
- This appliance contains batteries that cannot be replaced by users.
- This appliance contains batteries that are non-replaceable. When the battery is at end of life, the appliance shall be properly disposed of.

## 4 Product description

This product has been designed for the reduction of facial signs of ageing. It makes use of RF technology, which has an effect on the delicate regions of the face. It deeply penetrates the skin to reconstruct the collagen, tightens the tissue and reduces fine lines. You can treat the following parts of the face (see **Fig. B**):

- |                                    |                                    |
|------------------------------------|------------------------------------|
| <b>A</b> lines between eyebrows    | <b>D</b> wrinkles around the mouth |
| <b>B</b> crow's feet               | <b>E</b> nasolabial folds          |
| <b>C</b> fine lines under the eyes | <b>F</b> tear trough               |

## 5 Check the product and package contents

### NOTICE! RISK OF DAMAGE!

 If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.

1. Take the product out of the packaging.
2. Check to make sure that all parts are included (see **Fig. A**)
3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

## 6 Charging

### WARNING! RISK OF INJURY!

- Do not charge when the product is in use or if it is wet. Make sure it is dry before charging.
- Only charge the product with the supplied USB cable.
- Only charge the product at temperatures between 0 and 40 °C.

- Always unplug the USB cable after charging.

- Charge the battery completely before first use.
- The battery indicator light **7** will blink orange when the battery is low. Recharge the product.
- The product cannot be used during charging.
- When fully charged, the product can be used for 20-30 minutes. The duration of the product use may vary according to different energy levels.

1. Plug the smaller end of the USB cable **9** into the charging port **8** of the product.
2. Plug the larger end of the USB cable into the USB port of a compatible USB power source (e.g. USB power adapter or computer) for charging.  
Charging starts automatically and the battery indicator light will blink green.
3. When the product is fully charged the battery indicator light is constant green. Unplug the USB cable from the charging port and the USB power source.

## 7 Operation

### 7.1 Switching the product on/off and selecting an energy level

- The product has four energy levels. Level 1 is the lowest energy level and level 4 is the highest.
- With the control button **5**, you can turn the product on and off, and you can select an energy level.

1. Press and hold the control button for 2 seconds.  
Level 1 is activated. The energy level indicator light **6** lights up white.  
Press the control button again to move to the next level.  
After level 4, it goes back to level 1.

2. Press and hold the control button again for 2 seconds.  
The product is switched off.

## 7.2 Treatments

- Start with energy level 1, and change to a higher energy level as long as the temperature feels comfortable in the treatment area.
- Use the product for 3 to 5 minutes on each treatment area for maximum effectiveness.
- Use the product 2 to 3 times a week for the first eight weeks.
- After eight weeks, use the product weekly, to maintain the effect.

### For maximum effectiveness:

- Make sure you apply enough Slider Gel (or any other water-based gel) to cover the electrodes on the treatment area before use.
- Do not leave the product on one spot of the skin. Keep moving the product while covering the entire treatment area. This will prevent overheating of the skin and reduces the possibility of any adverse effect.
- Make sure the 3 electrodes **2** are always in full contact with the skin.
- If the energy level feels uncomfortable, reduce it to a level that is comfortable for you.
- Only change the energy level when the product is not in contact with skin.

## 7.3 What to expect

During a treatment session it is normal to experience or feel:

- A sensation of warmth
- Some redness or slight swelling: this should subside within a few hours, or up to one day.
- A buzzing sound of the RF current.

## 7.4 Possible side effects

If you follow the instructions fully and pay attention to all warnings and contraindications, complications and side effects are uncommon. However, every cosmetic treatment, including those

designed for home use, involves some degree of risk. Therefore it is important that you understand and accept the risks and complications that may occur:

- Excessive edema (swelling) and erythema (redness) are rare reactions to the treatment that, in most cases, subside within a few days. In case these signs appear, immediately use a damp cool cloth to cool the area for at least 10 minutes and apply an anti-burn cream.
- In rare cases the user may experience dilation of the superficial blood capillaries. This may result in a temporary redness for a few days. Also, some of the capillaries may present red spots that will also reduce with time. If these signs appear, we recommend that you cool the area with chilled water immediately.

**i** Consult your physician if you experience any persistent symptoms including pain, skin redness, swelling or a burning sensation.

**i** If unexpected side effects, such as itching, occur within 48 hours of your first treatment, seek medical advice before continuing treatments.

## 7.5 Get started!

1. Make sure that your skin is clean and dry, and that you have removed all jewelry.
2. Remove the product and the other components from the box.
3. Remove the protective film that covers the electrodes **2**.
4. Start the product: Press the control button **5**. The product starts at energy level I automatically, and I energy level indicator light will turn on.

**i** Check whether the Slider Gel has expired.

5. Apply a thin layer of Slider Gel over the treatment surface of the product.

**i** Slider Gel (or any other water-based gel) facilitates smooth movement of the treatment surface over the skin and enables conductivity. To ensure safe and efficient treatment, always start by applying Slider Gel and apply more gel if needed.

- Using energy level 1, place the product on the area you want to treat. Make sure that it is in good contact with the skin. You will feel a mild sensation of warmth.
- Move the product in gentle massage-like motions on the area you wish to treat. Progress slowly and carefully. Keep moving the product ensuring that the 3 electrodes stay in full contact with the skin as you treat the area.
- With the product constantly in motion, treat **each** area for at least 3-5 minutes.
- Increase the energy level gradually. Using the highest energy level will get the best results, but it's important that you set the energy level that is comfortable for the area you are treating.

**i** You will feel the skin temperature increase as you increase the energy level. It usually takes about 1 minute for the product to reach the required temperature. During the treatment your skin will get a little red. Don't worry; this is normal. Stop immediately if you feel any discomfort.

- Try to move the product slowly at first. If the treatment area is hot and uncomfortable, try to move the product faster.
- Clean the treated area: remove the Slider Gel residue off the treatment surface, and clean the skin with warm water. Dry the skin and apply moisturizing facial products.

## 8 Cleaning

### NOTICE! RISK OF SHORT CIRCUIT!

- Do not immerse the product in water or other liquids.
- Make sure that no water or other liquid gets inside the housing.

## NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners, sharp or metal cleaning utensils such as knives, hard scrapers or similar items. These could damage the surfaces.
1. Turn off the product. Gently wipe off the Slider Gel of the treatment surface **3** with clean, damp paper (such as toilet paper or a wipe). Make sure you remove all gel residue to prevent it from solidifying or damaging the product. You can also use cleaning alcohol (up to 65%).
  2. Wipe the product and the treatment surface with a slightly damp cloth.
  3. Dry off the product with a soft, dry cloth.

## 9 Storage

After cleaning, we recommend you keep the product in the original packaging. Make sure that the product is completely dry before storage.

## 10 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Malfunction	Possible cause	Solution
The product cannot be switched on.	The battery is flat.	Charge the product (see section "Charging").
The battery indicator <b>7</b> blinks orange.	The battery is low.	
Product does not get warm although it is switched on.	Insufficient contact between the treatment surface <b>3</b> and the skin or not enough Slider Gel is applied to the treatment surface.	Apply enough Slider Gel to the treatment surface and make sure it is in full contact with the skin.

Malfunction	Possible cause	Solution
One or several lights 6 blink continuously.	A malfunction has occurred.	Switch the product off completely and then switch it on again.

## 11 Technical data

Model:	FAC01
Technology:	RF & LED light energy technology
LED light energy:	630 nm $\pm$ 5 nm
RF:	1 MHz $\pm$ 30%; max. 13 W
Operation and safety:	<p><b>Thermal sensor and heat stabilizer:</b> These measure the skin temperature, stabilize the heat energy and prevent overheating. When the skin temperature reaches 43°C, the stabilizer will automatically stop emitting heat, until the skin temperature lowers enough for the stabilizer to resume heating.</p> <p><b>Movement sensor:</b> The sensor will only emit RF heat when you are moving the product across your skin.</p>
Product weight:	approx. 51 g
Operating temperature:	10 °C - 35 °C
Storage temperature:	-20 °C - +60 °C
Relative air humidity during operation:	30 - 75%
Relative air humidity during storage:	10 - 75%
Ambient pressure during operation:	500 - 1060 hPa
Ambient pressure during storage:	700 - 1060 hPa
Input voltage:	5 V/1 A $\equiv$
Battery capacity:	500 mAh
Charging connection:	USB type C
Charging time:	approx. 3 hours
Operating time:	20 to 30 minutes with fully charged battery

## 12 Declaration of conformity

**CE** You can request the EU declaration of conformity from the address provided at the end of this manual.

## 13 Disposal

### 13.1 Disposing of packaging

Dispose of packaging according to type. Sort the paperboard and cardboard as waste paper and the film as recyclable material.

### 13.2 Dispose of old electrical product

Applicable in the European Union and other European countries with systems for separate collections of recyclable waste.

**Old electrical products may not be disposed of as household waste!**

If the product cannot be used anymore, each consumer **is obligated by law to dispose of old products**

**separately from household waste**, e.g. at a municipal/city district collection point. It is thereby ensured that old products are professionally recycled and negative effects on the environment are avoided. Therefore, electrical products are marked with the above symbol.

**Batteries and rechargeable batteries must not go in domestic waste!**

As a consumer you are legally obliged to dispose of all batteries, irrespective of whether they contain harmful substances\* or not, at a collection point in your community/district or at a retail outlet, so that they can be disposed of in an environmentally friendly manner.

\* marked with: Cd = Cadmium, Hg = Mercury, Pb = Lead, Li = Lithium.

## 14 Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

## 15 Customer service

For more information about the products of Silk'n visit your regional Silk'n website: [www.silkn.eu](http://www.silkn.eu). If the product is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: [www.silkn.eu](http://www.silkn.eu).

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Copyright 2023 © Silk'n Beauty Ltd. Tous droits réservés.
Copyright 2023 © Silkn Beauty Ltd. Alle rechten voorbehouden.
Copyright 2023 © Silk'n Beauty Ltd. Todos os direitos reservados.
Copyright 2023 © Silk'n Beauty Ltd. Με την επιφύλαξη παντός δικαιώματος.
Copyright 2023 © Silkn Beauty Ltd. Tutti i diritti riservati.
Copyright 2023 © Silkn Beauty Ltd. Alle rettigheder forbeholdes.