Silk'n®

Beautiful Technology

BRIGHT


ND-FC04

- EN USER MANUAL
- DE BEDIENUNGSANLEITUNG
- FR MODE D'EMPLOI
- NL GEBRUIKSAANWIJZING
- ES MANUAL DE USUARIO
- PT MANUAL DO UTILIZADOR

- GR ΟΔΗΓΙΕΣ ΧΡΗΣΗΣ
- TR KULLANIM KILAVUZU
- IT MANUALE D'USO
- NO BRUKERHÅNDBOK
- SV ANVÄNDARMANUAL
- DA BRUGERMANUAL

EN
DE
FR
NL
ES
PT
GR
TR
IT
NO
SV
DA

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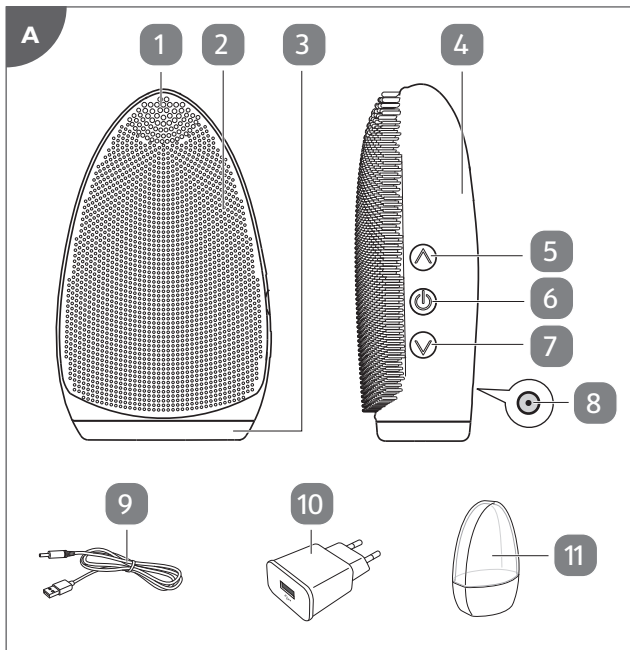Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silk'n Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silk'n Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.

Invention Works B.V.
Galileïstraat 17, 3029 AL Rotterdam, The Netherlands
www.silkn.eu – info@silkn.eu

Table of contents

1	Package contents/device parts _____	5
2	General information _____	5
2.1	Reading and storing the user manual _____	5
2.2	Explanation of symbols _____	5
3	Safety _____	7
3.1	Proper use _____	7
3.2	Safety instructions _____	7
4	Checking the facial brush and package contents _____	12
5	Charging _____	12
6	Operation _____	13
6.1	Using different areas of the facial brush _____	13
6.2	Operation modes and speed settings _____	14
6.3	Switching on/off and selecting an operation mode _____	14
6.4	Selecting a speed setting _____	14
6.5	Starting and ending the treatment _____	14
7	Treatment steps _____	15
8	Cleaning _____	16
9	Storage _____	16
10	Troubleshooting _____	17
11	Technical data _____	17
12	Disposal _____	19
12.1	Disposing of the packaging _____	19
12.2	Disposing of the facial brush _____	19
13	Warranty _____	20
14	Customer service _____	20

I Package contents/device parts

- | | |
|----------------------------------|-----------------------|
| 1 Cleansing area (thicker studs) | 7 Decrease button (▼) |
| 2 Cleansing area (thinner studs) | 8 Charging port |
| 3 LED light | 9 USB cable |
| 4 Massaging area | 10 USB wall adapter |
| 5 Increase button (▲) | 11 Travel case |
| 6 Control button (⏻) | |

2 General information

2.1 Reading and storing the user manual

This user manual accompanies the silicone facial cleansing brush – “Bright” (hereafter referred to as the “facial brush”), and contains important information on setup and handling.

Before using the facial brush, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the facial brush. The user manual is based on the standards and rules in force in the European Union. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the facial brush on to third parties.

2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the facial brush, USB wall adapter or packaging.

WARNING!

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

NOTICE!

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

This symbol identifies electrical devices that belong to the protection class III (device).

This symbol identifies electrical devices that belong to the protection class II (USB wall adapter).

The word "GS" stands for tested safety. Products marked with this symbol meet the requirements of the German Product Safety Act (ProdSG).

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

This symbol identifies voltage and the current rating for direct current.

The international efficiency mark indicates that the power supply meets the level VI requirements.

IPX7

The facial brush is protected against water immersion for up to 30 minutes at a depth of 1 m.

The USB wall adapter is equipped with a short-circuit proof safety transformer.

The USB wall adapter has a switching power supply.

This symbol indicates the polarity of the pins in the USB wall adapter.

IP20

The USB wall adapter is protected against solid objects up to 12 mm.

ta

Ambient temperature (USB wall adapter)

3 Safety

3.1 Proper use

The facial brush is exclusively designed for facial cleansing and massaging. It is only intended for private use and not suitable for commercial purposes. The facial brush has been tested for water ingress to the extent covered by IPX7. This means the device is waterproof and suitable for use in the shower or while taking a bath, but do not leave it submerged in water (max. 1 m) or under running water for more than 30 minutes.

Only use the facial brush as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer or vendor cannot be held liable for damages or injury incurred through improper or incorrect use.

3.2 Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

- Only connect the facial brush if the line voltage of the socket corresponds to the data on the rating plate of the USB wall adapter.
- Only connect the USB wall adapter to an easily accessible socket so that you can quickly disconnect it from the mains supply in the event of a problem.
- Only use the facial brush with the supplied USB cable and USB wall adapter.
- The supplied USB cable and USB wall adapter can only be connected to the facial brush. Do not use for other purposes.
- Do not use the facial brush, USB cable or USB wall adapter if they are damaged or defective. The manufacturer or customer

service must replace the damaged or defective parts to avoid risks.

- Do not open the housing; instead, have a qualified professional perform repairs. Contact the customer service team for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection or incorrect operation.
- Only parts that comply with the original device data may be used for repairs. This facial brush contains electrical and mechanical parts which are essential for providing protection against sources of danger.
- Do not use the facial brush while charging.
- Do not immerse the USB cable or USB wall adapter in water or other liquids.
- Keep the USB cable and USB wall adapter away from water. They are not waterproof.
- Do not leave the facial brush submerged in water (max. 1 m) or under running water for more than 30 minutes.
- Always switch the facial brush off when you are not using it, when you intend to clean it, or in the event of a malfunction. Always unplug the USB cable and USB wall adapter when you are not charging the facial brush.

- Never touch the USB cable and USB wall adapter with wet or damp hands.
- Never insert any objects into the housing.
- Do not pull the USB wall adapter out of the socket by the USB cable; instead, always pull it out by the adapter itself. Never move, pull, or carry the facial brush by the USB cable.
- Keep the facial brush, USB cable and USB wall adapter away from open flames and hot surfaces.
- Lay the USB cable so that it does not pose a tripping hazard.
- Do not kink the USB cable or wrap it around the facial brush. Do not lay the USB cable over sharp edges.
- Only use the facial brush indoors.

WARNING! DANGER FOR CHILDREN AND PERSONS WITH IMPAIRED PHYSICAL, SENSORY OR MENTAL CAPACITIES, OR LACK OF EXPERIENCE AND KNOWLEDGE.

- The facial brush can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the facial brush in a safe way

and understand the hazards involved. Children must not play with the facial brush. Cleaning and user maintenance must not be performed by unsupervised children.

- Keep children under the age of 8 away from the facial brush, USB cable, USB wall adapter and travel case.
- Make sure that children do not play with the plastic wrapping. They may get caught in it when playing and suffocate.

 **WARNING! RISK OF INJURY!**

- Do not use the facial brush on the eye lid, inside the eye orbit, around areas that contain inflammatory acne, or areas that have rashes or wounds.
- Stop using the facial brush immediately if you feel any discomfort or develop any abnormal reactions such as itching, uncomfortable inflammation, or scratches.
- Do not focus the facial brush on one spot on your skin and do not use the facial brush more than 2 minutes at a time.
- Never leave the facial brush unattended while it is in operation.
- Consult your doctor before using the facial brush if:
 - you suffer from allergies, skin irritation or skin diseases.

- you are pregnant.
- you are using a pacemaker.
- you are not sure whether the facial brush is suitable for you.

NOTICE! RISK OF DAMAGE!

- Do not use clay-based, silicone-based or grainy cleansers, exfoliators or scrubs with the facial brush as they may damage the soft silicone touch-points of the facial brush.
- Avoid leaving the facial brush in direct sunlight and never expose it to extreme heat or boiling water.
- Do not place the facial brush on the edge of a surface to avoid dropping.
- Never fill the housing of the facial brush with liquid.
- Clean the facial brush after each use to avoid the accumulation of grease and residues of face wash products.
- Stop using the facial brush if its plastic parts exhibit cracks or breaks or are deformed. Only replace damaged parts with corresponding original spare parts.
- Do not operate the facial brush if it has been dropped or is damaged in any way. Have the facial brush checked and/or repaired by qualified technicians if necessary.

4 Checking the facial brush and package contents

NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the facial brush. Be very careful when opening it.
1. Take the facial brush out of the packaging.
 2. Check to make sure that all parts are included (see **Fig. A**).
 3. Check whether the facial brush or individual parts are damaged. If this is the case, do not use the facial brush. Contact our customer service.

5 Charging

WARNING! RISK OF INJURY!

- Do not charge when the facial brush is in use or if it is wet. Make sure it is dry before charging.
 - Only charge the facial brush with the supplied USB cable and USB wall adapter.
 - Only charge the facial brush at temperatures between 0 and 40 °C.
 - Always unplug the facial brush after charging has been completed.
- 
- Fully charge the facial brush before first use and when the battery power is low.
 - When the battery power is low, the LED light becomes red and starts to blink. Recharge the facial brush.
 - The facial brush cannot be used during charging.

1. Plug the smaller end of the USB cable **9** into the charging port **8**.
2. Plug the larger end of the USB cable into the USB port of the USB wall adapter **10**.
You may also plug the USB cable to a compatible USB power source (e.g. USB port of a computer) for charging.
3. Plug the USB wall adapter into a mains socket for charging. During charging, the LED light **3** will become red and start to blink.
4. The LED light will become solid white when it is fully charged. Unplug the USB wall adapter from the mains socket and disconnect the USB cable from the facial brush.

6 Operation

- Clean the facial brush before first use as described in the chapter “Cleaning”.
- Before cleansing an area of the face for the first time, perform a skin sensitivity test on a spot. Start your treatment only if there is no adverse reaction.
- The LED light **3** is solid white during operation and flashes once whenever a button is pressed.

6.1 Using different areas of the facial brush

Area	Best for
Cleansing area (thicker studs) 1	Precise cleansing and removal of dead skin cells and horny layers.
Cleansing area (thinner studs) 2	Gentle cleansing of large areas such as cheeks.
Massaging area 4	Massaging the skin for improved blood circulation and promoting the absorption of skin care products.

6.2 Operation modes and speed settings

The facial brush is equipped with a vibrating and a pulsating mode.

- The vibrating mode uses uniform vibrations for deep, gentle cleansing.
- The pulsating mode uses alternating strong and weak vibrations for a massage-like soothing effect.
- The “Bright” version has 7 speed settings for each mode.

6.3 Switching on/off and selecting an operation mode

1. Press the control button (⏻) **6** once. The facial brush is now switched on and in the vibrating mode.
You will hear one beep sound when the facial brush is switched on.
2. Press the control button (⏻) again. The facial brush is now in the pulsating mode.
3. Press the control button (⏻) a third time to switch off the facial brush.
You will hear two beep sounds when the facial brush is switched off.

6.4 Selecting a speed setting

1. Select an operation mode (see chapter “Switching on/off and selecting an operation mode”).
2. Press the increase button (⬆) **5** or decrease button (⬇) **7** in the selected operation mode to select a speed setting (1–7).
The default level is 4.

6.5 Starting and ending the treatment

- You may start your treatment after selecting an operation mode and a speed setting.
- Once the facial brush starts, you will hear a beep sound every 30 seconds, hinting you about the time that has passed. The facial brush will automatically switch off after 2 minutes. Please note that whenever you change the operation mode or speed setting during treatment, the timer will be reset to zero.
- You can also end the treatment anytime by pressing the control button (⏻) **6**. You may have to press the button more than once, depending on the product version and current operation mode.

7 Treatment steps

1. Make sure the facial brush has enough battery power before use (see chapter “Charging”).
2. Dampen your face and the facial brush with lukewarm water. Apply your regular face wash product on the face (see step 1 in **Fig. B**).
3. Turn on the facial brush, select the vibrating mode and a speed setting (see chapter “Operation”).
4. Cleanse the face by moving the cleansing area (1 or 2) over the face in circular movements. Use the area with thicker studs for areas that need precise cleansing (e.g. around the nose, eyes and chin) and use the area with thinner studs for the larger areas (e.g. forehead, cheeks and neck) (see step 2 in **Fig. B**). The facial brush will automatically switch off after 2 minutes.
5. Thoroughly rinse the face and facial brush with water and then wipe them dry with a clean, dry towel (see step 3 in **Fig. B**).
6. Apply your serum or day/night cream on the face.
7. Turn on the facial brush and select the pulsating mode and a speed setting (see chapter “Operation”).
8. Massage the face by moving the massaging area 4 over the face in circular movements to stimulate blood circulation and promote the absorption of your skin care product (see step 4 in **Fig. B**). The facial brush will automatically switch off after 2 minutes. If you want to end the treatment earlier, press the control button (5) 6. You may have to press the button more than once, depending on the product version and current operation mode.
9. Thoroughly clean the facial brush after use (see chapter “Cleaning”).

- You can divide the skin into 6 areas: forehead, nose, chin, left cheek, right cheek and neck. We recommend you not to treat each area for more than 20 seconds.
- This facial brush is designed to cleanse your skin in just 2 minutes. We recommend you not to exceed this time.
- Do not use this facial brush more than twice a day as this can irritate the skin and cause redness.

8 Cleaning

WARNING! RISK OF INJURY!

- Always unplug the facial brush before cleaning.

NOTICE! RISK OF DAMAGE!

- Do not leave the facial brush submerged in water (max. 1 m) or under running water for more than 30 minutes. Make sure that no water or other liquids penetrate the housing.
- Do not use any aggressive cleaners, brushes with metal or nylon bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers and the like.

Always clean the facial brush after each use.

1. Wash the facial brush with lukewarm water and soap.
2. Thoroughly rinse it with water.
3. Wipe it dry with a dry cloth.

To clean the travel case , use a damp cloth to wipe it off. Then let all parts dry completely.

9 Storage

All parts must be completely dry before being stored.

- Store all parts in a clean and dry place, avoiding heat and direct sunlight.
- Store the facial brush in the travel case  when traveling to avoid damage.

10 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the facial brush yourself.

Fault	Possible cause and solution
The facial brush does not start.	Make sure the facial brush is properly charged.
The facial brush does not charge.	Make sure the USB cable is securely connected between the facial brush and USB wall adapter, and the USB wall adapter is properly connected to a mains socket.

11 Technical data

Facial brush

Model:	ND-FC04
Rated voltage:	5 V
Rated power:	3 W
Battery capacity:	500 mAh
IP code:	IPX7 (protected against water immersion for up to 30 minutes at a depth of 1 m)
Protection class:	III
Speed settings:	7
Modes:	Vibrating and pulsating
Material:	ABS and silicone
Weight:	106 g
Dimensions (L × W × H):	110 × 35 × 70 mm

USB wall adapter

Power input:	100–240 V~ 50/60 Hz, 0.25 A
Power output:	5 V === 0.5 A 2.5 W
Protection class:	II

Model number: PS06H050K0500EU (EU)
PS06H050K0500BU (GB)

IP code: IP20

Published information	Value and precision	Unit
Name or trademark of the manufacturer, commercial register number and address:	Shenzhen Flypower Technology Co., Ltd. Building A2, Haosan Linpokeng 2nd Industry Zone, Nanpu Road, Shajing, Bao'an, Shenzhen, Guangdong, China	
Model identifier:	PS06H050K0500EU (EU) PS06H050K0500BU (GB)	
Input voltage:	100–240	V
Input alternating current frequency:	50/60	Hz
Output voltage:	5	V
Output current:	0.5	A
Output power:	2.5	W
Average efficiency during operation:	73.16	%
Efficiency at low load (10 %):	28.5	%
Power consumption at zero load:	0.1	W

12 Disposal

12.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

12.2 Disposing of the facial brush

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

Old appliances may not be disposed of in the household waste!

Should the facial brush no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in**

your city or county. This ensures that old appliances are recycled in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

Batteries and rechargeable batteries may not be disposed of with household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances* or not, to a collection point run by the communal authority or to a retailer, so that they can be disposed of in an

environmentally friendly manner.

Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

13 Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

14 Customer service

For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: www.silkn.eu.

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Copyright 2023 © Silkn Beauty Ltd. Todos os direitos reservados.
Copyright 2023 © Silkn Beauty Ltd. Με την επιφύλαξη παντός δικαιώματος.
Copyright 2023 © Silkn Beauty Ltd. Tutti i diritti riservati.
Copyright 2023 © Silk'n Beauty Ltd. Alle rettigheder forbeholdes.
© 2023 Silk'n Beauty Ltd. All rights reserved. Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.
Silk'n UK Ltd.