Silk'n®

LED FACE MASK PRO

FAC08

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
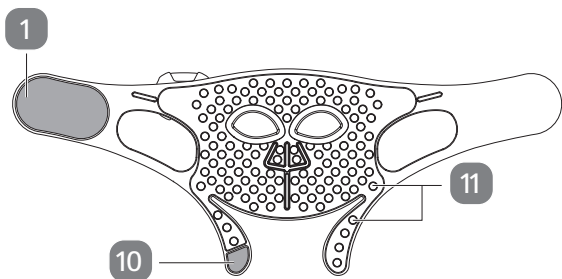
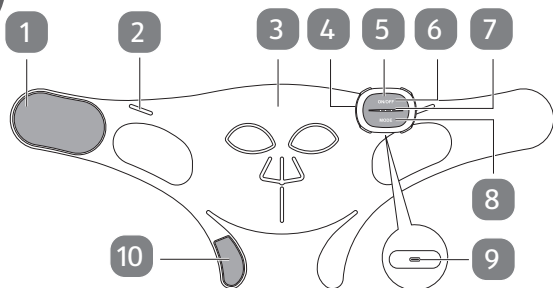
B

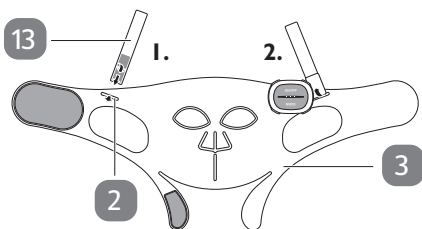Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silk'n Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silk'n Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silk'n Beauty Ltd.

Invention Works B.V.
Galileïstraat 17, 3029 AL Rotterdam, The Netherlands
www.silkn.eu – info@silkn.eu

Table of contents

I	Package contents/device parts	5
2	General information	5
2.1	Reading and storing the user manual	5
2.2	Explanation of symbols	5
3	Safety	6
3.1	Proper use	6
3.2	Safety instructions	6
4	Checking the product and package contents	10
5	Charging	10
6	Assembly	11
7	Operation	11
7.1	Preparation	11
7.2	Carrying out the session	12
7.3	Selecting the operation modes	12
8	Regimen	13
9	Cleaning	14
10	Storage	15
11	Troubleshooting	15
12	Technical data	15
13	Disposal	16
13.1	Disposing of the packaging	16
13.2	Disposing of the product	16
14	Warranty	17
15	Customer service	17

I Package contents/device parts

- | | |
|------------------------|-----------------------|
| 1 Head strap, 2× | 8 Mode button |
| 2 Small opening, 2× | 9 Charging port |
| 3 LED face mask | 10 Chin strap, 2× |
| 4 Magnetic holder | 11 LED lights |
| 5 Controller | 12 USB charging cable |
| 6 On/off button | 13 Elastic strap, 2× |
| 7 Indicator lights, 3× | 14 Storage pouch |

2 General information

2.1 Reading and storing the user manual

This user manual accompanies this LED Face Mask Pro (hereafter referred to as the “product”), and contains important information on setup and handling.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product. The user manual is based on the standards and rules in force in the European Union and United Kingdom. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

WARNING!

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

NOTICE!

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

Devices marked with this symbol are only suitable for household use (dry indoor rooms).

Direct current.

3 Safety

3.1 Proper use

The product is exclusively designed to help improve your skin's glow, elasticity and texture, and reduce the visibility of wrinkles, fine lines and pigmentation. It is only intended for private use and not suitable for commercial purposes. Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer cannot be held liable for damages or injury incurred through improper or incorrect use.

3.2 Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

- Only use the product with the supplied USB charging cable.
- The product must only be supplied at safety extra-low voltage corresponding to the marking on the product.
- Do not use the product if it is damaged or if the USB charging cable is defective.
- Do not open the housing; instead, have a qualified professional perform repairs.

Contact the customer service team for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection or incorrect operation.

- Only parts that comply with the original device data may be used for repairs. This product contains electrical and mechanical parts which are essential for providing protection against sources of danger.
- Do not immerse the product or USB charging cable in water or other liquids. Keep them away from water.
- Always switch the product off when you are not using it, when you intend to clean it, or in the event of a malfunction. Always unplug the USB charging cable when you are not charging the product.
- Never touch the USB charging cable with wet or damp hands.
- Never insert any objects into the housing.
- Keep all parts of the product away from open flames and hot surfaces.
- Lay the USB charging cable so that it does not pose a tripping hazard.
- Do not pull the USB charging cable out of the USB power source by the cord; instead, always pull it out by the cable plug.
- Do not kink the USB charging cable and do not lay it over sharp edges.

- Only use the product indoors.
- Do not use the product while charging.

WARNING! RISK OF INJURY!

- The product cannot be used by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge unless they have been given supervision or instruction concerning the use of the product in a safe way and understand the hazards involved. Children must not play with the product. Cleaning and user maintenance must not be performed by children.
- Keep children away from the product and USB charging cable.
- The product cannot be used by pregnant women, children, people with epilepsy, thyroid patients, people using thyroid medications, people who are photoallergic or use light allergy medication.
- The product cannot be used on skin diseases, wound ulceration, trauma and post-operative wounds.
- Do not use the modes that include Blue LED lights (modes F2 and F6) for more than 10 minutes per session, and max. 3 sessions per week. All other modes can be used on a daily basis starting from the fourth week.

- Consult your physician if you are uncertain whether you can use the product.
- Stop using the product immediately if you feel any discomfort, develop any abnormal reactions (such as itching or uncomfortable inflammation) or in case your skin becomes warm to the point of discomfort.
- Under abusive conditions, liquid may be ejected from the battery; avoid contact. Liquid ejected from the battery may cause irritation or burns.
- If battery fluid comes into contact with skin, eyes or mucous membrane, wash and rinse the affected areas immediately with clean water. Consult a physician if necessary.
- The battery must not be disassembled, short-circuited or thrown into a fire.
- Do not try to remove the battery from the product. This product contains batteries that are non-replaceable.
- Do not look directly at the light source part of the product.

NOTICE! RISK OF DAMAGE!

- Keep the product away from open fire, hot surfaces, high temperatures or the effects of the weather (e.g. rain).
- Stop using the product if its parts exhibit cracks or breaks or are deformed. Only

replace damaged parts with corresponding original spare parts.

4 Checking the product and package contents

NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.
1. Take the product out of the packaging.
 2. Check to make sure that all parts are included (see **Fig. A**).
 3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

5 Charging

- Fully charge the product before first use and when the battery power is low. The charging time is around 3 hours.
 - When the battery power is low, the product emits 5 beep sounds and the indicator lights flash red before the product automatically switches off. Recharge the product.
 - If the product is to be stored for an extended period, the battery level should be checked and the product should be charged regularly.
 - The product cannot be used during charging.
1. Remove the controller **5** from the magnetic holder **4**.
 2. Plug the smaller end of the USB charging cable **12** into the charging port **9** of the controller.
 3. Plug the larger end of the USB charging cable into the USB port of a compatible USB power source (not included) (e.g. USB wall adapter, 5 V DC, 1 A).

4. If a USB wall adapter is used, plug it into a mains socket.
The indicator lights **7** flash green during charging. Charging is done when the indicator lights illuminate green continuously.
If the voltage of the USB power source is too high, the indicator lights flash red rapidly and charging will not start. Use a suitable USB power source.
5. When charging is done, unplug the USB wall adaptor from the mains socket and disconnect the USB charging cable from the charging port.
6. Attach the controller back to the magnetic holder.
Make sure the controller is properly placed in the magnetic holder. Otherwise, the product will not operate.

6 Assembly

1. Open the hook-and-loop fasteners on the elastic straps **13**.
2. Guide the elastic straps through the small openings **2** (see step 1 in **Fig. B**).
3. Close the hook-and-loop fasteners to secure the elastic straps on the LED face mask **3** (see step 2 in **Fig. B**).

7 Operation

WARNING! RISK OF INJURY!

- Do not use the modes that include Blue LED lights (modes F2 and F6) for more than 10 minutes per session, and max. 3 sessions per week. All other modes can be used on a daily basis starting from the fourth week.

7.1 Preparation

1. Fully charge the product and attach the elastic straps **13** (see chapter "Charging" and "Assembly").
2. Thoroughly clean and dry your face. Make sure it is free from any make-up and dirt as they can reduce the amount of light entering the skin.
3. Sit or lie down and fit the LED face mask **3** comfortably on your face.

4. Use the hook-and-loop fasteners on the head straps **11**, chin straps **10** and elastic straps to secure the mask on your face. The elastic straps should be closed at the top of your head.

7.2 Carrying out the session

1. Press and hold the on/off button **6** to switch on the product. The product emits a beep sound and the indicator lights **7** illuminate white when switched on. The LED lights **11** turn on. By default, the "F1" operation mode is activated.
2. If desired, press the mode button **8** repeatedly to select another operation mode (F1–F7). The product emits a short beep sound when the operation mode is changed. See chapter "Selecting the operation modes" for more information.
3. The product automatically switches off after 10 minutes. You can also press and hold the on/off button to switch off the product at any time. The product emits 2 short beep sounds when switched off.
4. After the session, remove the mask and apply a serum and/or day/night cream on your face.
5. Thoroughly clean the product after use (see chapter "Cleaning").

7.3 Selecting the operation modes

- Press the mode button **8** repeatedly to cycle through the operation modes (F1–F7).
- The operation timer (10 minutes) is reset whenever you change to another operation mode.
- Description of the operation modes:

Mode	LED	Best for
F1	Red	Reduces fine lines, stimulates collagen production and helps to improve blood circulation.
F2	Blue	Reduces acne-related redness and imperfections for clearer skin.
F3	Amber	Improves skin hydration, soothes dry and sensitive skin, and helps reduce pigmentation and redness.

Mode	LED	Best for
F4	Near-Infrared + Deep Near-Infrared*	Helps reduce wrinkles, promote collagen and elastin production, and boosts blood circulation on a deeper level to improve skin tone, texture and resilience. It also helps maintain skin layers beyond the dermis.
F5	Red + Near-Infrared + Deep Near-Infrared (a red colour is visible)	Combines the benefits of F1 and F4: For younger-looking and tighter skin.
F6	Blue + Red + Near-Infrared + Deep Near-Infrared (a purple colour is visible)	Combines the benefits of F1, F2 and F4: For anti-ageing and the reduction of acne-related redness and imperfections to improve the skin's vitality, radiance and clarity.
F7	Amber + Red + Near-Infrared + Deep Near-Infrared (a dark orange colour is visible)	Combines the benefits of F1, F3 and F4: For calmer, younger-looking and tighter skin, and a more even skin tone.

*In mode F4, the LED lights  do not illuminate because Near-Infrared and Deep Near-Infrared are invisible lights. However, to indicate this mode is activated, 3 red LED lights in the mask illuminate: one on the forehead, one on the left cheek and one on the right cheek.

8 Regimen

WARNING! RISK OF INJURY!

- Do not use the modes that include Blue LED lights (modes F2 and F6) for more than 10 minutes per session, and max. 3 sessions

per week. All other modes can be used on a daily basis starting from the fourth week.

If you are using the product for the first time, follow the schedule below:

First week	Use once every 3 days, max. 10 minutes at a time
Second week	Use once every 2 days, max. 10 minutes at a time
Third week	Use once every 2 days, max. 10 minutes at a time
Fourth week	Use once every day, max. 20 minutes at a time

As soon as results show, adjust the frequency to once or twice per week.

9 Cleaning

WARNING! RISK OF ELECTRIC SHOCK!

- Always unplug the product from the power source before cleaning.

NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners, brushes with metal bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers, and the like.
 - Do not use hot water for cleaning.
1. Remove the controller  from the magnetic holder .
 2. Use a damp cloth or a water-based wipe to wipe off the LED face mask . Then wipe it dry with a dry cloth.
 3. Attach the controller back onto the magnetic holder after cleaning.

10 Storage

All parts must be completely dry before being stored.

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product so that it is not accessible to children.
- For best protection, store the product in the storage pouch **14**.

11 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not operate.	• Make sure the product is properly charged. • Make sure the product is not connected to a USB power source. The product does not work while charging. • Make sure the controller 5 is properly placed in the magnetic holder 4.
The product does not charge.	• Make sure the USB charging cable 12 is securely connected between the charging port 9 and the USB power source. • If a USB wall adapter is used for charging, make sure it is compatible with the product (5 V DC, 1 A).

12 Technical data

Model:	FAC08
Rated input:	5 V === 1 A
Working voltage:	3.7 V
Battery capacity:	2,000 mAh (lithium)
Number of LED lights:	120
IP code:	IPX5

Weight:	315 g
Dimensions (L × W × H):	254.7 × 674.8 × 44.7 mm
Material:	Silicone (mask) PC (remote control)

13 Disposal

13.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

13.2 Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

Old appliances may not be disposed of in the household waste!

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance with the regulations in force in your city or county.** This ensures that old appliances are recycled

in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

Batteries and rechargeable batteries may not be disposed of with household waste!

As the end user you are required by law to bring all batteries and rechargeable batteries, regardless whether they contain harmful substances* or not, to a collection point run by the communal authority or to a retailer, so that they can be disposed of in an

environmentally friendly manner.

Bring the entire product (including the rechargeable battery) to your collection point and make sure that it is in an uncharged state!

* labelled with: Cd = cadmium, Hg = mercury, Pb = lead

Note: To remove the rechargeable battery, wrap a towel around the controller and use a hammer to smash the housing. Make sure your eyes and hands are protected to avoid injury. Please note, the removal of the rechargeable battery is not reversible.

14 Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

15 Customer service

For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: www.silkn.eu.

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Silk'n UK Ltd.