Silk'n®

ULTRASOUND/ ELECTROPORATION ATTACHMENT HA2513

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
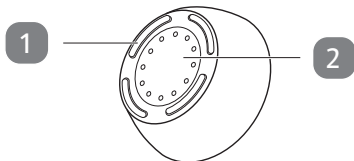
B

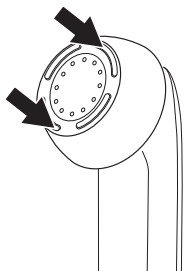Legal notice

Invention Works B.V. reserves the right to make changes to its products or specifications to improve performance, reliability, or manufacturability. Information furnished by Invention Works B.V. is believed to be accurate and reliable at the time of publication. However, Invention Works B.V. assumes no responsibility for its use. No license is granted by its implication or otherwise under any patent or patent rights of Invention Works B.V.

No part of this document may be reproduced or transmitted in any form or by any means, electronic or mechanical, for any purpose, without the express written permission of Invention Works B.V. Data is subject to change without notification.

Silkn Beauty Ltd. has patents and pending patent applications, trademarks, copyrights, or other intellectual property rights covering subject matter in this document. The furnishing of this document does not give you any license to these patents, trademarks, copyrights, or other intellectual property rights except as expressly provided in any written agreement from Silkn Beauty Ltd. Specifications are subject to change without notice.

Silk'n and the Silk'n logo are registered trademarks of Silkn Beauty Ltd.

Invention Works B.V.

Galileïstraat 17, 3029 AL Rotterdam, The Netherlands

www.silkn.eu – info@silkn.eu

Table of contents

I	Package contents/device parts	5
2	General information	5
2.1	Reading and storing the user manual	5
2.2	Explanation of symbols	5
3	Safety	6
3.1	Proper use	6
3.2	Safety instructions	6
3.3	Contraindications	8
4	Checking the product and package contents	10
5	Preparing for use	10
6	Operation	11
7	After use	11
8	Regimen	12
9	Cleaning	12
10	Storage	12
11	Troubleshooting	13
12	Technical data	13
13	Disposal	13
13.1	Disposing of the packaging	13
13.2	Disposing of the product	14
14	Warranty	14
15	Customer service	14

I Package contents/device parts

- 1 Electroporation electrodes, 4× 2 Ultrasound plate

2 General information

2.1 Reading and storing the user manual

This user manual accompanies this Ultrasound/ Electroporation attachment (hereafter referred to as the “product”), and contains important information on setup and handling. The product must only be used with the FaceTite MultiPlatform (H2501) device.

Before using the product, read the user manual carefully. This particularly applies to the safety instructions. Failure to do so may result in personal injury or damage to the product. The user manual is based on the standards and rules in force in the European Union and United Kingdom. When abroad, you must also observe country-specific guidelines and laws. Store the user manual for further use. Make sure to include this user manual when passing the product on to third parties.

2.2 Explanation of symbols

The following symbols and signal words are used in this user manual, on the product or packaging.

WARNING!

This signal symbol/word designates a hazard with moderate risk, which may result in death or severe injury if not avoided.

NOTICE!

This signal word warns of possible damage to property.

This signal symbol provides you with useful additional information on handling and use.

Read the user manual.

Declaration of Conformity: Products labelled with this symbol meet all applicable provisions of the European Economic Area.

3 Safety

3.1 Proper use

The product is exclusively designed to promote the absorption of your skin care products. It is only intended for private use and not suitable for commercial purposes.

Only use the product as described in this user manual. Any other use is considered improper and may result in damage to property or persons. The manufacturer cannot be held liable for damages or injury incurred through improper or incorrect use.

3.2 Safety instructions

WARNING! RISK OF ELECTRIC SHOCK!

- Only use the product with the FaceTite MultiPlatform (H250I) device.
- Carefully read the user manual of the FaceTite MultiPlatform (H250I) device before using the product.
- Do not use the product for other purposes.
- Do not open the housing; instead, have a qualified professional perform repairs. Contact the customer service team for this. Liability and warranty claims are waived in the event of repairs performed by the user, improper connection or incorrect operation.
- Only parts that comply with the original device data may be used for repairs. This product contains electrical and mechanical parts which are essential for providing

protection against sources of danger.

- Do not immerse the product in water or other liquids. Keep it away from water.
- Never insert any objects into the housing.
- Keep all parts of the product away from open flames and hot surfaces.
- Only use the product indoors.

WARNING! RISK OF INJURY!

- This product is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or instruction concerning use of the product by a person responsible for their safety.
- Keep children away from the product. Children should be supervised to ensure that they do not play with the product.
- Cleaning and user maintenance shall not be made by children without supervision.
- Make sure that children do not play with the plastic wrapping. They may get caught in it when playing and suffocate.

NOTICE! RISK OF DAMAGE!

- Keep the product away from open fire, hot surfaces, high temperatures or the effects of the weather (e.g. rain).

- Stop using the product if its parts exhibit cracks or breaks or are deformed. Only replace damaged parts with corresponding original spare parts.

3.3 Contraindications

Do not use the product if you:

- Have a pacemaker, internal defibrillator, or any other active electrical implant anywhere in the body.
- Suffer from poorly controlled endocrine disorders, such as diabetes and thyroid conditions.
- Have a history of skin cancer, or current condition of any other type of cancer or pre-malignant moles.
- Have received radiation therapy or chemotherapy treatments within the past 3 months, or have areas of potential skin malignancies.
- Have an impaired immune system due to immuno-suppressive diseases such as AIDS and HIV, or are using immuno-suppressive medications.
- Have used Isotretinoin (Accutane®) within 6 months prior to using this product.
- Have a disease related to photosensitivity, such as porphyria, polymorphic light eruption, solar urticaria, lupus, etc.

- Suffer from epilepsy.
- Are pregnant or nursing.

Do not use on areas that have:

- A permanent implant such as metal plates and screws (except oral ones), silicone implants, or areas with injected fills or toxins.
- Any active condition in the areas using this product, such as sores, psoriasis, eczema, rashes, open lacerations, abrasions, or lesions.
- purulent skin with inflammation, or history of keloid scars, abnormal wound healing or skin vulnerability.
- Any piercings or other metallic (e.g. titanium) materials.
- Had any surgery within 3 months prior to using this product.

Seek medical advice before continued use if you experience an unexpected side effect, such as itchy skin, within the first 48 hours of using this product.

Consult your physician if you have:

- Couperose or rosacea.
- A history of diseases stimulated by heat, such as recurrent herpes simplex in the areas using this product. If so, you may use this product following a prophylactic regimen.

- Used physician-prescribed drugs on your skin within the past 6 months.
- Skin diseases, keloids, or abnormal wound healing.
- Any other conditions which, in the opinion of your physician or dermatologist, would make it unsafe to use this product.

4 Checking the product and package contents

NOTICE! RISK OF DAMAGE!

- If you are not cautious when opening the packaging with a sharp knife or other pointed objects, you may quickly damage the product. Be very careful when opening it.
1. Take the product out of the packaging.
 2. Check to make sure that all parts are included (see **Fig. A**).
 3. Check whether the product or individual parts are damaged. If this is the case, do not use the product. Contact our customer service.

5 Preparing for use

1. For installation, gently press the product onto the FaceTite MultiPlatform device until it locks into place (see **Fig. B**).
2. Thoroughly clean and dry your skin. Make sure it is free from any make-up or dirt.
3. Apply your serum or day/night cream on your face. Make sure you use a serum/cream that is suitable to be used in combination with ultrasound or electroporation technology.

6 Operation

WARNING! RISK OF INJURY!

- Do not use the product for more than 5 minutes per session.
- The product is not operating when connected to the power supply. It can only be used cordless.
- Refer to the user manual of the FaceTite Multiplatform device to turn on/off and select the energy levels (1–4). The electroporation electrodes  and ultrasound plate  activate once the device is turned on.
- Move the device over the skin in circular movements during use.
- Always start using the product from energy level 1 and increase the energy level as long as you feel comfortable.
- Use the highest energy level for the best results. However, it is important that you feel comfortable with the selected energy level. Lower the energy level once you feel any discomfort.
- For each use, use the product for at least 5 minutes on each large area (e.g. cheeks, neck) and at least 3 minutes on each small area (e.g. nasolabial folds (smile lines) and areas between the eyebrows).

Do not leave the product on one spot of your skin. Keep moving the product to cover the entire area. This prevents overheating and reduces the possibility of any adverse effects.

7 After use

1. When done, turn off the device and remove the product from the device.
For removal, gently pull the product outwards until it detaches from the device.
2. Thoroughly clean the product after use (see chapter “Cleaning”).

8 Regimen

- Always start using the product from energy level 1 and increase the energy level as long as you feel comfortable.
- For each use, use the product for at least 5 minutes on each large area (e.g. cheeks, neck) and at least 3 minutes on each small area (e.g. nasolabial folds (smile lines) and areas between the eyebrows).
- Use the product twice a week for 10 weeks. You will begin to see visible results after 5 weeks/10 uses.
- After that, use once every 2 weeks to maintain results.
- To prevent temporary results, it is important to use the product consistently.

9 Cleaning

NOTICE! RISK OF DAMAGE!

- Do not use any aggressive cleaners, brushes with metal bristles, or sharp or metallic cleaning utensils such as knives, hard scrapers, and the like.
 - Do not use hot water for cleaning.
1. Clean the product with a damp cloth. Make sure there is no remaining skin care product on the product.
 2. Dry the product thoroughly with a dry, soft and clean cloth.

10 Storage

All parts must be completely dry before being stored.

- Store the product in a clean and dry place, avoiding heat and direct sunlight.
- Store the product so that it is not accessible to children.

11 Troubleshooting

Some problems may be caused by minor faults that you can fix yourself. To do so, follow the instructions in the following table. If it is still not possible to resolve the problem, contact customer service. Do not repair the product yourself.

Fault	Possible cause and solution
The product does not work.	Make sure the product has been properly installed onto the FaceTite MultiPlatform device. When properly installed, the battery/attachment indicator light of the device illuminates green.
Energy level indicator light 1 and 4 are blinking.	Make sure the product is not connected to the power supply in the mains socket. The product can only be operated cordless.

12 Technical data

Model:	HA25I3
Frequency	
Ultrasonic:	3 ± 0.5 MHz
Electroporation:	35 / 30 / 25 / 20 KHz
Weight:	86 ± 10 g
Dimensions (L × W × H):	$66.5 \times 66.5 \times 31$ mm

13 Disposal

13.1 Disposing of the packaging

Sort the packaging before you dispose of it. Dispose of paperboard and cardboard with the recycled paper service and wrappings with the appropriate collection service.

13.2 Disposing of the product

(Applicable in the European Union and other European countries with separate collection systems of recyclable materials)

Old appliances may not be disposed of in the household waste!

Should the product no longer be capable of being used at some point in time, **dispose of it in accordance**

with the regulations in force in your city or county. This ensures that old appliances are recycled

in a professional manner and also rules out negative consequences for the environment. For this reason, electrical equipment is marked with the symbol shown here.

14 Warranty

This product is covered by a 2-year warranty according to European regulations and law. The extent of the warranty on this product is limited to technical defects caused by faulty production processes. In the event you want to claim warranty please be sure to contact our customer service for instructions. They may be able to solve your problem without needing to return the product to the store or our service center. Our customer service is always pleased to help you!

15 Customer service

For more information about the products of Silk'n visit your regional Silk'n website: www.silkn.eu. If the device is damaged, defective, in need of repair or you need our assistance, please contact your nearest Silk'n service center. This manual is also available as a PDF download from: www.silkn.eu.

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Silk'n UK Ltd.